



GEMSTAR[®]
Infusion Pump
System Operating Manual

**For use with the
Hospira GemStar[®] Pump
List 13000-04-07/08
List 13100-04-07/08**

Hospira, Inc., Lake Forest, IL 60045
430-00897-A02 (Rev. 02/05)

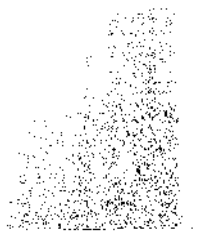

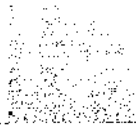

0
1
2
3
4
5
6
7
8
9
10
11
12
13
14
15
16
17
18
19
20
21
22
23
24
25
26
27
28
29
30
31
32
33
34
35
36
37
38
39
40
41
42
43
44
45
46
47
48
49
50
51
52
53
54
55
56
57
58
59
60
61
62
63
64
65
66
67
68
69
70
71
72
73
74
75
76
77
78
79
80
81
82
83
84
85
86
87
88
89
90
91
92
93
94
95
96
97
98
99

Hospira GemStar® Infusion Pump System Operating Manual

**For use with the
Hospira GemStar Pump**

List 13000-04-07/08

List 13100-04-07/08

**PLEASE
Read this entire manual before
using the Hospira GemStar Pump**

This manual is designed for use by healthcare professionals, caregivers, and patients. The Hospira Technical Support Operations hotline is available 24 hours a day (in the USA) to provide consultation and technical assistance regarding the Hospira GemStar Pump.

In the USA: **Hospira Technical Support Operations
1-800-241-4002**

Outside the USA: **Contact your local Hospira
sales office**

**To order additional copies of this manual
(List 13981) call Hospira Customer Care:
1-877-9HOSPIRA
(1-877-946-7747)**

Change History

Part Number	Description of Change	Pages Changed
430-00897-001 (Rev. 9/02)	Original issue	N/A
430-00897-A01 (Rev. 12/03)	Updated Cover	Front and Back
	Updated Change History	ii
	Updated Battery Pack Information	19
430-00897-002 (Rev. 10/04)	Second Release	Inclusive
430-00897-A02 (Rev. 02/05)	Updated Cover	Front and Back
	Updated Change History	ii
	Aligned Note and Tip Icons	21, 61, 69, 98, 97, 100–102, 106–107, 111, 122, 125, 144, 149, 178, 190
	Corrected typographical error in note	178
	Reformatted dashed line with scissors icon	121–126
	Added UL/CSA information	171
	Updated CSA information	Back

Note: Change page manuals assembled by Hospira include a change page identifier in the part number on the cover page and on the bottom of each changed page. This change page manual is identified as 430-00897-A02.

Contents

1. Introduction	1
GemStar Pump Layout	2
GemStar Pump Components	4
Therapies	5
GemStar Pump Features	6
Menu System	7
Operating Modes	9
Indications for Use	10
Contraindications for Use	11
Warnings and Cautions	11
2. Basic Operation	17
Overview	17
Selecting a Power Source	18
Connecting to the GemStar Docking Station	23
Powering On the Pump	24
Accessing the Programming Menu	26
Programming the Pump	27
Starting a New Container	29
GemStar Pump Set	29
Loading the Cassette	32
Releasing the Cassette	33
Priming a GemStar Pump Set	34
3. Total Parenteral Nutrition (TPN)	47
Therapy Features and Specifications	47
TPN Programming	48
Changing a TPN Program	49
TPN Auto Taper Down	50
TPN Programming Worksheet	51

Contents

4. Pain Management	53
Therapy Features and Specifications	53
Pain Management Programming	54
Bolus Delivery	56
Loading Dose Delivery	57
Changing a Pain Management Program	58
Pain Management Programming Worksheet	60
5. Intermittent	61
Therapy Features and Specifications	61
Intermittent Programming	62
Intermittent Delivery Interruptions	64
Changing an Intermittent Program	66
Intermittent Programming Worksheet	67
6. Continuous	69
Therapy Features and Specifications	69
Continuous Programming	70
Rate Titration	71
Piggybacking	71
Changing a Continuous Program	74
Continuous Programming Worksheet	75
7. Weight Dosed	77
Therapy Features and Specifications	77
Weight Dosed Programming	78
Auto KVO	79
Rate Titration	79
VTBI Titration	80
Changing a Weight Dosed Program	80
Weight Dosed Programming Worksheet	81

Contents

8. mL/hr Only	83
Therapy Features and Specifications	83
mL/hr Only Programming	84
Auto KVO	85
Rate Titration	85
VTBI Titration	86
Piggybacking	86
Changing a mL/hr Only Program	89
mL/hr Only Programming Worksheet	90
9. Variable Time	91
Therapy Features and Specifications	91
Variable Time Programming	92
Variable Time Delivery Interruptions	93
Changing a Variable Time Program	94
Variable Time Programming Worksheet	95
10. Options Menu	97
Overview	97
Tips for Using the Options Menu	97
Review Program	99
Histories	99
Keypad Lock	104
Change Pump Settings	106
Set Clock	108
Print	109
Speed Protocol	110
Next Dose	111
Connecting to the Data Port	112
Transferring Data Using a Computer	117
11. Clinician Instructions	121
Clinician Locking Sequence Instructions	122
Program and Deliver a Clinician Activated Loading Dose	124
Shift Totals and History	125

Contents

12. Troubleshooting	127
Alarm Messages	128
Alert Messages	135
Other Display Messages	138
13. Maintenance	139
Pump Storage	139
Cleaning and Disinfecting	139
Repair	141
14. Operation Test	143
Equipment Required	143
Test Setup	144
Performing the Operation Test	144
Printing Test Results	149
Operation Test Checklist	151
15. Optional System Components	153
Connecting the Bolus Cord	156
Using the GemStar Lockboxes	156
Using the Hospira Carrying Cases	159
Using the GemStar Docking Station	165
16. Pump Sets	167
17. GemStar Pump Specifications	171
Occlusion Information	174
18. Warranty/Technical Assistance	185
Index	191
IEC Symbols	inside back cover

Chapter 1

Introduction

The GemStar® Pump is a small and lightweight, single-channel infusion device designed for use in the home, in the hospital, or anywhere electronic infusion is required. The GemStar Pump can be powered by AC mains adaptor, rechargeable battery pack, docking station, or two disposable AA alkaline or lithium batteries. When powered by batteries, the GemStar Pump is ideal for ambulatory patients.

The following conventions are used throughout this manual:

Convention	Use	Example
[ALL CAPS]	Pump keys	Press [YES/ENTER] to continue.
ALL CAPS	Display messages	The pump displays: ENTER TO CONFIRM.

WARNING: A WARNING message contains special safety information and must be observed at all times. Failure to observe a warning message could be life threatening.

CAUTION: A CAUTION contains information that could prevent irreversible product damage or hardware failure. Failure to observe a caution could result in serious patient or user injury.



Note: A NOTE contains information or examples to help explain instructions.




Tip: A TIP contains information that makes the pump easier to use.



Note: Examples are provided only to help explain a pump function and may not be clinically relevant. Graphics are for illustrative purposes only and may not reflect the actual product.

1. Introduction

GemStar Pump Layout

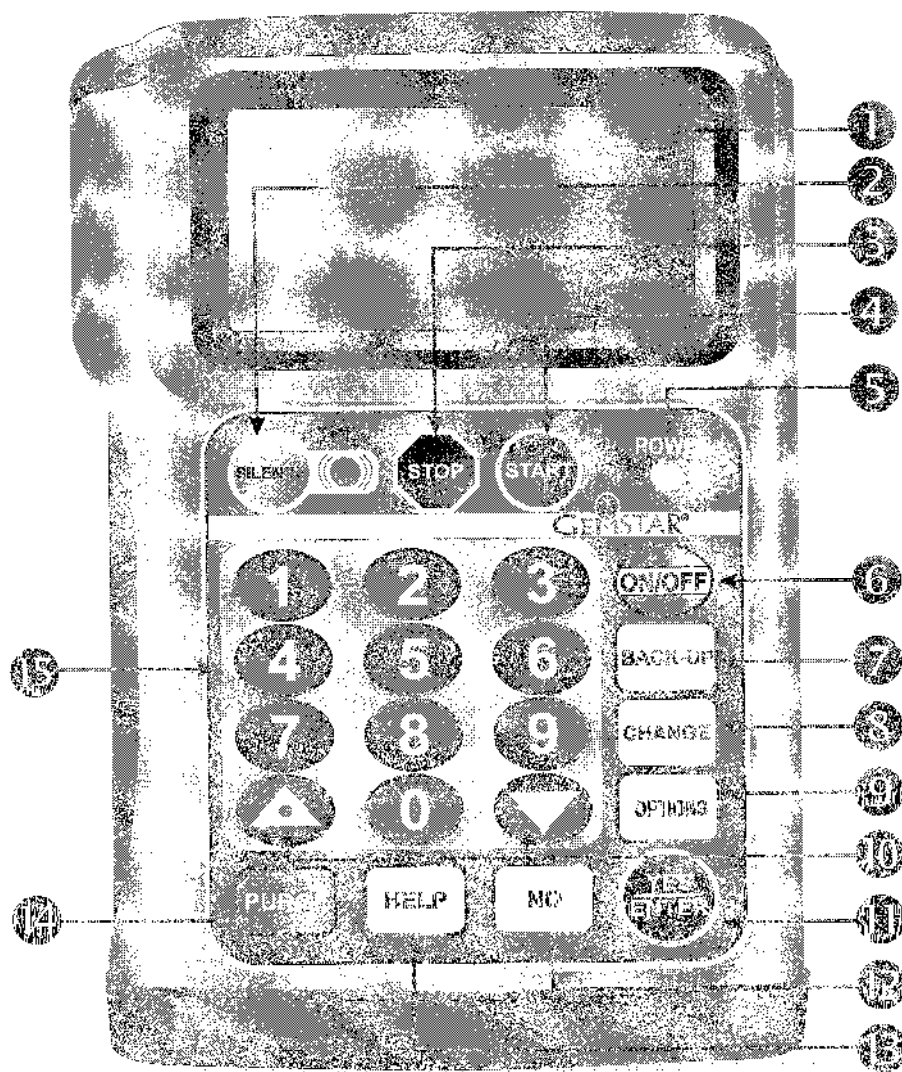

Shown approximate size


1. Introduction

1	DISPLAY	• Displays the pump's status.
2		• Silences an alarm.
3		• Stops infusion.
4		• Starts infusion.
5	 Power Indicator	• Lights continuously when the pump is powered by AC mains. • Flashes when the pump is powered by external batteries (e.g., battery pack). • Not lit when using disposable AA batteries.
6		• Powers the pump on or off. To power-off the pump, press and hold until the display clears.
7		• Returns to previous steps during programming. • Exits Help, Change, and Options menus.
8		• Displays the following menu items from the STOP or RUN mode*: - New Container Change Program (or Titrate VTBI) - New Program Piggybacking • Resets a numeric value to zero.
9		• Displays the following Options Menu items from the STOP or RUN mode*: 1 Review Program 5 Set Clock 2 Histories 6 Print 3 Keypad Lock 7 Speed Protocol 4 Pump Settings 8 Next Dose
10	 	• Scrolls through display messages. • Toggles between AM and PM in time-based entries. • UP arrow enters a decimal point in numeric entries. • DOWN arrow toggles between units of measure when entering concentrations.
11		• Responds "YES" to display questions. • Accepts entries and goes to the next step.
12		• Responds "NO" to display questions.
13		• Displays information during programming steps and alarm conditions. • Displays Shift Totals, except during programming steps and alarm conditions.
14		• Clears an AIR-IN-LINE alarm. • Primes a Gemstar Pump Set.
15	[0] - [9]	• Selects menu items and enters numeric values. If two numbers appear before a menu item (e.g., 02), press both number keys to select that item.

* Not all menu items are available with all modes.

1. Introduction

GemStar Pump Components

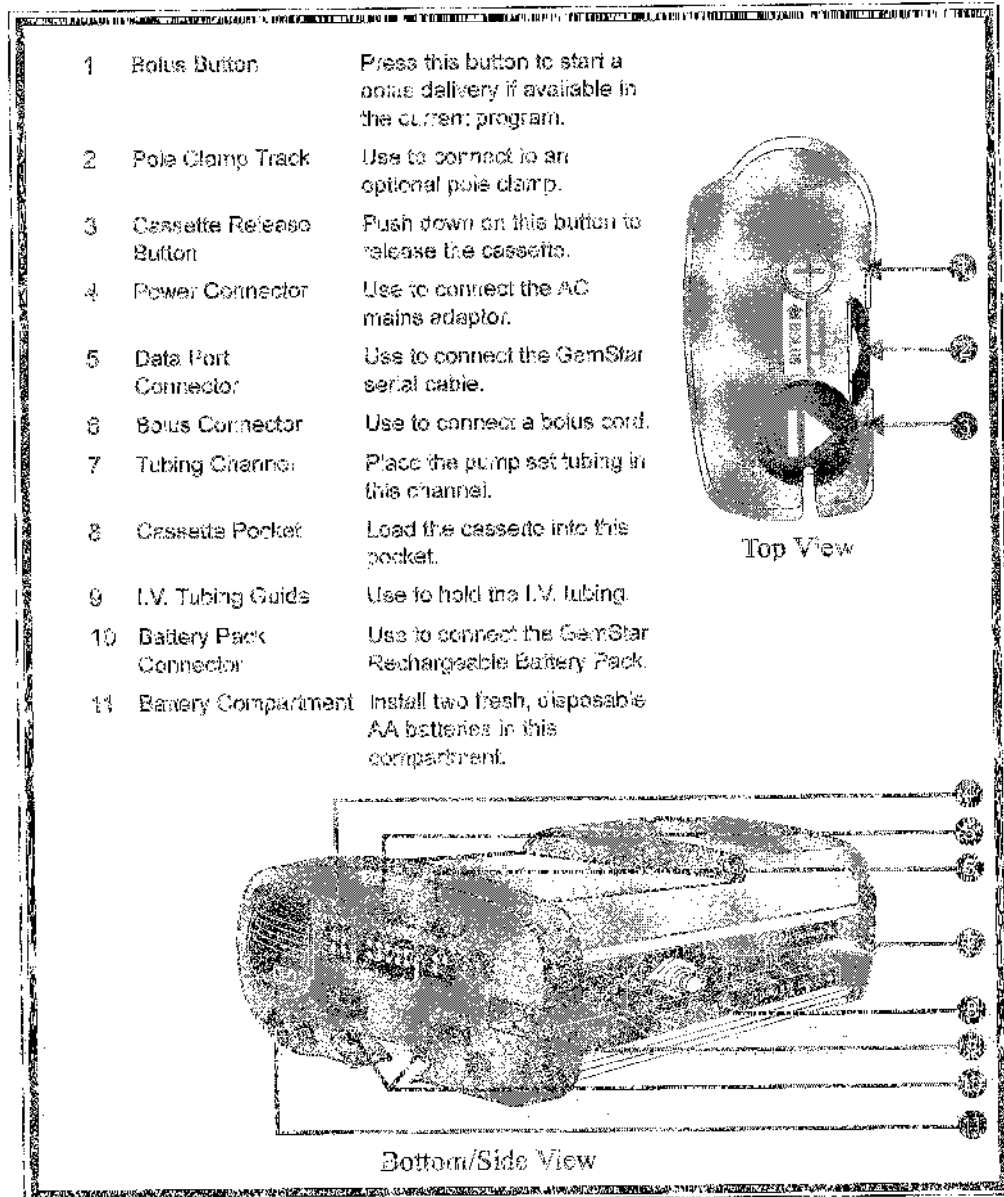

1. Introduction

Therapies

Seven types of therapy are available for infusion of medications and fluids. The following table lists possible applications and available features of each therapy.

Therapy	Possible Applications	AUTO KVO	TITRATE RATE	PIGGY-BACKING	BOLUS	LOADING DOSE	FASTER UP/DOWN	TITRATE VTEB	OPTIONAL KVO
TPN	Total Parenteral Nutrition fluid delivery						*		*
Pain Management	Infusion of analgesics or anesthetics for patient-controlled pain management (PCA)				*	*			
Intermittent	Interval-based therapies such as antibiotics								*
Continuous	Infusion of mg, mcg, or mL/hr with a selectable KVO (Keep Vein Open) option		*	*					*
Weight Dosed	Critical care, pediatrics	*	*					*	
mL/hr Only	Simple rate programming	*	*	*				*	
Variable Time	Chronotherapy				*				*
* = Available									

Note: Not all therapies are available in all pump configurations.

1. Introduction

GemStar Pump Features

Feature	Description	Reference Page
• <i>Adjustable Settings</i>	Adjusts the distal occlusion pressure limits, air sensor sensitivity, real-time clock, and audible alarm volume.	Chapter 10 Options Menu
• <i>Alerts and Alarms</i>	Signals audible and visual alarms when attention is required.	Chapter 12 Troubleshooting
• <i>Audible Keypad Response</i>	Sounds a single beep when any key is pressed. If an invalid key is pressed, the pump sounds a flutter tone.	Chapter 1 Introduction, <i>GemStar Pump Layout</i>
• <i>Keypad Lock</i>	Restricts access to various options. Four lock levels are available.	Chapter 10 Options Menu and Chapter 11 Clinician Instructions
• <i>Next Container</i>	Repeats the current program without reprogramming the pump.	Chapter 2 Basic Operation
• <i>Piggybacking</i>	Allows delivery from a secondary container during a Continuous or mL/hr Only program.	Chapter 6 Continuous and Chapter 8 mL/hr Only
• <i>Power-on Self-test</i>	Performs a self-test each time the pump is powered on.	Chapter 2 Basic Operation
• <i>Operation Test</i>	Performs user-assisted tests to verify the pump is operating properly.	Chapter 14 Operation Test
• <i>Program History</i>	Maintains a history event log, with a date and time stamp of each event for each program. Review the history on the display, or download it to a printer or computer.	Chapter 10 Options Menu
• <i>Speed Protocols</i>	Stores up to nine frequently used programs in the pump's memory for quick access.	Chapter 10 Options Menu

Menu System

The GemStar Pump uses a menu system that provides step-by-step instructions through all of the pump's functions. Refer to the following table for a list of functions in each menu.

Use the keypad to select menu items, respond to display messages, and enter numeric values. For example, press [2] to select PAIN MANAGEMENT from the Programming Menu.



Note: Menu items may vary depending on the pump's configuration.



Tip: A menu item does not have to be on the display to be selected. For example, when the Programming Menu first displays, only the first four therapies are shown. Instead of pressing the down arrow to see the other menu items, press [7] to go directly to Variable Time Programming.

1. Introduction

To display this menu	Do this	Menu Items Available
Programming Menu	1 Press CHANGE to display the Change menu. 2 Select NEW PROGRAM.	1 TPN 2 PAIN MANAGEMENT 3 INTERMITTENT 4 CONTINUOUS 5 WEIGHT DOSED 6 mL/HR ONLY 7 VARIABLE TIME <i>Note: If a program is not currently stored in pump memory at power-on, the Programming Menu automatically displays after the self-test.</i>
Change Menu	Press CHANGE .	1 NEW CONTAINER 2 NEW PROGRAM 3 CHANGE PROGRAM 4 PIGGYBACKING <i>Note: If the current therapy is Weight Dosed or mL/hr Only and the pump is in RUN mode:</i> 3 ELITERATE VEST will replace 3 CHANGE PROGRAM on the menu display. <i>Note: Piggybacking is available only in mL/hr Only and Continuous therapies.</i>
Options Menu	Press OPTIONS .	1 REVIEW PROGRAM 2 HISTORIES 3 KEYPAD LOCK 4 PUMP SETTINGS 5 SET CLOCK 6 PRINT 7 SPEED PROTOCOL 8 NEXT DOSE

1. Introduction

 For more information on:

Programming Menu

Change Menu

Options Menu

Refer to the following pages:

Chapter 2 Basic Operation, page 17.

For menu item [1], refer to *Starting a New Container* in Chapter 2 Basic Operation, on page 29. For menu items [2], [3], and [4], refer to the appropriate therapy chapter.

Chapter 10 Options Menu, page 97.

Operating Modes

The GemStar Pump has four operating modes:

Mode	Description
Help Mode	The pump enters the HELP mode when  is pressed during programming or an alarm condition. The HELP mode displays additional information about the current condition.
Programming Mode	The pump enters the PROGRAMMING mode after the pump is powered on or when  is pressed.
Run Mode	After  is pressed, the pump enters the RUN mode, even when not infusing.
Stop Mode	When programming is complete, the pump enters the STOP mode until  is pressed. When the programmed therapy is complete, and during certain alarm conditions, the pump automatically enters the STOP mode. Press  to enter the STOP mode at any time, except during a CHECK CASSETTE alarm.

1. Introduction

Indications for Use

Physicians or certified, licensed healthcare professionals should always oversee infusions. Healthcare professionals should always supervise pump users and instruct users in pump operation and troubleshooting. Instruction should include how to prevent intravenous (I.V.) related complications, including how to prevent accidental infusion of air.

The pump is suitable for intravenous, arterial, subcutaneous, short-term epidural infusion and parenteral administration of general I.V. fluids, medications, nutritional fluids, and blood/blood products. The epidural route is recommended to provide anesthesia or administer analgesia for periods up to 96 hours.

Note: Epidural administration of anesthetics is recommended in the continuous mode only. Epidural administration of analgesics may be delivered by continuous, bolus, or continuous with bolus.

For epidural use, the administration of drugs is restricted to those anesthetic and analgesic drugs approved for continuous epidural administration: e.g., Chloroprocaine Hydrochloride USP, Lidocaine Hydrochloride USP, and Morphine Sulfate Injection USP (Preservative Free).

WARNING: Delayed respiratory depression following continuous epidural administration of preservative-free morphine sulfate has been reported.

For epidural administration, the following is recommended:

- Nylon or Teflon® catheter
- Pump sets without Y-sites
- Sackers for the tubing indicating ongoing administration

Note: Facilities practicing epidural anesthesia/analgesia must be staffed and equipped to manage cardio-pulmonary resuscitation (CPR). Supplies should include oxygen, naloxone, and other appropriate resuscitative drugs and equipment. Monitor patients continuously (e.g., oximetry) during epidural administration. Frequently observe patient for side effects (for up to 24 hours) following completion of drug administration by the epidural route.

Contraindications for Use

Persons who do not have the mental and physical capability or the emotional stability to properly use the GemStar Pump should not be operate this pump.

DO NOT use the following with the GemStar Pump System:

- ♦ Drugs not compatible with silicone rubber or PVC plastic
- ♦ Drugs not stable under infusion conditions

Warnings and Cautions

Pay attention to the following warnings and cautions when operating the GemStar Pump. These warnings and cautions are repeated in the manual when relevant to the procedure being discussed. Pay attention to all alert messages.

General Warnings and Cautions

- ♦ Federal (USA) law restricts this device to sale by or on the order of a physician or other licensed practitioner.
- ♦ Manual references to specific values are approximate only, unless indicated otherwise. Air-in-line sensitivity values are approximate only.
- ♦ When a Proximal Occlusion alarm activates, DO NOT power-off the pump.
- ♦ Closely supervise patients who are likely to be adversely affected by unintended operations and failures, including interrupted medication or fluid infusion. Ensure provisions for immediate corrective action are available.

Regarding Drugs Used, Pump Sets, and Containers

- ♦ USE ONLY GemStar Pump Sets with the GemStar Pump. Use of unauthorized sets may result in injury to the patient or damage to the pump.
- ♦ NEVER use drugs that are incompatible with silicone rubber or PVC plastic.
- ♦ To reduce the loss of potency for drugs known to be absorbed by PVC plastic and silicone, begin infusion as soon as practical after priming the pump set. Use of high flow rates during infusion will minimize drug absorption.

1. Introduction

- DO NOT use medications that are unstable under infusion conditions.
- ALWAYS use connections with Luer Lock™ fittings.
- To prevent contamination, use aseptic technique with all fluid-path connections. Remove the protective coverings as assembly progresses.
- ALWAYS place the flow stop in the closed position when the set is not in the pump.
- Make sure the cassette is properly installed. When properly installed, all four cassette latches are visibly holding the cassette securely in the cassette pocket.
- Ensure the tubing is behind the tubing guide (at the distal end).
- Arrange tubing, cords, and cables to minimize the risk of patient strangulation or entanglement.
- NEVER use vented fluid containers (e.g., glass or rigid plastic) unless suspended from a pole. When using a vented drip set, ALWAYS use an air-eliminating filter OR set the air sensitivity to either ON or 2 mL.
- Stop infusion if signs or symptoms of infiltration occur.
- GemStar Pump Sets contain an anti-siphon valve to reduce the potential for gravity flow at a typical 36-inch head height.
- Failure to use the anti-siphon valve may result in unrestricted flow.
- ALWAYS close the slide clamps before removing the cassette from the pump.
- When using the pump for secondary delivery (piggybacking), ensure the fluids being infused are chemically and physically compatible.
- When infusing short-half-life drugs (i.e. those dosed with units of mcg/kg/min) at very low rates (0.8 mL/hr or less) the rate may be too slow for that drug. Consider using a lower drug concentration with these drugs so a pump rate above 0.8 mL/hr can be used.
- When using pre-filled vials, the system may underdeliver at low flow rates.

1. Introduction

Regarding Air-in-Line and Infusion

- ♦ To reduce the risk of infusing air, use an air-eliminating filter when the air sensitivity is set to OFF.
- ♦ ALWAYS prime the pump set to remove all air from the cassette, tubing, and injection site prior to connecting to the patient.
- ♦ ALWAYS disconnect the pump set from the patient prior to priming or purging.

Regarding Pump Operation

- ♦ If the pump does not perform as stated in this manual, remove the pump from service IMMEDIATELY.
- ♦ If the pump does not display the self-test or if beeps do not sound at power-on, check all power connections and install two fresh, disposable AA batteries. If the problem continues, contact Hospira Technical Support Operations before using the pump.
- ♦ Ensure a cassette is installed in the pump before connecting to a printer or computer. **If the pump is connected to a patient:** BEFORE connecting the pump to a computer or printer, press [STOP], close the slide clamp, and eject the cassette. Leave the pump set connected to the patient. Install another cassette in the pump before connecting to a printer or computer.
- ♦ Disconnect the pump from the patient BEFORE performing the Operation Test.
- ♦ ALWAYS connect to a grounded AC mains power source when using AC mains power. If the quality of the grounding source is in doubt, operate the pump only with batteries.
- ♦ USE ONLY the AC mains adaptors specifically labeled for use with the GemStar Pump.
- ♦ USE ONLY the AC mains adaptor specifically labeled for use with the GemStar Pump to charge the battery pack.
- ♦ The GemStar AC Mains Adaptor is for use with the GemStar Pump or GemStar Rechargeable Battery Pack only. DO NOT use the GemStar AC mains adaptor with other products.
- ♦ During charging, the battery pack is warm. If the battery pack becomes hot to the touch, IMMEDIATELY unplug the AC mains adaptor and contact Hospira Technical Support Operations.

1. Introduction

- When removing the cord from the bottom of the pump, firmly grasp the connector and pull straight out. DO NOT twist or bend the cord or connector.
- To ensure proper pump operation, ALWAYS replace BOTH disposable batteries with fresh, disposable AA batteries when a change is required.
- ALWAYS fully close the battery door when using disposable batteries.
- DO NOT use rechargeable batteries in the battery compartment.
- Ensure two fresh, disposable AA batteries are installed for backup power when using other power sources (e.g. docking station, battery pack). The pump will continue to operate on backup power if the external power source fails.
- ALWAYS avoid sources of high-intensity electromagnetic radiation (e.g., radio transmitters, MRI scanners, microwave ovens, X-ray machines, and CAT scanners).
- Use of radio-frequency emitting devices, such as cellular telephones, 2-way radios, and Electrical Surgical Devices (ESU), in close proximity to this device may affect its operation.
- Possible explosion hazard exists if the pump is used in the presence of flammable anesthetics. NEVER use the pump in the presence of flammable or explosive vapors.
- Non-hazardous, low-level electrical potentials are commonly observed when fluids are administered using infusion devices. These potentials are well within accepted safety standards, but may create artifacts on voltage-sensing equipment, such as ECG, EMG, and EEG machines. These artifacts vary at a rate that is associated with the infusion rate. If the monitoring machine is not operating correctly or has loose or defective connections to its sensing electrodes, these artifacts may be accentuated so as to simulate actual physiological signals. To determine if the abnormality in the monitoring equipment is caused by the infusion device instead of some other source in the environment, set the infusion device so that it is temporarily not delivering fluid (if clinically acceptable). Disappearance of the abnormality indicates that it was probably caused by electronic noise generated by the infusion device. Proper setup and maintenance of the monitoring equipment should eliminate the artifact. Refer to the appropriate monitoring system documentation for setup and maintenance instructions.

1. Introduction

- ♦ The GemStar pump, when connected to the GemStar AC adaptor, has been tested and found to comply with EMC/EMI limits in accordance with IEC/EN 60601-1-2 (2001). When connected to other power sources and/or optional system components, the GemStar system complies with EMC/EMI limits in accordance with IEC/EN 60601-1-2 (1993). These limits are designed to provide reasonable protection against harmful interference in a typical medical installation. The equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to other devices in the vicinity. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference with other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - ♦ Reorient or relocate the receiving device
 - ♦ Increase the separation between the equipment
 - ♦ Connect the equipment to an outlet on a circuit different from that to which the other device is connected
 - ♦ Consult the manufacturer or field service technician for help
 - ♦ The use of portable and mobile RF equipment may have an impact on this and other pieces of medical equipment

Regarding Handling and Maintenance

- ♦ Use proper care during unpacking, installation, and operation of the pump. If the pump is inadvertently mishandled, check connections and programmed data to ensure there is no damage. Refer to Chapter 14 Operation Test for further information.
- ♦ ALWAYS avoid dropping or hitting the pump. If the pump is dropped or hit, ALWAYS check connections and programmed data.
- ♦ NEVER use sharp objects (e.g., fingernails, pens, pencils, or other probes) to program or clean the pump.
- ♦ Use care not to damage the silicone seals around the sensor bodies in the cassette pocket.
- ♦ To avoid mechanical or electronic damage, NEVER submerge pump in water or other fluids and avoid fluid spills. If the pump becomes wet, dry it immediately with a dry, lint-free cloth. Check connections and programmed data.

1. introduction

- Some cleaning and sanitizing compounds may slowly degrade components made from some plastic materials. DO NOT use compounds containing combinations of isopropyl alcohol and dimethyl benzyl ammonium chloride (i.e., alcohol wipes).
- DO NOT sterilize by heat, steam, ethylene oxide (ETO), or radiation. Apply disinfectants to the outside surface of the pump only. DO NOT use abrasive cleaners or materials on the pump. Using abrasive cleaners or cleaning solutions not recommended by Hospira may result in product damage.
- Please comply with local disposal and recycling regulations as appropriate for disposable batteries, rechargeable battery packs, medical electronic components, and pump sets.

Chapter 2

Basic Operation

The minimum elements required for using the GemStar Pump are:

- ♦ GemStar Pump
- ♦ Appropriate power source
- ♦ Appropriate GemStar Pump Set
- ♦ Appropriate fluid container
- ♦ Patient-access device

Overview

To start a therapy:

- 1 Determine the program settings.** Refer to the appropriate therapy chapter for a worksheet that may be copied and used to record program settings.
- 2 Select a power source.** Refer to *Selecting a Power Source* on page 18.
- 3 Prime the pump set.** Refer to *GemStar Pump Set* on page 29.
- 4 Load the cassette into the pump.** Refer to *Loading the Cassette* on page 32.
- 5 Power-on the pump and go to the Programming Menu.** Refer to *Powering On the Pump* on page 24.
- 6 Select the appropriate optional system components, such as a bolus cord.**
- 7 Select the therapy type and enter a program.** Refer to the appropriate therapy chapter for detailed instructions.
- 8 Press [START] to begin delivery.**

2. Basic Operation

Selecting a Power Source

The GemStar Pump may be powered by one or more of the following:

Power Source	Description	Ambulatory	Hospital	Home
GemStar AC Mains Adaptor (USA)	Wall-mount adaptor plugs into a 110 VAC mains supply.		*	*
GemStar AC Mains Adaptor (Universal)	Tabletop adaptor plugs into a 100–240 VAC mains supply.		*	*
GemStar Docking Station	Plugs into a 100–240 VAC Volt mains supply while mounted on an I.V. pole. The docking station's internal batteries power the GemStar Pump during periods of transport.	*	*	*
GemStar Battery Pack	The battery pack's internal batteries power the GemStar Pump during periods of transport or when use of AC power is not desirable.	*	*	*
Disposable Batteries	Two disposable AA alkaline or lithium batteries (available from local sources)	*	*	*

CAUTION: To ensure proper pump operation, ALWAYS replace BOTH batteries with fresh, disposable AA batteries when a change is required.


DO NOT use rechargeable batteries in the battery compartment.

Ensure two fresh, disposable AA batteries are installed for backup power when using other power sources (e.g. docking station, battery pack). The pump will continue to operate on backup power if the external power source fails.

2. Basic Operation

Installing Disposable Batteries

- 1** Flip open the tab of the battery door on the bottom of the pump.
- 2** With the pump keypad facing up, hold onto the tab and turn the battery door to the left (counterclockwise) until it is released from the battery compartment.
- 3** Pull the door free and remove any existing batteries.
- 4** Install two fresh, disposable AA batteries into the battery compartment. Ensure the batteries are inserted as shown.
- 5** Return the battery door to the battery compartment.
- 6** Turn the battery door to the right (clockwise). Ensure the battery door is secured under the slot.
- 7** Flip closed the tab of the battery door to lock in place.

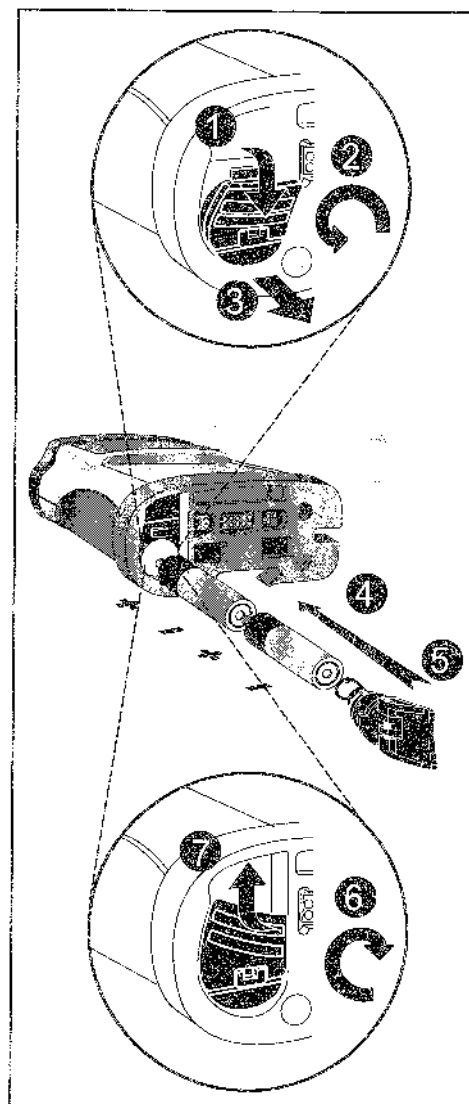

2. Basic Operation


Connecting the GemStar AC Mains Adaptor

- 1 Connect the cord of the AC mains adaptor to the connector labeled "3 VDC" on the bottom of the pump.
- 2 Plug the AC mains adaptor into a standard wall outlet.

CAUTION: ALWAYS connect to a grounded AC mains power source when using the AC mains adaptor. USE ONLY AC mains adaptors specifically labeled for use with the GemStar Pump.

When removing the cord from the bottom of the pump, firmly grasp the connector and pull straight out. DO NOT twist or bend the cord or connector.

Note: The green LED  on the pump remains lit when connected to AC mains power.

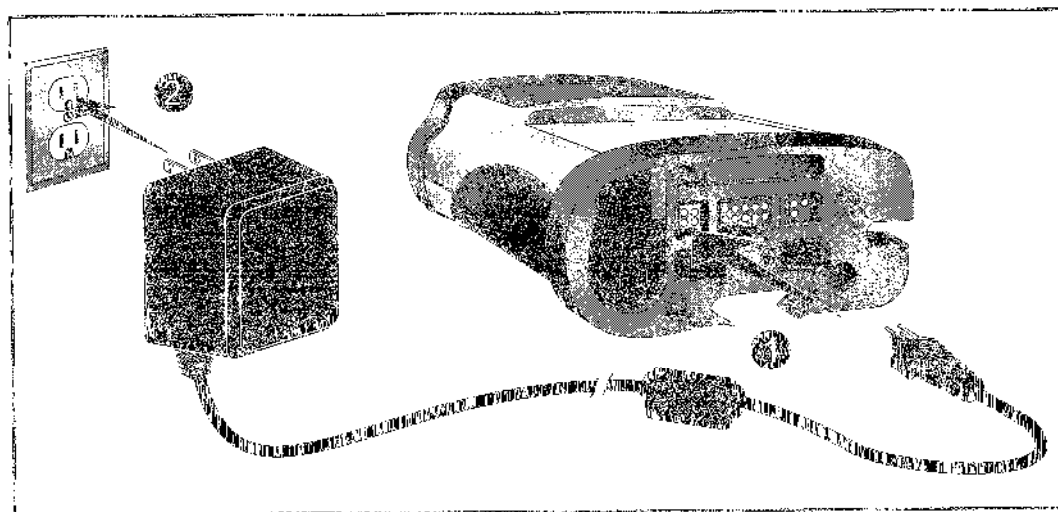

2. Basic Operation

Using the GemStar Rechargeable Battery Pack

Use the GemStar AC Mains Adaptor (List 13026-01 [USA] or 13072-01 [Universal]) to recharge the GemStar Battery Pack.

CAUTION: ALWAYS connect to a grounded AC mains outlet when using AC mains adaptor.


Use only AC mains adaptors specifically labeled for use with the GemStar Pump to charge the battery pack.

During charging, the battery pack is warm. If the battery pack becomes hot to the touch, IMMEDIATELY unplug AC mains power and contact Hospira Technical Support Operations.

Note: A charged battery pack will last approximately 24 hours when the pump is infusing at rates less than 125 mL/hr.

To connect the battery pack:

- 1 Line up the connectors on the top of the battery pack with the holes on the bottom of the pump.
- 2 Snap the battery pack into place.

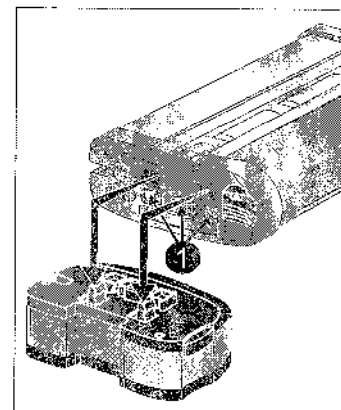


To release the battery pack:

- 1 Press the release button on the battery pack.
- 2 Pull the battery pack straight out from the pump.

Note: When the pump is connected to the battery pack, the pump's green LED flashes and the battery pack light is off.

When the pump is connected to the battery pack at power-on, the pump displays USING EXT BATT.

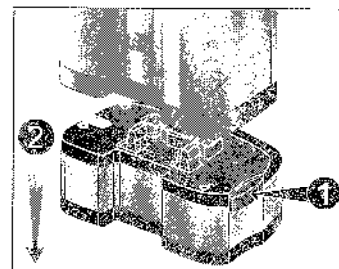


2. Basic Operation

To recharge the battery pack:

- 1 Plug the pin connector of the AC mains adaptor into the port on the bottom of the battery pack.
- 2 Plug the AC mains adaptor into a standard wall outlet.
- 3 Leave the battery pack connected to the AC mains adaptor for at least 3 hours to recharge the battery pack.

CAUTION: When removing the cord from the bottom of the pump, firmly grasp the connector and pull straight out. **DO NOT** twist or bend the cord or connector.


Note: The battery pack light is yellow while the battery pack is charging. When the battery pack is charged, the battery pack light is green.

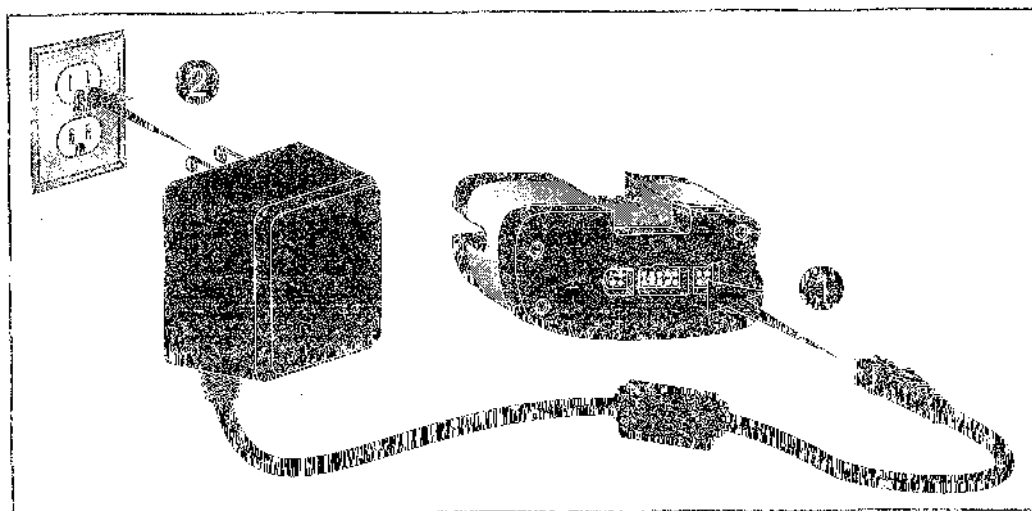

Connecting to the GemStar Docking Station

Connecting the Pump

- 1 Mount the docking station to the I.V. pole and tighten the knob.
- 2 Line up the connectors on the bottom of the pump with the connectors on the docking station.
- 3 Slide the pump into the docking station until it snaps into place.
- 4 Connect the docking station to AC mains power.

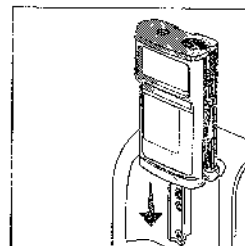


Note: The docking station light is yellow while the docking station is charging. When the docking station is charged, the light is green.

Disconnecting the Pump

To remove the pump:

- 1 Press the release button on the top of the docking station.
- 2 Slide the pump out.

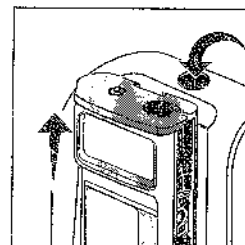

2. Basic Operation

Powering On the Pump

CAUTION: If the pump does not display the self-test or if beeps do not sound at power-on, check all power connections and install two fresh, disposable AA batteries. If the problem continues, contact Hospira Technical Support Operations before using the pump.

Note: If a cassette is not properly loaded, the pump will alarm after the power-on test is complete.

- 1 Press [ON/OFF] to power-on the pump. The pump completes a self-test and displays the time and date. No response is required.

GEMSTAR 7
UNIT SELF-TEST
TIME IS 10:00 AM
FRI, JUN 11, 04

*Note: The self-test takes approximately 30 seconds to complete. If the pump is locked, lock symbols appear at the beginning and end of the second line:
Ⓜ UNIT SELF-TEST Ⓜ.*

- 2 The pump checks for available power sources. Respond to the display message as shown in the following table:

Display Message	Power Source(s)	User Response
No message displays	AC power and disposable batteries.	No response required, the pump goes to the next step.
LOW BATTERIES ENTER TO CONFIRM	AC power and weak (low voltage) disposable batteries; or weak disposable batteries only.	Press  to confirm. Power-off the pump and install two fresh, disposable AA batteries as a backup power source.
NO INTERNAL BATT NO POWER LOSS ALARM AVAILABLE ENTER TO CONFIRM	No disposable batteries or dead disposable batteries; and external power.	Press  to confirm. Power-off the pump and install two fresh, disposable AA batteries as a backup power source.

2. Basic Operation

Display Message	Power Source(s)	User Response
USING BATTERIES ENTER TO CONFIRM	Disposable batteries only.	Press  to confirm. If an external power source is in use, ensure the connections are secure.
USING EXT BATT ENTER TO CONFIRM	External batteries (docking station battery or battery pack) and disposable batteries.	Press  to confirm. If an AC mains adaptor is in use, ensure the connections are secure.

Note: To power-off the pump, press and hold  until the display clears. The pump cannot be powered off during a CHECK CASSETTE alarm (except during PROGRAMMING mode). Refer to Chapter 12 Troubleshooting for information on clearing a CHECK CASSETTE alarm.



For more information on:

Setting the clock
Locking the keypad

Refer to the following pages:

Set Clock in Chapter 10 Options Menu, on page 108.
Keypad Lock in Chapter 10 Options Menu, on page 104.

2. Basic Operation

Accessing the Programming Menu

After power-on, the pump checks its memory for a stored program. Respond to the display message as shown in the following table:

Display Message	Condition	User Response
The pump displays the Programming Menu.	Pump memory is clear.	> To enter a new program: Select a therapy from the Programming Menu and enter a new program. Refer to the appropriate therapy chapter for instructions.
1 RESUME PROGRAM 2 CLR PROG&SHIFT 3 CLR PROG, SHIFT AND HISTORY	A program is currently stored in pump memory.	> To continue the current program, press [1]. > To reset the current program: 1 Press [1]. 2 Press [CHANGE] to display the Change Menu and select NEW CONTAINER. > To enter a new program and clear shift totals: Press [2]. The pump clears the program and shift totals from memory and enters PROGRAMMING mode. > To enter a new program, clear shift totals, and clear history: Press [3]. The pump clears the program, shift totals, and history from memory and enters PROGRAMMING mode.



Tip: To enter a new program when the pump is already powered on:

- 1 Press [STOP] to place the pump in STOP mode.
- 2 Press [CHANGE] to display the Change Menu.
- 3 Select NEW PROGRAM. The pump clears the current program and displays the Programming Menu.

Programming the Pump

Programming Overview

The pump's menu system provides step-by-step guidance through all of the programming steps. Refer to the appropriate therapy chapter for detailed programming instructions.

Programming worksheets are available at the end of each therapy chapter. Each worksheet lists the minimum and maximum value (program range) allowed for that therapy. If a number outside of the program range is entered, the pump displays the minimum or maximum value allowed.

- 1** Select the desired therapy from the Programming Menu.
- 2** Respond to display messages and enter requested values.

Note: Program entries can be changed any time before the Program Review is complete. Press [BACK-UP] to return to previous program displays and make changes as desired.

- 3** Review the program, if required. Press  to begin review. Press [YES/ENTER] when done.
- 4** SAVING PROGRAM displays. The pump stores the program in memory and then enters STOP mode.

 **Tip:** If you enter an incorrect number during programming, press [CHANGE] BEFORE pressing [YES/ENTER] to reset the number to zero. Then enter the correct number.

2. Basic Operation

Programming Tips

The following table explains how to use the keypad during programming:

Key	Function
[0] - [9]	Press the number keys to select menu items and to enter numeric values. <i>Note: Some menu items require pressing two number keys to select that item.</i> <i>For example, when selecting a unit of measure for a weight dosed program, the menu may display as:</i> 01 mcg/kg/mg 02 mcg/kg/hr 03 mcg/kg/min 04 mg/kg/min <i>To select mg/kg/min, press [0][4].</i>
	Press  to accept entries or to continue when done viewing a display.
	Press  during programming and alarm conditions to display information about the current condition.
 	Press  or  to scroll through display messages. Press  to insert a decimal point in a numeric value.
	When entering a numeric value during programming, press  before pressing  to reset the value to zero.
	Press  to return to a previous display, to exit the Options Menu, or to exit the HELP mode.

Starting a New Container

Note: The NEW CONTAINER function is not available when the pump is in FULL Lock. Refer to Chapter 10 Options Menu, page 97 for more information on locking the keypad.

To repeat the current program with a new container:

- 1 From the STOP mode, press [CHANGE] to display the Change Menu.
- 2 Select NEW CONTAINER from the Change Menu. The pump displays NEW CONTAINER for several seconds and clears the amount infused from memory.
- 3 Review the program. Press  to begin review.
- 4 Press [START] to begin delivery.

GemStar Pump Set

The GemStar Pump Set is a sterile, single-use, disposable pump set. It is designed for one-time use only and should be properly discarded after use.

Fluid path and areas beneath undisturbed protective set covers are sterile and nonpyrogenic in the intact unit package.

STERILE 


Contact a Hospira representative for more information on selecting the appropriate pump set configuration. Additional disposable components, such as air-eliminating filters and extension sets, may be added to the line as required. To use a GemStar Pump Set, follow the instructions included with the set.

Note: GemStar Pump Sets contain an anti-siphon valve to reduce the potential for gravity flow at a typical 36-inch head height.

False occlusion alarms and fluid delivery inconsistencies may occur when pumping viscous fluids in microbore tubing at rates greater than 500 mL/hr.

2. Basic Operation

A typical GemStar Pump Set is shown below:

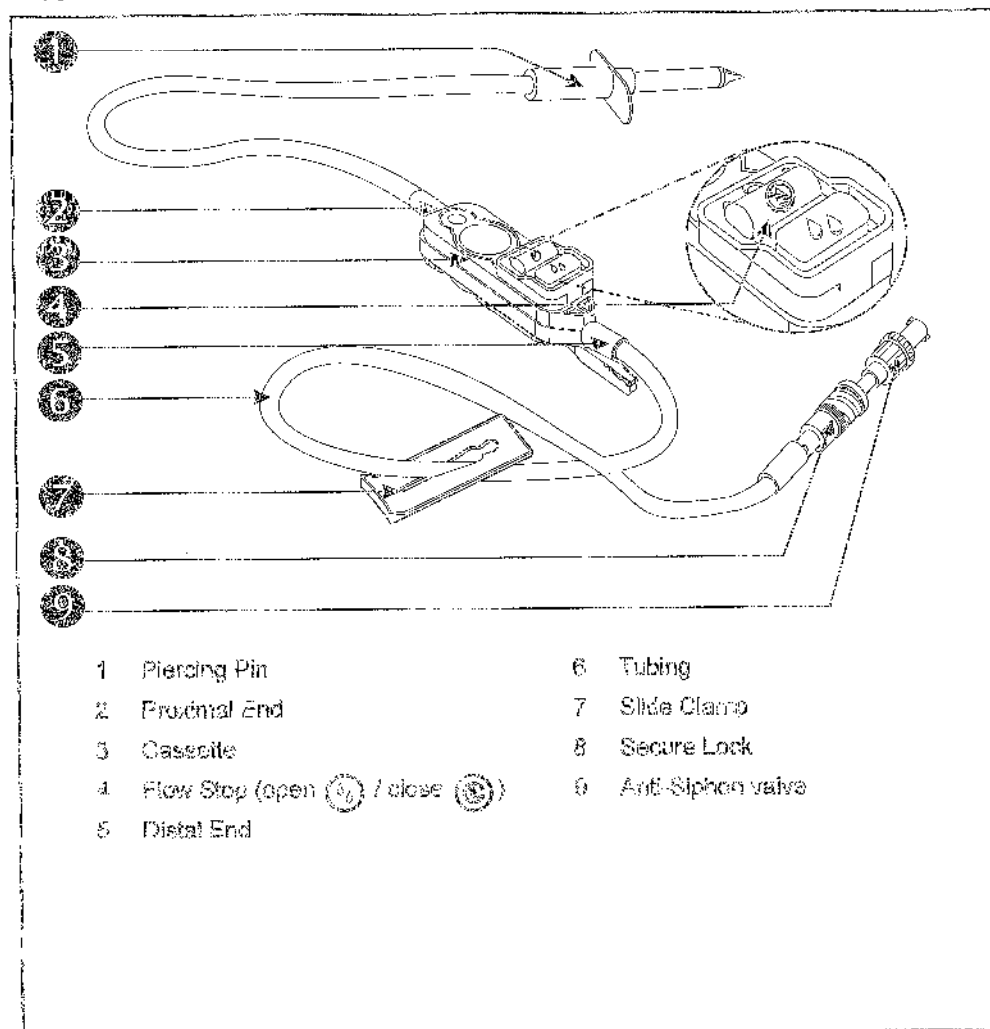

2. Basic Operation

WARNING: ALWAYS prime the pump set to remove air from the cassette, tubing, and injection sites prior to connecting to the patient. ALWAYS disconnect the pump set from the patient prior to priming or purging. Arrange tubing, cords, and cables to minimize the risk of patient strangulation or entanglement.


Failure to use the anti-siphon valve may result in unrestricted flow.

CAUTION: USE ONLY GemStar Pump Sets with the GemStar Pump. Use of unauthorized sets may result in injury to the patient or damage to the pump.

To prevent contamination, use aseptic technique with all fluid-path connections. Remove protective coverings as assembly progresses.

Opening the Pump Set


- 1 Open the GemStar Pump Set package and remove the contents.
- 2 Remove the protective covers from the fluid container administration port and the pump set.

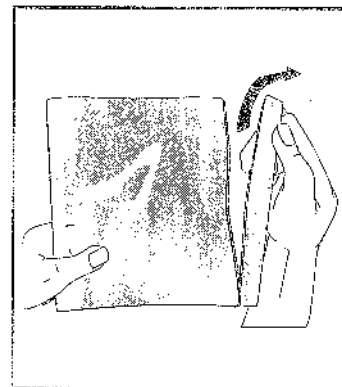


2. Basic Operation

Loading the Cassette


- 1 Identify the four black cassette latches inside the cassette pocket.

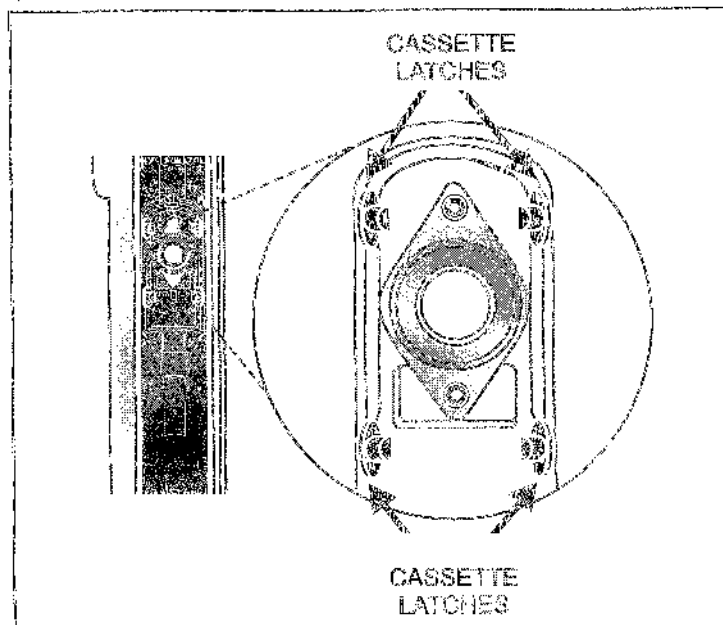

- 2 Line up the cassette as shown (with the flow stop button facing the pump).

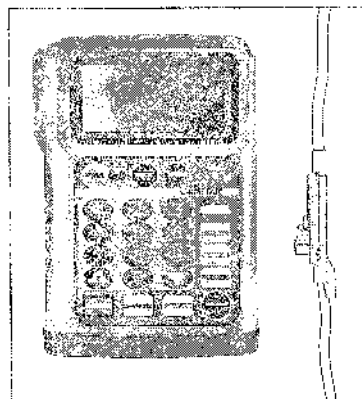


- 3 Insert the cassette into the cassette pocket by pushing firmly along the entire cassette until it is firmly seated.

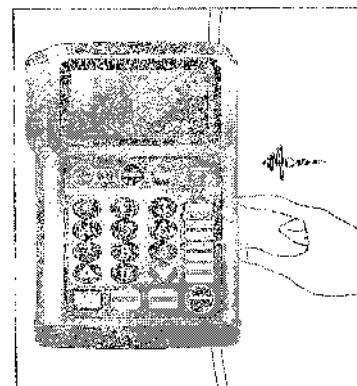


2. Basic Operation

- 4 Verify that all four cassette latches are securely holding the cassette.

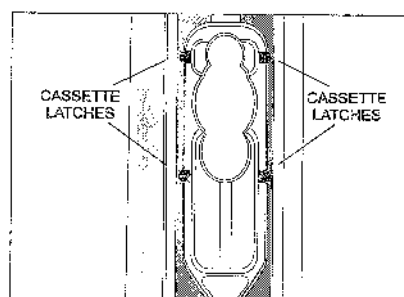

CAUTION: Ensure the cassette is properly installed. When properly installed, all four cassette latches are visibly holding the cassette securely in the cassette pocket.

- 5 Place the tubing in the tubing channel.

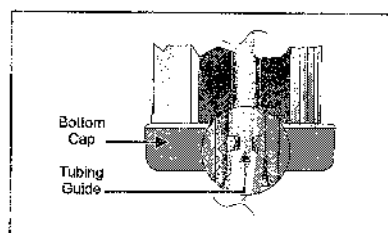


Note: Ensure the tubing is behind the tubing guide (at the distal end).


Note: The pump cannot be powered off during a CHECK CASSETTE alarm (except during PROGRAMMING mode). Refer to Chapter 12 Troubleshooting for information on clearing a CHECK CASSETTE alarm.

Releasing the Cassette

- 1 Close the CAIR® (roller) clamp or slide clamp on the distal line.

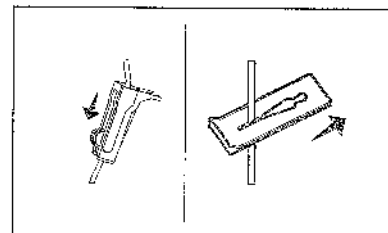

- 2 Push down on the cassette release (large black button) on the top of the pump.

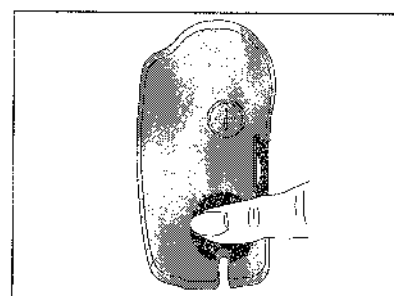

- 3 Remove the cassette.

Note: When the cassette is released from the pump, the cassette flow stop automatically closes, preventing free flow.


2. Basic Operation

Priming a GemStar Pump Set


Priming Pump Sets WITHOUT a Drip Chamber

WARNING: Failure to use the anti-siphon valve may result in unrestricted flow.


- 1 Push down on the flow stop (A) to close the cassette.

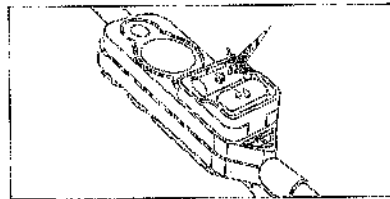


- 2 Hold the container with the administration port toward the ceiling and insert the piercing pin.

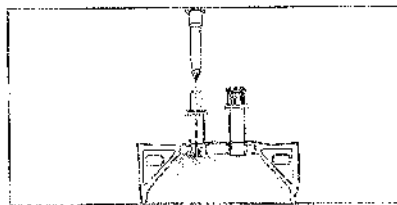


- 3 Push down on the flow stop (A) to open (priming position) the cassette.

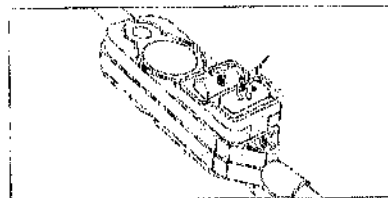


- 4 Invert the cassette approximately 45°.

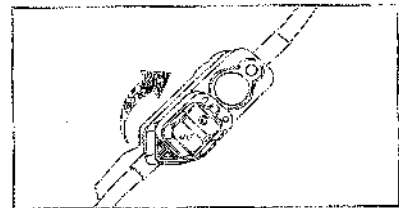

- 5 While holding the rotated cassette, squeeze or roll the container to remove air from the container and cassette.

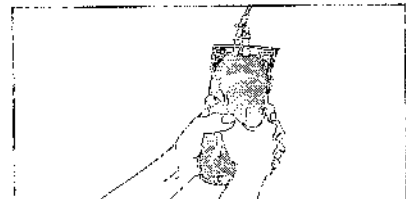


- 6 Return the cassette to the upright position. Continue priming until fluid fills the tubing.

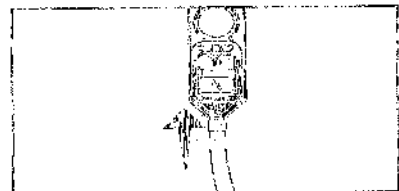


2. Basic Operation


- 7 Push down on the flow stop  to close the cassette.

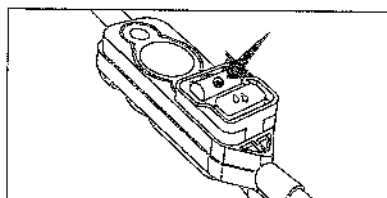

Note: Ensure the tubing is behind the tubing guide (at the distal end).

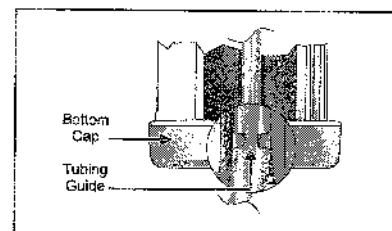


- 8 Load the cassette.

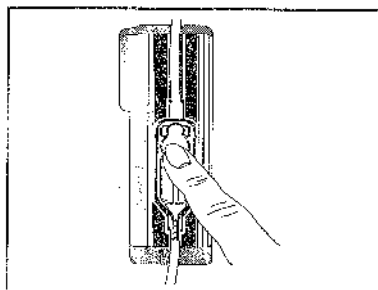

CAUTION: When priming is complete, ensure no fluid flows at the distal end of the pump set. If flow is observed, **DO NOT** use the pump set.

- 9 Connect the pump set to the patient-access device.

CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

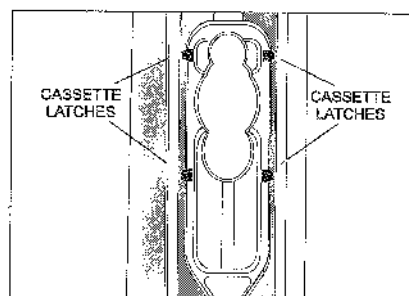


2 Basic Operation

Priming Pump Sets **WITHOUT** a Drip Chamber Using the **PURGE** Key


WARNING: Failure to use the anti-siphon valve may result in unrestricted flow.

Note: The PURGE function is available only if there is a program in the pump.


- 1 Push down on the flow stop (Ⓢ) to close the cassette.

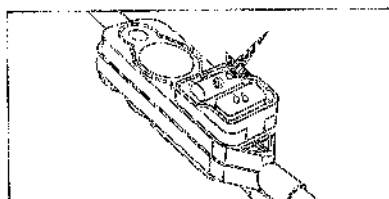


- 2 Hold the container with the administration port toward the ceiling and insert the piercing pin.

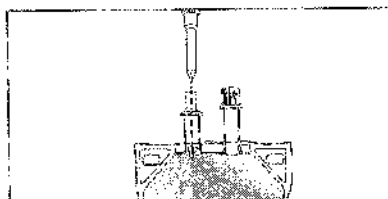


- 3 Load the cassette.

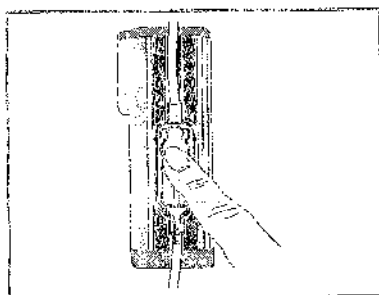

CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

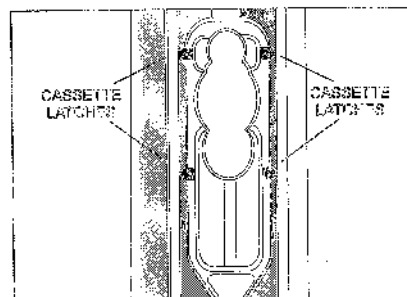


Note: Ensure the tubing is behind the tubing guide (at the distal end).

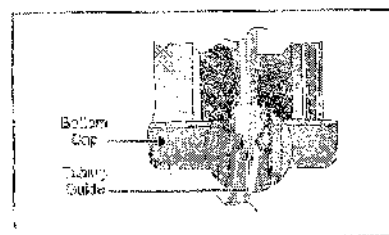

- 4 Program desired therapy.

2. Basic Operation

- 5 Press [PURGE].

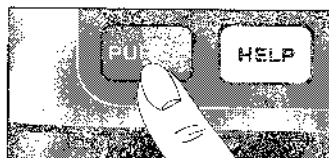

- 6 The pump displays PRIME THE SET? Press [YES/ENTER].
- 7 While holding the container with the administration port toward the ceiling, press and hold [PURGE].

*Note: Purge rate is 250 mL/hr.
The pump purges for up to 2 minutes each time you press and hold [PURGE].*

Continue using the [PURGE] key until fluid fills the tubing. Ensure all air is removed from the container, cassette, and tubing.

CAUTION: When priming is complete, ensure no fluid flows at the distal end of the pump set. If flow is observed, **DO NOT** use the pump set.

- 8 Connect the pump set to the patient-access device.


2. Basic Operation

Priming Pump Sets WITH a Drip Chamber


WARNING: Failure to use the anti-siphon valve may result in unrestricted flow.

When using a vented chamber, **ALWAYS** use an air-eliminating filter OR set the air sensitivity to either ON or 2 mL.


- 1 Push down on the flow stop (⑤) to close the cassette.

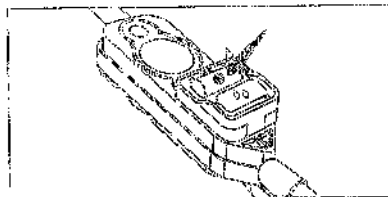


- 4 Push down on the flow stop (⑤) to open (priming position) the cassette.

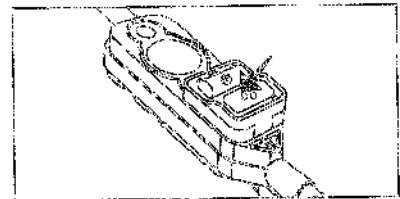


- 2 Insert the piercing pin into the container. Suspend the container on an I.V. Pole.

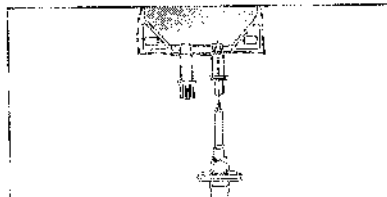


- 5 Ensure the CAIR® (roller) clamp is open.

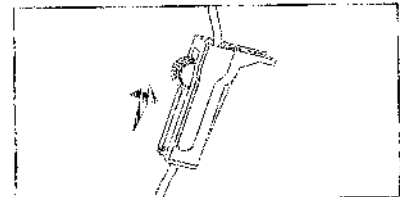

- 3 Squeeze the drip chamber until filled halfway.

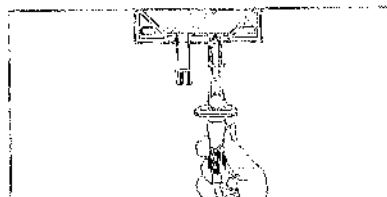


- 6 Invert the cassette approximately 45°.

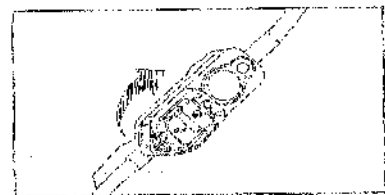


2. Basic Operation


- 7** While holding the rotated cassette, remove air from the container and cassette.

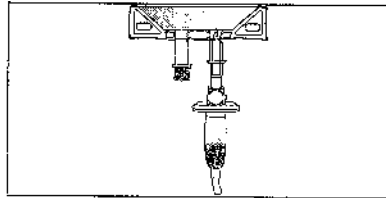


- 8** Return the cassette to the upright position. Continue priming until fluid fills the tubing.

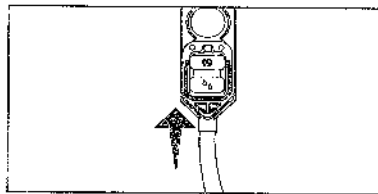


- 9** Push down on the flow stop (C) to close the cassette.

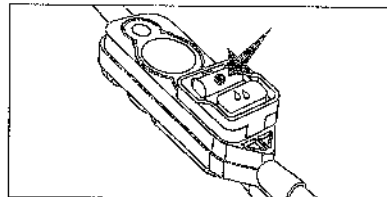


- 10** Load the cassette.

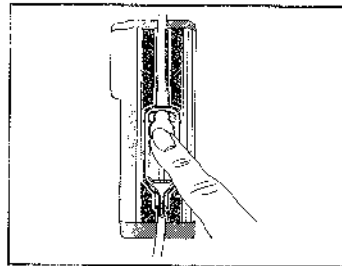

CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

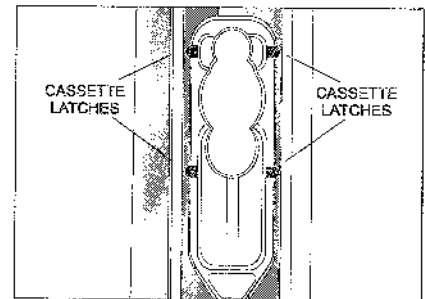

Note: Ensure the tubing is behind the tubing guide (at the distal end).

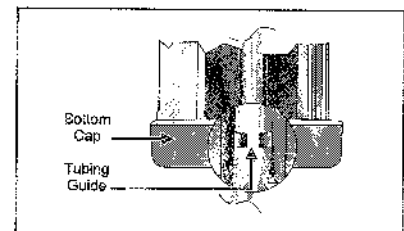

CAUTION: When priming is complete, ensure no fluid flows at the distal end of the pump set. If flow is observed, **DO NOT** use the pump set.

- 11** Connect the pump set to the patient-access device.

2. Basic Operation


Priming Pump Sets WITH a Drip Chamber Using the PURGE Key

WARNING: Failure to use the anti-siphon valve may result in unrestricted flow.


When using a vented set, ALWAYS use an air-eliminating filter OR set the air sensitivity to either ON or 2 mL.

Note: The Purge function is available only if there is a program in the pump.


- 1 Push down on the flow stop (A) to close the cassette.

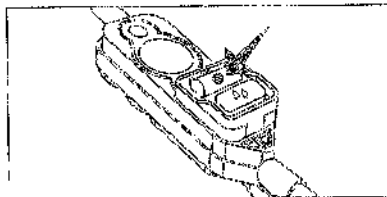


- 2 Insert the piercing pin into the container. Suspend the container on an I.V. Pole.

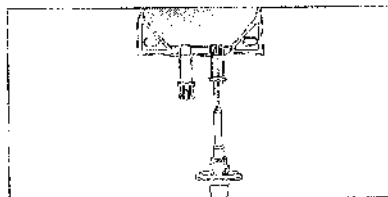


- 3 Squeeze the drip chamber until filled halfway.

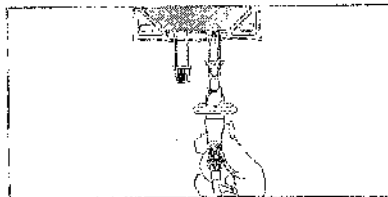

- 4 Push down on the flow stop (B) to open (priming position) the cassette.

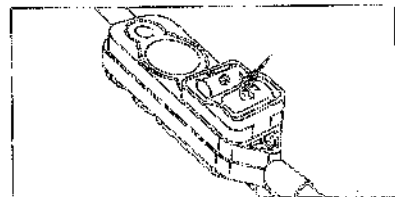


- 5 Ensure the CAIR® (roller) clamp is open.

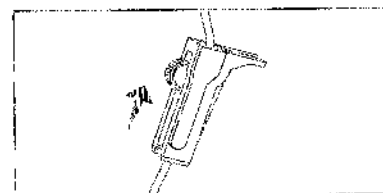


2. Basic Operation


- 6 Load the cassette.

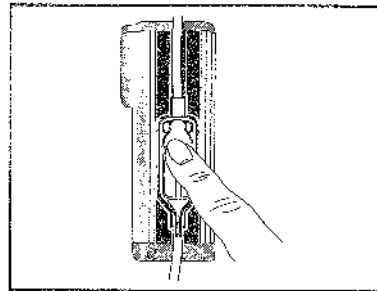

CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

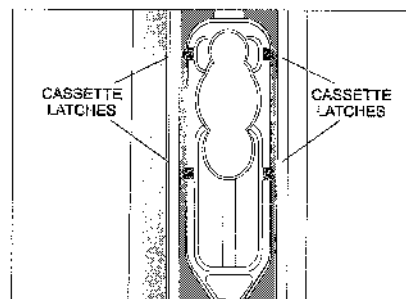


Note: Ensure the tubing is behind the tubing guide (at the distal end).

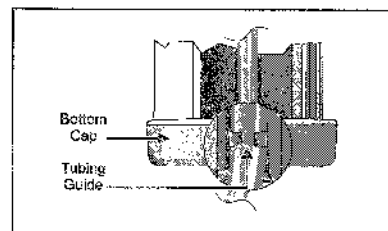

- 7 Program the desired therapy.

- 8 Press [PURGE].

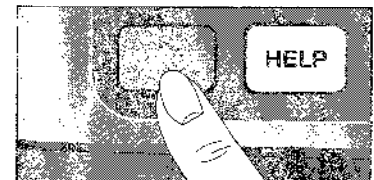

- 9 The pump displays PRIME THE SET? Press [YES/ENTER].

- 10 Press and hold [PURGE].

Note: Purge rate is 250 mL/hr. The pump purges for up to 2 minutes each time you press and hold [PURGE].

Continue using the [PURGE] key until fluid fills the tubing. Ensure all air is removed from the container, cassette, and tubing.

CAUTION: When priming is complete, ensure no fluid flows at the distal end of the pump set. If flow is observed, DO NOT use the pump set.

- 11 Connect the pump set to the patient-access device.


2. Basic Operation

Priming Pump Sets Using Pre-filled Vials


WARNING: Failure to use the anti-siphon valve may result in unrestricted flow.

CAUTION: When using pre-filled vials, the system may underdeliver at low flow rates.


- 1 If using a pump set with the Y-extension, use the slide clamp to close off the KVO (Keep Vein Open) side of the tubing.

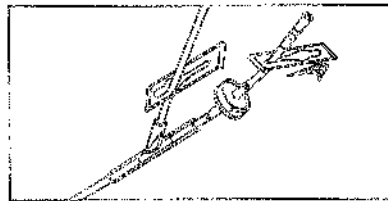

- 2 Open the upper and lower slide clamps.

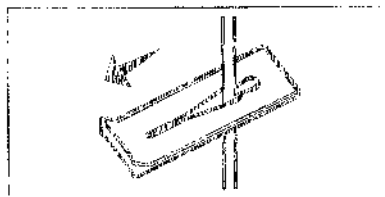


- 3 Push down on the flow stop (A) to close the cassette.

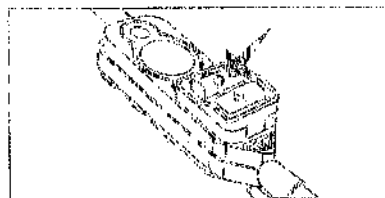

- 4 Connect the pump set to a Hospira pre-filled vial and injector.




- 5 Push down on the flow stop (A) to open (priming position) the cassette.

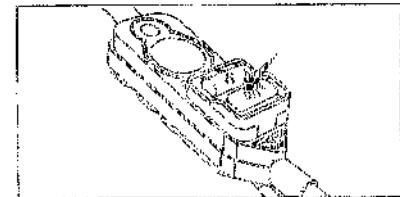


2. Basic Operation


- 6 Hold the vial and set upright (vertically) on a flat surface with the pump set extending from the top.
- 7 Slowly push down on the injector until fluid flows at the distal end of the tubing, which indicates that all air is cleared from the vial, injector, cassette, and tubing.

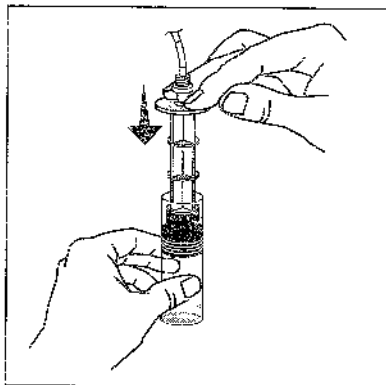


- 8 Push down on the flow stop (Ⓢ) to close the cassette.

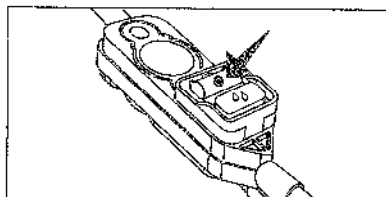


- 9 If the set includes a KVO Y-extension, close the slide clamp below the pre-filled vial and open the slide clamp on the KVO side. Prime the KVO Y-extension.

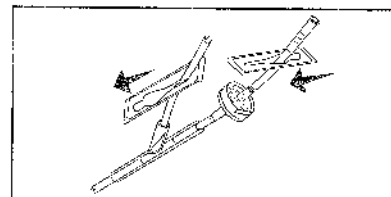

- 10 Load the cassette.

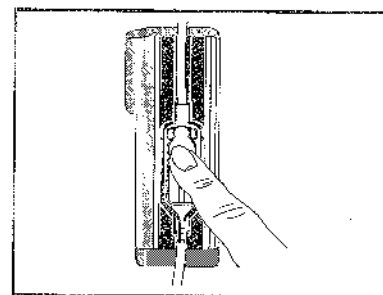


CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

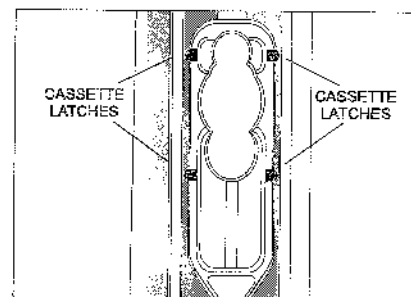

2. Basic Operation



Note: Ensure the tubing is behind the tubing guide (at the distal end).

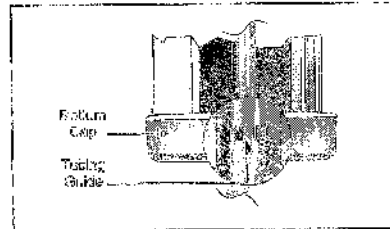

CAUTION: When priming is complete, ensure no fluid flows at the distal end of the pump set. If flow is observed, **DO NOT** use the pump set.


- 11 Connect the pump set to the patient-access device.

2. Basic Operation


Priming Pump Sets with Extensions

WARNING: Failure to use the anti-siphon valve may result in unrestricted flow.


- 1 Push down on the flow stop  to close the cassette.

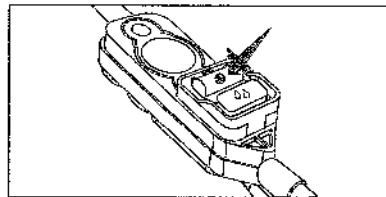


- 2 Hold the container with the administrative port toward the ceiling and insert the piercing pin.

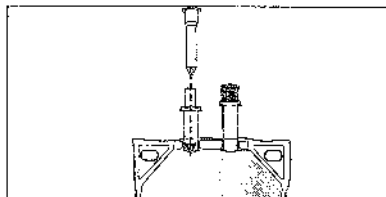


- 3 Push down on the flow stop  to open (priming position) the cassette.

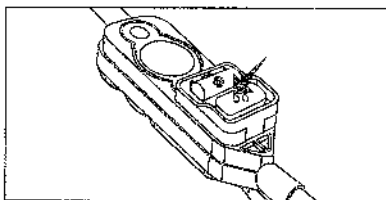


- 4 Invert the cassette approximately 45°.

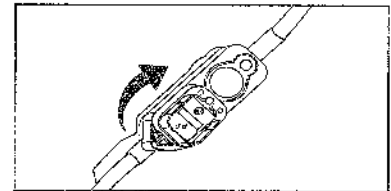

- 5 While holding the rotated cassette, squeeze or roll the container to remove air from the container and cassette.

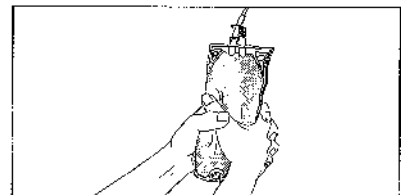


- 6 Return the cassette to the upright position. Continue priming until fluid fills the tubing.

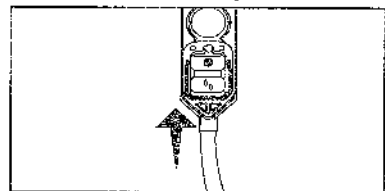


2. Basic Operation


- 7 Push down on the flow stop  to close the cassette.

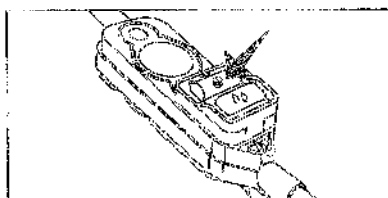


- 8 Connect the extension set to the pump set.
- 9 Load the cassette.

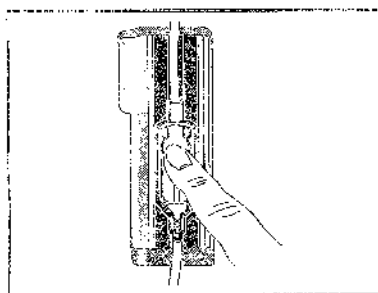


CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

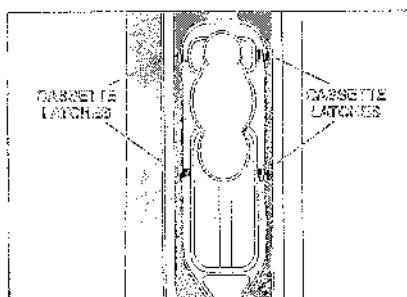

Note: Ensure the tubing is behind the tubing guide (at the distal end).

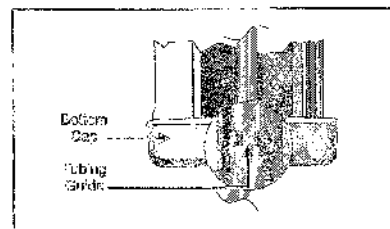

- 10 Program the desired therapy.
- 11 Press [PURGE].

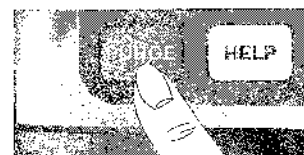

- 12 The pump displays PRIME THE SET? Press [YES/ENTER].
- 13 Press and hold [PURGE].

Note: Purge rate is 250 mL/min. The pump purges for up to 2 minutes each time you press and hold [PURGE].

Continue using the [PURGE] key until fluid fills the tubing. Ensure all air is removed from the container, cassette, and tubing.

CAUTION: When priming is complete, ensure no fluid flows at the distal end of the pump set. If flow is observed, **DO NOT** use the pump set.

- 14 Connect the pump set to the patient-access device.

Total Parenteral Nutrition (TPN)

Therapy Features and Specifications

This therapy allows programming of Total Parenteral Nutrition (TPN) protocols.

- ♦ Taper Up and Taper Down may be programmed to gradually increase and decrease the infusion rate at the beginning and end of a therapy. If Taper Down is selected during programming, the pump allows auto tapering if the continuous delivery is stopped before the programmed volume has been infused.
- ♦ Set an optional KVO (Keep Vein Open) during programming.



Tip: The Programming Worksheet at the end of this chapter lists the minimum and maximum values allowed for this protocol. Use copies of the programming worksheet to record program values before programming the pump.

TPN Programming Tips

The following table explains how to use the keypad during TPN programming:

Key	Function	Key	Function
[0] - [9]	Selects menu items and enters numeric values.	 	Scrolls through program review and help messages. The up arrow inserts a decimal point in a numeric value.
	Accepts entries or continues when done viewing a display.		Clears a numeric value (resets to zero) during programming.
	Displays information during programming and alarm conditions.		Returns to a previous display, exits the Options Menu, and exits the HELP mode.

3. Total Parenteral Nutrition (TPN)

TPN Programming

1 Select TPN from the Programming Menu

2 Select the TPN delivery mode.

3 Enter the total VTBI (Volume To Be Infused).

Note: Maximum VTBI is 9600 mL minus any taper amount.

4 Enter the Taper Up and Taper Down times (if applicable).

Note: Time is entered in hours:minutes (00:00). For one hour thirty minutes, press [1] [3] [0]; not [1] [0].

5 Enter the total time to infuse.

6 Enter the KVO rate. If no KVO is desired, press [YES/ENTER] while the value displayed is zero.

Note: KVO runs after Taper Down until the total container is infused.

7 Enter the container size (VTBI + KVO).

Note: If KVO is programmed, increase the container size to accommodate the desired KVO volume.

8 Select the air sensitivity, if required.

9 Review the program. Press  to begin review.

10 If using the PURGE key to prime the pump set, press [PURGE] now.

WARNING: ALWAYS disconnect the pump set from the patient before priming or purging.

3. Total Parenteral Nutrition (TPN)

Changing a TPN Program

Can be changed before  is pressed	Can be changed at any time	Cannot be changed
• VTBI • Taper times • Total TPN time	• KVO rate • Increase container size only • Air sensitivity	• Delivery mode

 *Note: You cannot change a program when the pump is in FULL or CONTAINER Lock.*

To change a program:

- 1 From the STOP mode, press [CHANGE] and select CHANGE PROGRAM.
- 2 As each message displays, make the desired changes or press [YES/ENTER] to accept the current setting.
- 3 Review the program. Press  to begin review.
- 4 Press [START] to begin infusion.

3. Total Parenteral Nutrition (TPN)

TPN Auto Taper Down

Auto Taper Down is available only if a Taper Down time was entered during programming. Auto Taper Down may be activated when the continuous delivery portion of a TPN with Taper Down protocol is stopped. Auto Taper Down is not available while the pump is in an alarm condition.

- 1 Press [STOP] during TPN delivery.

The pump displays DO YOU WANT AUTO TAPERING?

- Press [NO] to return to the STOP mode. Then press [START] to continue therapy.
- Press [YES/ENTER] to set Auto Taper Down.

- 2 Enter the Auto Taper Down time.

The maximum Auto Taper Down time is three hours OR two times the remaining volume divided by the continuous infusion rate, whichever is less.

Note: Time is entered in hours:minutes (00:00). For one hour thirty minutes, press [1] [3] [0]; not [2] [0].

- 3 Press [YES/ENTER] to return to the RUN mode, then press [START] to begin Auto Taper Down.

3. Total Parenteral Nutrition (TPN)

GemStar® TPN Programming Worksheet

GIVE: _____

Parameters	Selections and/or Values to be Entered	Programmable Ranges
Therapy	TPN	N/A
Delivery Mode	☐ Continuous ☐ Continuous + Taper Down ☐ Continuous + Taper Up ☐ Continuous + Taper Up and Down	N/A
VTBI (Volume To Be Infused)	_____ mL	1 - 9600 mL*
Taper Up	☐ No ☐ Yes _____ hrs:min	1 min - 3 hrs
Taper Down	☐ No ☐ Yes _____ hrs:min	1 min - 3 hrs**
Total Time to Infuse	_____ hrs:min	1 min - 24 hrs
KVO (Keep Vein Open)	☐ No KVO Rate: 0.0 ☐ Yes KVO Rate: _____	1.0 - 5.0 mL/hr
Container Size (VTBI + KVO)	_____ mL	9999 mL max
Air Sensitivity	☐ On Alarms at approx. 0.5 mL of air ☐ 2 mL Alarms at approx. 2 mL of air ☐ Off No alarm	N/A

* Maximum VTBI is 9600 mL minus any taper amount.

** Programmable ranges for taper down are 1 minute to 3 hours OR 2 times the remaining volume divided by the continuous delivery rate, whichever is less.

Notes

Pain Management

Therapy Features and Specifications

This therapy allows programming of a continuous only, bolus only, or continuous with bolus delivery.

- ♦ In Continuous + Bolus protocols, you can enter zero for either the continuous rate or the bolus dose when entering or changing a program. The continuous rate and bolus dose cannot be zero at the same time.
- ♦ The default bolus delivery rate is 125 mL/hr.
- ♦ A loading dose may be set during programming, which can be delivered after programming and priming are complete. The default loading dose rate is 125 mL/hr.
- ♦ Priming using [PURGE] is only allowed before [START] is pressed or after one of the following events occur:
 - ♦ Complete Review of New Program
 - ♦ New Container
 - ♦ Air-in-Line Alarm

CAUTION: When using pre-filled vials, the system may underdeliver at low flow rates.

 *Note: New Container is allowed only when at least 1/2 of the container volume has been infused or 1 hour has elapsed since starting the infusion.*

The use of a lockbox is recommended with pain management protocols. Refer to Chapter 15 Optional System Components, page 153, for instructions on using a lockbox with the GemStar Pump.



Tip: The Programming Worksheet at the end of this chapter lists the minimum and maximum values allowed for this protocol. Use copies of the programming worksheet to record program values before programming the pump.

4. Pain Management

Pain Management Programming Tips

The following table explains how to use the keypad during pain management programming:

Key	Function	Key	Function
[0] - [9]	Selects menu items and enters numeric values.	 	Scrolls through program review and help messages. The up arrow inserts a decimal point in a numeric value. The down arrow toggles units of measure when entering concentrations.
	Accepts entries or continues when done viewing a display.		Clears a numeric value (resets to zero) during programming.
	Displays information during programming and alarm conditions.		Returns to a previous display, exits the Options Menu, and exits the HELP mode.

Pain Management Programming

- 1 Select **PAIN MANAGEMENT** from the Programming Menu.
- 2 Select the delivery mode.

*Note: If Continuous + Bolus is selected, the continuous delivery rate or bolus dose can be set to zero through the **CHANGE** menu without reprogramming the pump.*

- 3 Select the unit of measure and concentration as required.

4. Pain Management

- 4 Set the continuous delivery rate.

Note: The pump only delivers in increments of 0.1 mL. Values entered in a Pain Management therapy may be rounded to the nearest 0.1 mL. For example, at a concentration of 15 mg/mL, a rate of 5 mg rounds to 4.5 mg.

$$5 \text{ mg} = 0.333 \text{ mL}$$

$$4.5 \text{ mg} = 0.3 \text{ mL}$$

When a value is rounded, the pump sounds four quick beeps, displays ROUNDING for a few seconds, and then displays the rounded value. Press [YES/ENTER] to accept the rounded value and continue to the next step.

- 5 Program a loading dose as required.
- 6 Program a bolus dose with lockout time and dose limit as required.

Regarding a 1-Hour or 4-Hour Dose Limit: The hour dose limit is the maximum volume (continuous plus bolus or bolus only amount) that can be delivered over the selected period. When the programmed dose limit is reached, the continuous delivery stops and a new bolus is not allowed. As the oldest delivery amounts age out of the 1-hour or 4-hour dose limit record, the pump accepts bolus requests or resumes continuous delivery (if programmed).

Note: If a 1-hour or 4-hour dose limit is programmed, a bolus in progress will stop when the amount infused reaches the dose limit.

Loading doses are included in a 1-hour or 4-hour dose limit. A dose limit in effect could prevent the delivery of a loading dose.

- 7 If programmed in units of mcg or mg, the line displaying the units flashes. Enter the container size in proper unit measure.

Note: To enter the total in mL when programming a concentration, press  to move to the mL line. Enter the mL value when mL flashes on the display.

- 8 Select the air sensitivity, if required.
- 9 Review the program. Press  to begin review.
- 10 If using the PURGE key to prime the pump set, press [PURGE] now.

WARNING: ALWAYS disconnect the pump set from the patient before priming or purging.

4. Pain Management

Bolus Delivery

Set bolus doses and delivery limits during programming. The default bolus dose delivery rate in a pain management protocol is 125 ml/hr.

A bolus cord is available for patient bolus requests. Refer to *Connecting the Bolus Cord* in Chapter 15 Optional System Components, on page 156, for more information.

Delivering a Bolus Dose

From the RUN mode, press either the bolus button [+] on the top of the pump or the button on the end of the bolus cord to begin bolus dose delivery. The pump flashes **BOLUS DELIVERY** on the display and the amount infused accrues as the delivery progresses.

When the bolus dose is complete, the pump does the following:

- + Records the bolus request and delivery amount to the event and bolus history logs
- + Adds the delivered boluses to the total bolus amount in the program amounts
- + Resets the bolus lockout time (if applicable) and ignores other bolus requests until the lockout time has elapsed
- + Continues programmed infusion



Note: ALL bolus requests are recorded to the event and bolus history logs.

If the bolus delivery does not start, it could be locked out by one of the following conditions:

- + Bolus lockout period
- + Loading dose delivery in progress
- + Bolus dose delivery in progress
- + Bolus per-hour limit
- + 1-hour or 4-hour dose limit

4. Pain Management

Bolus Dose Interruptions

Bolus deliveries may be interrupted by:

- ♦ Pressing [STOP]
- ♦ An alarm condition

To continue bolus delivery after an interruption:

- 1** Press [START]. The pump displays COMPLETE BOLUS NOW?
- 2** Respond as follows:

To do this:	Press this key:	Bolus lockout time is set to:
Deliver the remaining bolus amount		The time of delivery completion
Clear the undelivered bolus amount		The time the bolus was interrupted

Loading Dose Delivery

You can set a loading dose during programming, if desired. A loading dose can be delivered after programming and priming are complete.

Note: Loading dose amounts are included in a 1-hour or 4-hour dose limit.

Delivering a Loading Dose

- 1** After priming and programming, press [START].
- 2** The pump displays DELIVER THE LOADING DOSE NOW? Press [YES/ENTER] to begin loading dose.

When the loading dose is complete, the pump does the following:

- ♦ Automatically begins programmed infusion
- ♦ Records the dose amount to the history log
- ♦ Resets the bolus lockout time (if applicable) and ignores other bolus requests until the lockout time has elapsed

4. Pain Management

Loading Dose Interruptions

Loading dose deliveries can be interrupted by:

- Pressing [STOP]
- An alarm condition

To continue a loading dose:

- 1 Resolve the interruption.
- 2 Press [START]. The pump displays **COMPLETE THE LOADING DOSE?**
- 3 Respond as follows:

To do this	Press this key
Deliver the remaining loading dose amount	
Clear the remaining loading dose amount	

Changing a Pain Management Program

Can be changed at any time		Cannot be changed
• Delivery rate	• Container size	• Delivery mode
• Bolus	• Air sensitivity	• Unit of measure
• Bolus lockout	• Dose limit amount	• Concentration
		• Dose limit type

Note: You cannot change a program when the pump is in FULL or CONTAINER Lock. With a RATE CHANGE Lock, the continuous delivery rate and bolus amount can be changed within the set ranges.

Although the delivery mode cannot be changed, if Continuous + Bolus is programmed, the continuous rate or bolus amount can be changed to zero. The minimum values set with a RATE CHANGE Lock could restrict this function.

4. Pain Management

To change a program:

- 1** From the STOP mode, press [CHANGE] and select CHANGE PROGRAM.
- 2** As each message displays, make the desired changes or press [YES/ENTER] to accept the current setting.
- 3** Review the program. Press  to begin review.
- 4** Press [START] to begin delivery.

4 Pain Management

GemStar®

Pain Management Programming Worksheet

GIVE: _____

Parameters	Selections and/or Values to be Entered	Programmable Ranges
Therapy	Pain Management	N/A
Delivery Mode	☐ Continuous ☐ Bolus Only ☐ Continuous + Bolus	P/A
Unit of Measure	☐ Concentration: _____ mg/mL ☐ Concentration: _____ mcg/mL ☐ mL	0.1 - 100 mg/mL 1 - 1000 mcg/mL
Continuous Rate	_____	0.1 - 2500 mg/hr 1 - 25000 mcg/hr 0.1 - 25.0 mL/hr
Loading Dose	☐ No ☐ Yes Amount: _____	0.1 - 2500 mg 1 - 25000 mcg 0.1 - 25.0 mL
Bolus Dose	_____	0 - 2500 mg 0 - 25000 mcg 0 - 25.0 mL
Bolus Lockout	_____ minutes	5 - 999 min
Optional Dose Limit	☐ 4-hour limit Amount: _____ ☐ 1-hour limit Amount: _____ ☐ # boluses/hour Number: _____ ☐ No dose limit	4X rate - 400 mL, 40000 mg, 400000 mcg* Rate - 100 mL, 10000 mg, 100000 mcg* 1 - 12 *For continuous plus bolus delivery.
<i>Note: For Bolus Only, the 1-hour or 4-hour dose limit minimum is 1 bolus dose.</i>		
Container Size	_____	0.1 - 99999 mg 1 - 999999 mcg 0.1 - 9999 mL
Air Sensitivity	☐ On Alarms at approx. 0.5 mL of air ☐ 2 mL Alarms at approx. 2 mL of air ☐ Off No alarm	N/A

Chapter 5

Intermittent

Therapy Features and Specifications

This therapy allows programming of multiple doses at regular intervals.

- Optional KVO (Keep Vein Open) can be set during programming.
- The pump can be disconnected from the patient between dose deliveries. A call back alarm can be set to alert the patient or caregiver to reconnect the pump before the next scheduled delivery begins. The pump must remain powered on for the call back alarm to function. At the set time before each dose, the pump produces an audible and visual alarm and flashes START on the display.



Tip: The Programming Worksheet at the end of this chapter lists the minimum and maximum values allowed for this protocol. Use copies of the programming worksheet to record program values before programming the pump.

Intermittent Programming Tips

The following table explains how to use the keypad during intermittent programming:

Key	Function	Key	Function
[0] - [9]	Selects menu items and enters numeric values.	 	Scrolls through program review and help messages. Toggles between AM and PM in time-based entries. The up arrow inserts a decimal point in numeric values.
	Accepts entries or continues when done viewing a display.		Clears a numeric value (resets to zero) during programming.
	Displays information during programming and alarm conditions.		Returns to a previous display, exits the Options Menu, and exits the HELP mode.

5. Intermittent

Intermittent Programming

- 1 Select **INTERMITTENT** from the Programming Menu.
- 2 Enter the dose amount in mL.
- 3 Enter the infusion time for each dose in hours:minutes (00:00).
For example, press [4] [0] if each dose is to deliver for 40 minutes.
- 4 Enter the frequency of dose deliveries.
For example, press [4] [0] [0] to start a dose every four hours.
- 5 Enter the total number of doses in the container.
To determine the number of doses in the container, divide the volume in the container by the dose amount.
- 6 Enter the KVO rate in mL.
KVO delivers before a delayed start time, between doses, and after the last dose until the total container is infused.
If no KVO is desired, press [YES/ENTER] while the value displayed is zero.
- 7 Enter the container size.
The pump automatically calculates the minimum container size required to deliver the number of programmed doses and KVO between doses.

Note: If KVO is desired before a delayed start, increase the container size to include the desired KVO volume.
- 8 Select the air sensitivity, if required.

Note: DO NOT press [YES/ENTER] after selecting the air sensitivity. If [YES/ENTER] is accidentally pressed and a call back alarm is not desired, press [BACK-UP] to return to the beginning of step 8 and press [NO].
- 9 The pump displays **SET CALL BACK ALARM FOR NEXT DOSE?**
 - Press [NO] if you do not want to set the call back alarm.
 - Press [YES/ENTER] to set the call back alarm.Enter the time before the next dose that the alarm should sound. For example, press [1] [0] to set the alarm to sound 10 minutes before the start of the next scheduled delivery.

5. Intermittent

Note: When programming the call back alarm, be sure to allow enough time for the patient or caregiver to reconnect the pump to the patient and to press [START] BEFORE the scheduled dose start time.

10 The pump displays the current time.

- If the displayed time is correct, press [YES/ENTER].
- To set the clock, press [OPTIONS].

11 The pump displays the current time and START NOW?

- Press [YES/ENTER] to begin delivery when the pump enters the RUN mode.
- Press [NO] to set a delayed start time. The pump displays SET START TIME. Enter the desired start time and press [YES/ENTER] when done.

Note: When entering the time in 12-hour format, press  or  to toggle between AM and PM.

12 Review the program. Press  to begin review.

13 If using the PURGE key to prime the pump set, press [PURGE] now.

WARNING: ALWAYS disconnect the pump set from the patient before priming or purging.

Note: The call back alarm works only when the pump is powered on AND in the STOP mode. DO NOT power-off the pump while it is disconnected from the patient; if you have set the call back alarm.

When the alarm sounds, press [SILENCE] to silence the alarm for three minutes. Then reconnect the pump to the patient.

Tip: To use a delayed start without specifying a start time:

- 1 Press [YES/ENTER].
- 2 Complete the program review.
- 3 Power-off the pump.
- 4 When the pump is powered on again, complete the program review, then press [START] to begin delivery.

5. Intermittent

Intermittent Delivery Interruptions


Missed Start Time

- If the programmed start time is missed, it can be changed through the Options Menu. Refer to *Next Dose* in Chapter 10 Options Menu, on page 111.
- If the pump remains in the STOP mode and the first dose is missed, press [START] to begin delivery immediately and all remaining doses will be shifted.
- If the pump is programmed with a delayed start time and is then powered off before [START] is pressed, the pump will skip to the NEXT occurrence of the programmed start time.

For example, the current time is 10:00 PM. The user sets the delayed start time for 10:30 PM and then powers off the pump. The user does not power on the pump until 10:40 PM, which is 10 minutes AFTER the scheduled start time. The pump does not begin delivery until 10:30 PM the FOLLOWING day.

Dose Interruptions

To understand how interruptions are resolved, compare the normal Intermittent delivery diagram below for a 1-hour dose every four hours, to the interruption diagrams shown in the examples:

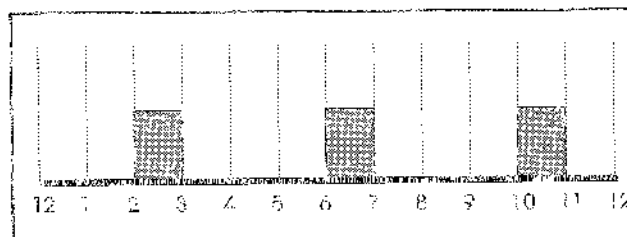


Dose Delivered -

- If an Intermittent dose is interrupted and the pump is placed in the RUN mode during the on or off cycle of the current dose, the dose continues delivery from the point where it was interrupted and all remaining doses are shifted.

For example, in the following diagram, an intermittent dose was started at 2:00. The dose was interrupted during the current dose cycle and was started again before the end of the current dose cycle.


5. Intermittent

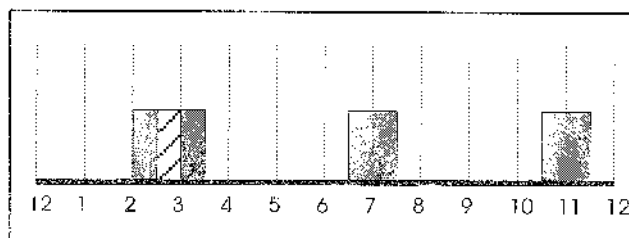

Dose Delivered -  Dose Interrupted - 

- If an Intermittent dose is interrupted and the pump is placed in the RUN mode after the off cycle of the current dose, the next dose begins delivery. All remaining doses are shifted if necessary.

For example, in the following diagram, an intermittent dose was started at 2:00. It was interrupted during the current delivery cycle at 2:30, but was not started again until after the current delivery cycle at 6:45:

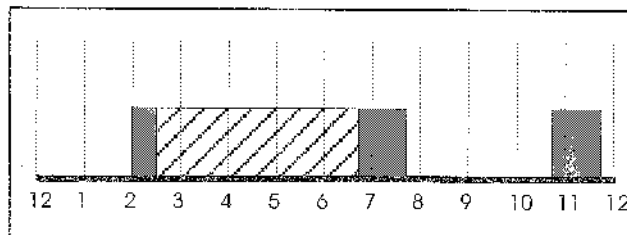

Dose Delivered -  Dose Interrupted - 

- When the dose is complete, the start time for the Next Dose displays. To change the start time refer to Chapter 10 Options Menu, page 97.
- Any missed Intermittent doses, except the first dose, are added to the end of the delivery schedule.

5. Intermittent

Changing an Intermittent Program

Can be changed at any time	Can be changed if no dose is in progress	Cannot be changed after  is pressed
• Dose parameters • KVO rate • Container size • Air sensitivity • Call back alarm	• Next dose start time	• Current time



Note: You cannot change a program when the pump is in FULL or CONTAINER Lock.

To change a program:

- 1 From the STOP mode, press [CHANGE] and select CHANGE PROGRAM.
- 2 As each message displays, make the desired changes or press [YES/ENTER] to accept the current setting.
- 3 Review the program. Press  to begin review.
- 4 Press [START] to begin delivery.

GemStar® Intermittent Programming Worksheet

GIVE: _____

Parameters	Selections and/or Values to be Entered	Programmable Ranges
Therapy	Intermittent	N/A
Dose Amount	_____	0.1 - 9600 mL
Infusion Time Per Dose	_____ hrs:min	1 min - 24 hrs
Dose Frequency	_____ hrs:min	(Infusion time per dose + 1 min) to 96 hrs
Number of Doses Per Container	_____	1 - 999
KVO (Keep Vein Open)	☐ No KVO Rate: 0.0 ☐ Yes KVO Rate: _____	0.1 - 5.0 mL/hr
Container Size*	_____	0.1 - 9999 mL
Air Sensitivity	☐ On Alarms at approx. 0.5 mL of air ☐ 2 mL Alarms at approx. 2 mL of air ☐ Off No alarm	N/A
Call Back Alarm	☐ No ☐ Yes Alarm: _____ hrs:min before start	1 min to time between doses
Start Time	☐ Start immediately after programming ☐ Start time: _____ Specify AM or PM if not 24-hour clock.	12:00 AM - 11:59 PM or 00:00 - 23:59

* The pump automatically calculates the minimum container size required to deliver the number of programmed doses and KVO between doses. If KVO is desired before a delayed start, increase the container size to include the desired KVO volume.

Notes

Continuous

Therapy Features and Specifications

This therapy allows programming of a continuous delivery.

- Optional KVO (Keep Vein Open) can be set during programming.
- The rate can be titrated while the pump is infusing.
- An optional secondary container can be “piggybacked” on the primary container during regular delivery using a GemStar Piggyback Pump Set. Refer to *Piggybacking* in this chapter for more information.



Tip: The Programming Worksheet at the end of this chapter lists the minimum and maximum values allowed for this protocol. Use copies of the programming worksheet to record program values before programming the pump.

Continuous Programming Tips

The following table explains how to use the keypad during continuous programming:

Key	Function	Key	Function
[0] - [9]	Selects menu items and enters numeric values.	 	Scrolls through program review and help messages. The up arrow inserts a decimal point in a numeric value.
	Accepts entries or continues when done viewing a display.		Clears a numeric value (resets to zero) during programming.
	Displays information during programming and alarm conditions.		Returns to a previous display, exits the Options Menu, and exits the HELP mode.

6. Continuous

Continuous Programming

- 1 Select CONTINUOUS from the Programming Menu.
- 2 Select the unit of measure and concentration as required.
- 3 Enter the rate.
- 4 Enter the VTBI (Volume to be infused).
- 5 Enter the KVO rate, if desired.

KVO runs after the programmed VTBI until the total container is infused. If no KVO is desired, press [YES/ENTER] while the value displayed is zero.

- 6 Enter the container size (VTBI + KVO).
VTBI is the default value on the CONTAINER display.

Note: If KVO is programmed, increase the container size to include the KVO delivery.

To enter the total in mL when programming a concentration, press  to move to the mL line. Enter the mL size when mL flashes on the display.

- 7 Select the air sensitivity, if required.
- 8 Review the program. Press  to begin review.
- 9 If using the PURGE key to prime the pump set, press [PURGE] now.

WARNING: ALWAYS disconnect the pump set from the patient before priming or purging.

Rate Titration

To titrate the rate while the pump is infusing:

- 1 In the RUN mode, press the number key for the new rate. Press [YES/ENTER].
For example: if the current rate is 100 mL/hr, press [5] [0] [YES/ENTER] to change the rate to 50 mL/hr.

Note: You have approximately five seconds after pressing a key to press another numeric key, the decimal key, or [YES/ENTER]. If no key is pressed within five seconds, the rate does not change and the delivery continues as previously programmed.

- 2 NEW RATE displays with the new value.
 - Press [YES/ENTER] within 10 seconds to accept the change and begin delivery at the new rate.
 - Press [NO] to continue delivery at the previously programmed rate.

Piggybacking

Piggybacking with the GemStar Pump requires a GemStar Primary Piggyback Pump Set with a proximal Y-site backcheck valve and a secondary piggyback set with an extension hook. **The secondary container must be higher than the primary container.** Refer to the instructions included with the piggyback pump set for more information.

WARNING: When using the pump for secondary delivery (piggybacking), ensure the fluids being infused are chemically and physically compatible.

The piggyback rate and piggyback VTBI display until the secondary container has been delivered. Then the primary infusion continues.

Note: The piggyback infusion amount is NOT added to the primary programmed infusion amount.


6. Continuous

Preparing for Secondary Delivery (Piggybacking)

- 1 Ensure a primary piggyback pump set with a proximal Y-site backcheck valve is in place.
- 2 Use the extension hook from the secondary set to suspend the primary container from the I.V. pole.
- 3 Prepare the secondary set according to the instructions included with the set.
- 4 Attach a needle or blunt cannula (if appropriate) and prime the secondary set.
- 5 Suspend the secondary container from the I.V. pole.

Note: When using a secondary container of 500 mL or greater, ensure the bottom of the secondary container is at least 7 inches (17.8 cm) above the fluid level in the primary container. Use additional extension hooks if necessary.

- 6 Connect the secondary set to the proximal Y-site (or proximal luer port).

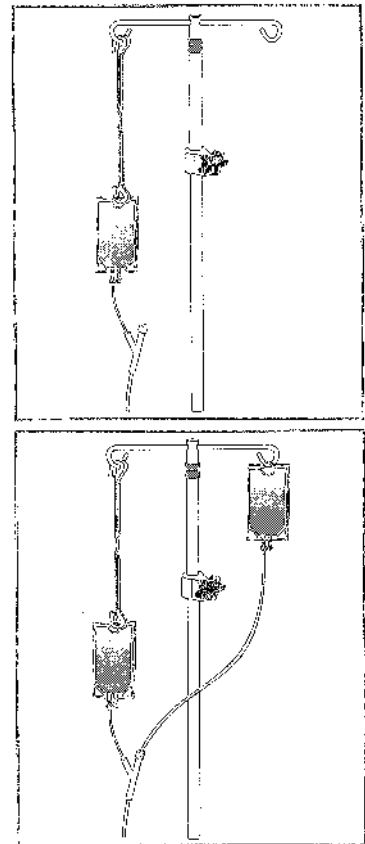

Programming a Secondary Delivery (Piggybacking)

Note: To distinguish the piggyback delivery from the primary delivery, PIGGYBACK displays and the piggyback VTBI counts down to completion. For example, a piggyback rate of 250 mL/hr with a remaining VTBI of 218.5 mL will display as shown.

PIGGYBACK
RATE 250.0 mL/hr
VTBI 218.5 mL



Tip: Prepare and set up containers BEFORE programming the pump for piggyback delivery.

- 1 Press [CHANGE] and select PIGGYBACKING from the Change Menu.
- 2 Enter the dose amount for the secondary container.
- 3 Enter the infusion time for the secondary container.
- 4 Open the slide clamp or CAIR® (roller) clamp on the tubing of the secondary container.
- 5 Press [YES/ENTER] to begin piggyback delivery.



Note: Primary flow stops until the secondary VTBI is delivered; then primary flow continues. Any remaining secondary fluid will be delivered at the primary rate. The maximum piggyback rate is 300 mL/hr. For secondary rates above 300 mL/hr, flow may occur from the primary container.

Changing a Piggyback Delivery

- 1 Press [CHANGE] and select PIGGYBACKING from the Change Menu.
- 2 Press [CHANGE] to reset the dose amount to zero.
- 3 Enter the new piggyback dose amount and press [YES/ENTER].
- 4 Set the infusion time for the secondary container.
- 5 Press [YES/ENTER] to begin piggyback delivery.



Note: You must know the amount that has already been delivered. When the piggyback dose amount is changed, the pump begins a new VTBI countdown.

8. Continuous

Canceling a Piggyback Delivery

- 1 Press [CHANGE] and select **PIGGYBACKING** from the Change Menu.
- 2 Press [CHANGE] to reset the dose amount to zero.
- 3 Press [YES/ENTER] to cancel piggyback delivery.

Changing a Continuous Program

Can be changed at any time		Cannot be changed
• Delivery rate	• Container size	• Unit of measure
• VTBI	• Air sensitivity	• Concentration
• KVO rate	• Piggybacking	



Note: You cannot change a program when the pump is in FULL or CONTAINER Lock.

To change a program:

- 1 From the STOP mode, press [CHANGE] and select **CHANGE PROGRAM**.
- 2 As each message displays, make the desired changes or press [YES/ENTER] to accept the current setting.
- 3 Review the program. Press **⏮** to begin review.
- 4 Press [START] to begin delivery.

6. Continuous

GemStar® Continuous Programming Worksheet

GIVE: _____

Parameter	Selections and/or Values to be Entered	Programmable Ranges
Therapy	Continuous	N/A
Unit of Measure	☐ Concentration: _____ mg/mL ☐ Concentration: _____ mcg/mL ☐ mL	0.1 - 1000 mg/mL 1 - 1000 mcg/mL
Rate	_____	0.1 - 9999 mg/hr 1 - 999900 mcg/hr 0.1 - 1000 mL/hr
VTBI (Volume To Be Infused)	_____	0.1 - 99999 mg 1 - 9999000 mcg 0.1 - 9999 mL
KVO (Keep Vein Open)	☐ No KVO Rate: 0.0 ☐ Yes KVO Rate: _____	0.1 - 5.0 mL/hr
Container Size (VTBI + KVO)	_____	0.1 - 99999 mg 1 - 9999000 mcg 0.1 - 9999 mL
Air Sensitivity	☐ On Alarms at approx. 0.5 mL of air ☐ 2 mL Alarms at approx. 2 mL of air ☐ Off No alarm	N/A
Piggyback Dose Amount	_____	1 - 7200 mL
Piggyback Infusion Time	_____ hrs:min	1 min or dose amount / 300 mL/hr to 24 hrs or dose amount / 0.1 mL/hr

Notes

Weight Dosed

Therapy Features and Specifications

This therapy allows programming of delivery rate as a function of body weight.

- ♦ The air sensitivity in a Weight Dosed program is automatically set to ON. This setting can be changed through the Options Menu after programming is complete. Refer to *Air Sensitivity* in Chapter 10 Options Menu, on page 107.
- ♦ The pump automatically delivers KVO (Keep Vein Open) after the programmed VTBI (Volume To Be Infused) has completed.
- ♦ The rate can be titrated while the pump is infusing.



Tip: The Programming Worksheet at the end of this chapter lists the minimum and maximum values allowed for this protocol. Use copies of the programming worksheet to record program values before programming the pump.

WARNING: When using a vented drip set, ALWAYS use an air-eliminating filter OR set the air sensitivity to either ON or 2 mL.

CAUTION: When infusing short-half-life drugs (i.e., those dosed with units of mcg/kg/min) at very low rates (0.8 mL/hr or less) the rate may be too slow for that drug. Consider using a lower drug concentration with these drugs so a pump rate above 0.8 mL/hr can be used.

7. Weight Dosed

Weight Dosed Programming Tips

The following table explains how to use the keypad during weight dosed programming:

Key	Function	Key	Function
[0] - [9]	Selects menu items and enters numeric values.	 	Scrolls through program review and help messages. The up arrow inserts a decimal point in a numeric value.
	Accepts entries or continues when done viewing a display.		Clears a numeric value (resets to zero) during programming.
	Displays information during programming and alarm conditions.		Returns to a previous display, exits the Options Menu, exits the HELP mode.

Weight Dosed Programming

- 1 Select WEIGHT DOSED from the Programming Menu.
- 2 Select the unit of measure.

Note: If two numbers appear before a menu time (e.g., 02), press both number keys to select that item.

- 3 Enter the patient weight, if required.
- 4 Enter the concentration:
 - Enter the value for the unit of measure when it flashes on the display.
 - Press [YES/ENTER] to confirm.

Tip: Press  to toggle between units of measure.

- 5 Enter the rate.

Note: The pump automatically displays the mL/hr equivalent even when another unit of measure is selected.

- 6 Enter the VTB!
- 7 Review the program. Press  to begin review.

7. Weight Dosed

- 8 If using the PURGE key to prime the pump set, press [PURGE] now.

WARNING: ALWAYS disconnect the pump set from the patient before priming or purging.

Auto KVO

When the programmed VTBI has completed delivery, Auto KVO delivery begins automatically and the empty container alarm sounds.

Note: Auto KVO delivers even while the empty container alarm sounds. Press [SILENCE] to quiet the alarm for 2 minutes.

If the weight dosed rate is:	Auto KVO delivers at:
Greater than or equal to 1 mL/hr	1 mL/hr
Less than 1 mL/hr	The weight dosed rate

Rate Titration

To titrate the rate while the pump is infusing:

- 1 In the RUN mode, press the number keys for the new rate. Press [YES/ENTER].

For example: if the current rate is 4.0 mg/kg/hr, press [5] [YES/ENTER] to change the rate to 5.0 mg/kg/hr.

Press  to enter the rate in mL/hr. When using units other than mL/hr, the pump automatically adjusts the mL/hr value.

Note: You have approximately five seconds after pressing a key to press another numeric key, the decimal key, or [YES/ENTER]. If no key is pressed within five seconds, the rate does not change and the delivery continues as previously programmed.

- 2 NEW RATE displays with the new value.
- Press [YES/ENTER] within 10 seconds to accept the change and begin delivery at the new rate.
 - Press [NO] to continue delivery at the previously programmed rate.

7. Weight Dosed

VTBI Titration

To titrate the VTBI while the pump is infusing:

- 1 Press [CHANGE] and select TITRATE VTBI.
- 2 Press the number keys for the new VTBI. Press [YES/ENTER].

For example: if the current VTBI is 50 mL, press [1][0][0][YES/ENTER] to change the VTBI to 100 mL.

Note: You have approximately five seconds after pressing a key to press another numeric key, the decimal key, or [YES/ENTER]. If no key is pressed within five seconds, the VTBI does not change and the delivery continues as previously programmed.

- 3 NEW VTBI displays with the new rate.
 - Press [YES/ENTER] within 10 seconds to accept the change and begin delivery with the new VTBI.
 - Press [NO] to continue delivery at the previously programmed VTBI.

Changing a Weight Dosed Program

Can be changed at any time	Cannot be changed
• Delivery rate • VTBI • Patient weight	• Unit of measure • Concentration (amount of drug)

Note: You cannot change a program when the pump is in FULL or CONTAINER Lock.

To change a program:

- 1 From the STOP mode, press [CHANGE] and select CHANGE PROGRAM.
- 2 As each message displays, make the desired changes or press [YES/ENTER] to accept the current setting.
- 3 Review the program. Press  to begin review.
- 4 Press [START] to begin delivery.

7. Weight Dosed

GemStar® Weight Dosed Programming Worksheet

GIVE: _____

Parameter	Selections and/or Values to be Entered	Programmable Ranges
Therapy	Weight Dosed	N/A
Unit of Measure**	☐ mcg/kg/min ☐ mcg/hr ☐ mcg/kg/hr ☐ grams/hr ☐ mcg/min ☐ ng/kg/min ☐ mg/kg/min ☐ U/kg/hr ☐ mg/kg/hr ☐ U/min ☐ mg/min ☐ U/hr ☐ mL/hr ☐ mU/min ☐ mg/hr ☐ mEq/hr	N/A
Patient Weight	(When programming in kg/min or kg/hr) _____	2.0 - 200.0 kg
Concentration (2 entries are required. Use  to toggle)	(For all units of measure EXCEPT mL/hr) _____ mg _____ U _____ mEq _____ in _____ in _____ in _____ mL _____ mL _____ mL	0.1 - 9999000 mg 0.1 - 99999999 U 0.1 - 9999 mEq 0.1 - 9999 mL
Rate	_____	(see below *)
VTBI (Volume To Be Infused)	_____	0.1 - 9999 mL

* Programmable Ranges for Rate (not all units of measure are available in all pump configurations):

0.001 - 833 mg/kg/min	0.1 - 1000 mL/hr	0.001 - 99999999 U/min
0.001 - 500000 ng/kg/hr	0.001 - 999999 mg/hr	0.001 - 99999999 U/hr
0.001 - 16666 mg/min	0.001 - 9999999 mcg/hr	0.001 - 999999999 mU/min
0.001 - 166666 mcg/kg/min	0.001 - 999 grams/hr	0.001 - 999 mEq/hr
0.001 - 9999999 mcg/kg/hr	0.009 - 999 ng/kg/hr	
0.002 - 166666 mcg/min	0.001 - 99999999 U/kg/hr	

** If two numbers appear before a menu item (e.g., 02), press both number keys to select that item.

CAUTION: When infusing short-half-life drugs (i.e., those dosed with units of mcg/kg/min) at very low rates (0.8 mL/hr or less) the rate may be too slow for that drug. Consider using a lower drug concentration so a pump rate above 0.8 mL/hr can be used.

Notes

Chapter 8

mL/hr Only

Therapy Features and Specifications

This therapy allows quick programming of delivery rate in milliliters (mL) only.

- ♦ The air sensitivity is automatically set to ON. This setting can be changed through the Options Menu after programming is complete. Refer to *Air Sensitivity* in Chapter 10 Options Menu, on page 107.
- ♦ The pump automatically delivers KVO (Keep Vein Open) after the programmed VTBI (Volume To Be Infused) is complete.
- ♦ The rate can be titrated while the pump is infusing.
- ♦ An optional secondary container can be “piggybacked” onto the primary container during regular delivery using a GemStar Piggyback Pump Set. Refer to *Piggybacking* in this chapter for more information.
- ♦ A program review is not required and is available only through the Options Menu.



Tip: The Programming Worksheet at the end of this chapter lists the minimum and maximum values allowed for this protocol. Use copies of the programming worksheet to record program values before programming the pump.

WARNING: When using a vented set, **ALWAYS** use an air-eliminating filter **OR** set the air sensitivity to either ON or 2 mL.

8. mL/hr Only

mL/hr Programming Tips

The following table explains how to use the keypad during mL/hr programming:

Key	Function	Key	Function
[0] - [9]	Selects menu items and enters numeric values.	 	Scrolls through program review and help messages. The up arrow inserts a decimal point in a numeric value.
	Accepts entries or continues when done viewing a display.		Clears a numeric value (resets to zero) during programming.
	Displays information during programming and alarm conditions.		Returns to a previous display, exits the Options Menu, exits the HELP mode.

mL/hr Only Programming

- 1 Select mL/hr Only from the Programming Menu.
- 2 Enter the rate.
- 3 Enter the VTBI.
- 4 If using the PURGE key to prime the pump set, press [PURGE] now.

WARNING: ALWAYS disconnect the pump set from the patient before priming or purging.



Note: A program review is not required for mL/hr Only programs.

Auto KVO

When the programmed VTBI has completed delivery, Auto KVO delivery begins automatically and the empty container alarm sounds.

Note: Auto KVO delivers even while the empty container alarm sounds. Press [SILENCE] to quiet the alarm for 2 minutes.

If the weight dosed rate is:	Auto KVO delivers at:
Greater than or equal to 1 mL/hr	1 mL/hr
Less than 1 mL/hr	The programmed rate

Rate Titration

To titrate the rate while the pump is infusing:

- 1 In the RUN mode, press the number keys for the new rate. Press [YES/ENTER].
For example, if the current rate is 100 mL/hr, press [5] [0] [YES/ENTER] to change the rate to 50 mL/hr.

Note: You have approximately five seconds after pressing a key to press another numeric key, the decimal key, or [YES/ENTER]. If no key is pressed within five seconds, the rate does not change and the delivery continues as previously programmed.

- 2 NEW RATE displays with the new value.
 - Press [YES/ENTER] within 10 seconds to accept the change and begin delivery at the new rate.
 - Press [NO] to continue delivery at the previously programmed rate.

8. mL/hr Only

VTBI Titration

To titrate the VTBI while the pump is infusing:

- 1 Press [CHANGE] and select TITRATE VTBI.
- 2 Press the number keys for the new VTBI. Press [YES/ENTER].

For example: if the current VTBI is 50 mL, press [1][0][0][YES/ENTER] to change the VTBI to 100 mL.

Note: You have approximately five seconds after pressing a key to press another numeric key, the decimal key, or [YES/ENTER]. If no key is pressed within five seconds, the VTBI does not change and the delivery continues as previously programmed.

- 3 NEW VTBI displays with the new rate.
 - Press [YES/ENTER] within 10 seconds to accept the change and begin delivery with the new VTBI.
 - Press [NO] to continue delivery at the previously programmed VTBI.

Piggybacking

Piggybacking with the GemStar Pump requires a GemStar Primary Piggyback Pump Set with a proximal Y-site backcheck valve and a secondary piggyback set with an extension hook. The secondary container must be higher than the primary container. Refer to the instructions included with the piggyback pump set for more information.

WARNING: When using the pump for secondary delivery (piggybacking), ensure the fluids being infused are chemically and physically compatible.

The piggyback rate and piggyback VTBI display until the secondary container has been delivered. Then the primary infusion continues.



Note: The piggyback infusion amount is NOT added to the primary programmed infusion amount.

Preparing for Secondary Delivery (Piggybacking)

- 1** Ensure a primary piggyback pump set with a proximal Y-site backcheck valve is in place.
- 2** Use the extension hook from the secondary set to suspend the primary container from the I.V. pole.
- 3** Prepare the secondary set according to the instructions included with the set.
- 4** Attach a needle or blunt cannula (if appropriate) and prime the secondary pump set.
- 5** Suspend the secondary container from the I.V. pole.

Note: When using a secondary container of 500 mL or greater, ensure the bottom of the secondary container is at least 7 inches (17.8 cm) above the fluid level in the primary container. Use additional extension hooks if necessary.

- 6** Connect the secondary set to the proximal Y-site (or proximal luer port).

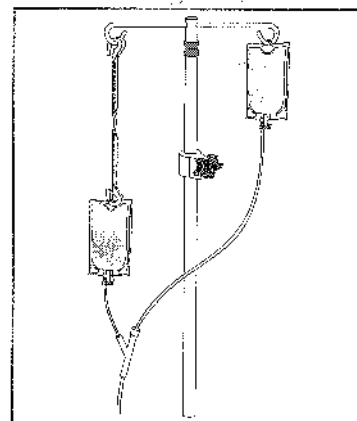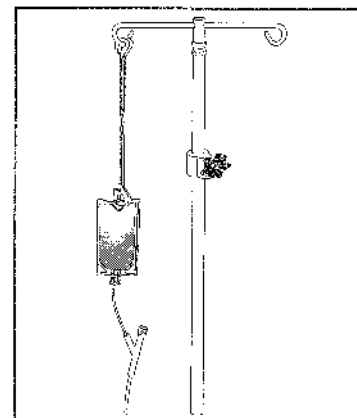

8. mL/hr Only

Programming a Secondary Delivery (Piggybacking)

Note: To distinguish the piggyback delivery from the primary delivery, PIGGYBACK displays and the piggyback VTBI counts down to completion. For example, a piggyback rate of 250 mL/hr with a remaining VTBI of 218.5 mL will display as shown.

PIGGYBACK	
RATE	250.0 mL/hr
VTBI	218.5 mL
▶▶▶▶▶	

Tip: Prepare and set up containers BEFORE programming the pump for piggyback delivery.

- 1 Press [CHANGE] and select PIGGYBACKING from the Change Menu.
- 2 Enter the dose amount for the secondary container.
- 3 Enter the infusion time for the secondary container.
- 4 Open the slide clamp or CAIR® (roller) clamp on the tubing of the secondary set.
- 5 Press [YES/ENTER] to begin piggyback delivery.

Note: Primary flow stops until the secondary VTBI is delivered; then primary flow continues. Any remaining secondary fluid will be delivered at the primary rate. The maximum piggyback rate is 300 mL/hr. For secondary rates above 300 mL/hr, flow may occur from the primary container.

Changing a Piggyback Delivery

- 1 Press [CHANGE] and select PIGGYBACKING from the Change Menu.
- 2 Press [CHANGE] to reset the dose amount to zero.
- 3 Enter the new piggyback dose amount and press [YES/ENTER].
- 4 Set the infusion time for the secondary container.
- 5 Press [YES/ENTER] to begin piggyback delivery.

Note: You must know the amount that has already been delivered. When the piggyback dose amount is changed, the pump begins a new VTBI countdown.

Canceling a Piggyback Delivery

- 1 Press [CHANGE] and select PIGGYBACKING from the Change Menu.
- 2 Press [CHANGE] to reset the dose amount to zero.
- 3 Press [YES/ENTER] to cancel piggyback delivery.

Changing a mL/hr Only Program

Can be changed at any time

- | | |
|------------------|----------------|
| • Delivery rate | • Piggybacking |
| • Container size | • VTBI |



Note: You cannot change a program when the pump is in FULL or CONTAINER Lock.

To change a program:

- 1 From the STOP mode, press [CHANGE] and select CHANGE PROGRAM.
- 2 As each message displays, make the desired changes or press [YES/ENTER] to accept the current setting.
- 3 Press [START] to begin delivery.



Note: A program review is not required and is available only through the Change or Options Menu.

6. mL/hr Only

GemStar® mL/hr Only Programming Worksheet

GIVE: _____

Parameter	Selections and/or Values to be Entered	Programmable Ranges
Therapy	mL/hr Only	N/A
Rate	_____	0.1 - 1000 mL/hr
VTBI (Volume To Be Infused)	_____	0.1 - 9999 mL
Piggyback Dose Amount	_____	1 - 7200 mL
Piggyback Infusion Time	_____ hrs:min	1 min or dose amount / 300 mL/hr to 24 hrs or dose amount / 0.1 mL/hr

Variable Time

Therapy Features and Specifications

This therapy allows programming of multiple doses at specified times.

- Optional KVO (Keep Vein Open) can be set during programming.
- The default bolus delivery rate is 400 mL/hr.



Tip: The Programming Worksheet at the end of this chapter lists the minimum and maximum values allowed for this protocol. Use copies of the programming worksheet to record program values before programming the pump.

Variable Time Programming Tips

The following table explains how to use the keypad during variable time programming:

Key	Function	Key	Function
[0] - [9]	Selects menu items and enters numeric values.	 	Scrolls through program review and help messages. Toggles between AM and PM in time-based entries. The up arrow inserts a decimal point in a numeric value.
	Accepts entries or continues when done viewing a display.		Clears a numeric value (resets to zero) during programming.
	Displays information during programming and alarm conditions.		Returns to a previous display, exits the Options Menu, and exits the HELP mode.

9. Variable Time

Variable Time Programming

- 1 Select VARIABLE TIME from the Programming Menu.
- 2 Select the phase type.
- 3 Select the unit of delivery.
- 4 Change the displayed time if necessary.
- 5 Enter the optional base rate with start and stop time.
- 6 If programming in Rates or Doses, go to step 7. If programming in Percentages, set the daily amount. Then go to step 8.

Note: If programming in percentages, the total daily amount of all phases programmed must equal 100%.

- 7 Enter the phase information for each phase. Up to 24 phases can be programmed.

Note: When entering time in 12-hour format, press  or  to toggle between AM and PM.

- 8 Enter the optional bolus dose and its lockout.
- 9 Enter the optional KVO.
- 10 Enter the container size.

The pump calculates the minimum container size required to deliver the phases, base rate, and KVO for one day or in a 24-hour period. If delivering the protocol for more than one day, multiply this amount by the number of days of therapy.

Note: Increase the container size to include bolus deliveries.

- 11 Select the air sensitivity, if required.

9. Variable Time

12 Review the program. Press  to begin review.

13 If using the PURGE key to prime the pump set, press [PURGE] now.

WARNING: ALWAYS disconnect the pump set from the patient before priming or purging.

Note: With variable time programming, the pump keeps track of where it is in the program each minute of the entire 24-hour period. For example, if Phase 1 is missed, and the pump is placed in RUN mode in the middle of Phase 2, it starts delivery according to the program parameters from the middle of Phase 2.

Variable Time Delivery Interruptions

If a phase is interrupted and the pump is placed in the RUN mode:	Then:
Before the end of the phase.	The phase runs to its scheduled stop time.
After the end of the phase, but before the beginning of the next start time.	The next phase begins delivery at the scheduled start time. All applicable Base Rate or KVO delivery infuses while waiting for the start of the next phase.
In the middle of a phase and the phase start was missed because of a program interruption or a missed program start.	The pump immediately begins delivery of the current phase at its current time-based point.

9. Variable Time

Changing a Variable Time Program

Can be changed at any time	Can ONLY be changed before  is pressed	Cannot be changed
• KVO rate • Container size • Air sensitivity • Bolus	• Current time • Base rate • Daily amount • Phase parameters	• Unit of measure • Concentration • Phase type

 *Note: You cannot change a program when the pump is in FULL or CONTAINER Lock.*

To change a program:

- 1 From the STOP mode, press [CHANGE] and select CHANGE PROGRAM.
- 2 As each program parameter displays, make the desired changes or press [YES/ENTER] to accept the current setting.
- 3 Review the program. Press  to begin review.
- 4 Press [START] to begin delivery.

 *Note: The container size cannot be made smaller than the daily total of phases, base rate, and KVO combined; or the current amount delivered, whichever is greater.*

GemStar® Variable Time Programming Worksheet

GIVE:.

Parameter	Selections and/or Values to Be Entered	Programmable Ranges
Select Phase Type	☐ Percentages ☐ Rates ☐ Doses	0.1 - 100%* 0.1 - 400 mL/hr 0.1 - 9600 mL
Unit of Measure	☐ Concentration: _____ mg/mL ☐ Concentration: _____ mcg/mL ☐ mL	0.1 - 1000 mg/mL 1 - 1000 mcg/mL
Current Time	Confirm or correct current time.	12:00 AM - 11:59 PM or 00:00 - 23:59
Base Rate	☐ No ☐ Yes Start time: _____ Stop time: _____ Rate: _____	12:00 AM - 11:59 PM or 00:00 - 23:59 0.1 - 400 mL/hr
Total Daily Amount for All Phases: _____ (required only when programming in percentages) _____ SEE ATTACHED WORKSHEET _____		
Bolus Dose	_____ _____	0.1 mL - 25.0 mL 0.1 mg - 2500 mg 1 mcg - 25000 mcg
Bolus Lockout?	_____ minutes	5 - 999 min
KVO (Keep Vein Open)	☐ No KVO Rate: 0.0 ☐ Yes KVO Rate: _____	0.0 - 5.0 mL/hr
Container Size (Total Daily Amount X Number of Days of Therapy)	_____	0.1 - 99999 mg 1 - 9999000 mcg 0.1 - 9999 mL
Air Sensitivity	☐ On Alarms at approx. 0.5 mL of air ☐ 2 mL Alarms at approx. 2 mL of air ☐ Off No alarm	N/A

*Sum of all phases must equal 100% when programmed in percentages.

9. Variable Time

Phone	Start Time	Stop Time	Percentage	Rate	Dose
1					
2					
3					
4					
5					
6					
7					
8					
9					
10					
11					
12					
13					
14					
15					
16					
17					
18					
19					
20					
21					
22					
23					
24					
TOTAL	Note: The sum of all phones must equal 100% when programmed in percentages.				

Options Menu

Overview

- 1** Press [OPTIONS] from RUN or STOP mode
- 2** Select the number key for the desired function.
- 3** Follow the display messages for the selected function. Refer to the appropriate description in this chapter for more information.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE

Note: Menu selections may vary depending on the configuration of the pump. For example, NEXT DOSE is only available when an Intermittent therapy is programmed.

When the pump exits OPTIONS, the pump returns to the mode it was in when [OPTIONS] was pressed.

Tips for Using the Options Menu

The following keys are helpful when using the Options Menu:

Keys	Functions
[1] - [8]	Press the corresponding number key to select a menu function.
	Press  to accept changes or to continue when done viewing a display.
 	Press the up or down arrow to scroll through long display messages or to toggle AM and PM in time-based entries.
	Press  to return to a previous display or to exit OPTIONS.

Tip: A menu item does not have to be on the display to be selected. For example, only items 1 through 4 are shown when the Options Menu is first displayed. To select the Print Menu without scrolling the display, press [6].

10. Options Menu

Select this item from the Options Menu:	To do this function:	Comments:
REVIEW PROGRAM	Display the current program settings.	N/A
HISTORIES	Display the history for the latest 400 events, shift or volume totals, bolus history (if available), or current configuration settings. Clear the history. Display the Operation Test.	N/A
KEYPAD LOCK	Lock and unlock the keypad to prevent unauthorized changes to the program.	The lock levels allow the clinician to restrict access to various levels of pump operation.
PUMP SETTINGS	Change the air sensor, occlusion sensor, and sound level settings.	N/A
SET CLOCK	Set the clock.	The clock can be displayed in 12-hour or 24-hour format.
PRINT	Print the history event log, bolus history (if available), current program settings, and speed protocols.	Files may be downloaded to a printer or personal computer.
SPEED PROTOCOLS	Assign, review/recall, and delete the stored protocols.	Allows up to 9 speed protocols to be saved in pump memory.
NEXT DOSE	Change the start time of the next dose when an interruption has caused a missed dose in an intermittent program.	NEXT DOSE appears on the menu only when the pump is in an intermittent program and changing the next dose is allowed.

10. Options Menu

Review Program

- 1 Press [OPTIONS] and select REVIEW PROGRAM. The pump displays PROGRAM REVIEW.
- 2 Press  to begin review.
- 3 Press [YES/ENTER] or [BACK-UP] when done to exit OPTIONS.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE

Histories

- 1 Press [OPTIONS] and select HISTORIES.
- 2 Select the number key for the desired function.
- 3 Complete the appropriate steps for the selected function.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE

10. Options Menu

Select this item from the Histories Menu:	To do this function:	Display includes:
HISTORY	Display history event log of the latest 400 events.	• Current date and time • Date and time the program was cleared, a new container was selected, or a speed protocol was recalled • History event log of the latest 400 events
SHIFT TOTALS	Display or clear shift totals.	• Time and date when shift was cleared • Volume infused in the shift • Any relevant program-specific information
VOLUME TOTALS	Display volume totals.	• Amount infused • VTBI • Any relevant program-specific information
CONFIGURATION	Display pump configuration.	• Dependent upon pump configuration
CLEAR HISTORY	Clear history event log and bolus history.	• CLEARING HISTORY displays for 3 to 4 seconds
BOLUS HISTORY	Display bolus history in hourly increments for the last 48 hours.	• Current date and time • Number of boluses delivered • Number of bolus demands • Amount delivered through bolus delivery
OPERATION TEST	Perform user-assisted Operation Test.	• Messages that guide the user through each section of the test



Note: Some menu items may not be available depending on the configuration of the pump.



Tip: Press and HOLD  or  to scroll through the history displays faster.

10. Options Menu

History

To display the history event log:

- 1 Press [OPTIONS] and select HISTORIES.
- 2 Select HISTORY from the Histories Menu.
- 3 Press  or  to scroll the display.
- 4 Press [BACK-UP] when done to exit OPTIONS.

Shift Totals

To display the shift totals:

- 1 Press [OPTIONS] and select HISTORIES.
- 2 Select SHIFT TOTALS from the Histories Menu.
- 3 The shift totals display. Press  or  to scroll the display.

VERSION N.NNN
S/N XXXXXXXX
TIME IS 09:00 AM
FRI, JUN 11, 04

NEW CONTAINER
08:15AM JUN 11

HISTORY CLEARED
11:00PM JUN 10
EVENT LOG
PM11:18 SET ON
AIR SENSITIVITY
PM11:18 SET CONT
SIZE 100 mL
PM11:18 SET RATE
1.0 mL/hr
PM11:18 SHIFT
CLEARED
END OF EVENT LOG

Sample History Event Log

To clear the shift totals:

- 1 Follow the steps for displaying shift totals.
- 2 At the end of the shift totals, the pump displays CLEAR TOTALS? YES OR NO:
 - Press [YES/ENTER] to clear the shift totals.
 - Press [NO] to exit OPTIONS without clearing the shift totals.

END OF
SHIFT TOTALS
CLEAR TOTALS?
YES OR NO 



Tip: Press [HELP] to display shift totals when the pump is not in PROGRAMMING mode or an alarm condition.

10. Options Menu

Volume Totals

To display the volume totals:

- 1 Press [OPTIONS] and select HISTORIES.
- 2 Select VOLUME TOTALS from the Histories Menu.
- 3 Press [BACK-UP] when done to exit OPTIONS.

To clear the volume totals:

- 1 Follow the steps for displaying volume totals.
- 2 At the end of the volume totals, the pump displays CLEAR TOTALS? YES OR NO:
 - Press [YES/ENTER] to clear the volume totals.
 - Press [NO] to exit OPTIONS without clearing the volume totals.

END OF
VOLUME TOTALS
CLEAR TOTALS?
YES OR NO

Configuration

To display the pump configuration:

- 1 Press [OPTIONS] and select HISTORIES.
- 2 Select CONFIGURATION from the Histories Menu.
- 3 Press [BACK-UP] when done to exit OPTIONS.

Clear History

To clear the history event log and bolus history:

- 1 Press [OPTIONS] and select HISTORIES.
- 2 Select CLEAR HISTORY from the Histories Menu.
- 3 The pump display CLEARING HISTORY for 3–4 seconds before exiting OPTIONS.



Tip: Press [BACK-UP] to exit OPTIONS immediately.

10. Options Menu

Bolus History

To display the bolus history:

- 1** Press [OPTIONS] and select HISTORIES.
- 2** Select BOLUS HISTORY from the Histories Menu.
- 3** Press [BACK-UP] when done to exit OPTIONS.

TIME IS 09:15 AM
FRI, JUN 11, 04
BOL DEL 30
BOL DEM 37
BOL AMT 30.0 mL
06:10A TO 09:14A

Sample Bolus History

Note: Bolus History is only available if the current program has a bolus programmed, or had a bolus programmed at one time. The Bolus History is maintained in hourly increments for the latest 48 hours. The Bolus History is cleared when a new program is entered, a Speed Protocol is recalled, or the History is cleared.

Operation Test

The Operation Test checks several pump functions to ensure the GemStar Pump is working properly. Refer to Chapter 14 Operation Test on page 143 for instructions.

10. Options Menu

Keypad Lock

- 1 Press [OPTIONS] and select KEYPAD LOCK from the Options Menu.
- 2 Select the desired lock level.

The following lock levels, listed from lowest to highest restrictions, are available:

Lock Level	Abbreviation
Therapy Lock	T
Rate Change Lock	R
Container Lock	C
Full Lock	F


1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE

1	THERAPY LOCK 
2	RATE CHANGE
3	CONTAINER 
4	FULL LOCK

The GemStar Pump allows two lock levels at the same time. Therapy Lock can be used with either Full Lock, Container Lock, or Rate Change Lock.

When two lock levels are applied, if a function is available under one lock level but not the other, the more stringent restrictions apply. The table on the following page lists the functions available for each lock level.

When the pump is locked under a particular level, a small lock symbol () appears next to the applied lock level as shown above. When the pump is in RUN or STOP mode, the lock symbol appears in the lower left corner of the display and the lock level abbreviation appears in the lower right corner.

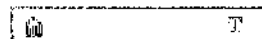

Clinician: refer to Chapter 11 Clinician Instructions on page 121 for detailed instructions.

10. Options Menu

GemStar® Keypad Lock Level Functions						
U = Unlock						C = Contain Lock
T = Therapy Lock						P = Full Lock
R = Rate Change Lock						* = Function is Available
Function	U	T	R	C	P	Comments
Lock/Unlock Keypad	*	*	*	*	*	Rate Change Lock is available only if the current program is Pain Management.
Power Pump On or Off	*	*	*	*	*	
Start or Stop Infusion	*	*	*	*	*	
Access Display, Help, Print, and Program Review Functions	*	*	*	*	*	
Deliver a Bolus	*	*	*	*	*	If programmed.
Silence an Alarm	*	*	*	*	*	
Select a New Container	*	*	*	*	*	
Purge	*	*	*	*	*	
Change a Program	*	*	*	*	*	Under Rate Change Lock, only the delivery rate or bolus dose volume can be changed and must be within the preset values.
Review/Recall Speed Protocol	*	*	*	*	*	Under Therapy Lock, the recalled protocol must be the same type that was programmed when the pump was locked.
Assign or Delete Speed Protocol	*	*	*	*	*	
Set Air Sensitivity, Distal Occlusion Level, Sound Level, and Clock	*	*	*	*	*	
Clear Shift Totals	*	*	*	*	*	
Program and Deliver a Loading Dose (Outside of Programming Mode)	*	*	*	*	*	
Enter a New Program	*	*	*	*	*	Under Therapy Lock, the new program must be of the same type that was programmed when the pump was locked.
Peggybacking	*	*	*	*	*	
VTB Titration	*	*	*	*	*	
Rate Titration	*	*	*	*	*	
Operation Test	*	*	*	*	*	
Clear History	*	*	*	*	*	
Select Another Therapy	*	*	*	*	*	

10. Options Menu

Change Pump Settings

- 1 Press [OPTIONS] and select PUMP SETTINGS.
- 2 Press the number key for the desired pump setting.
- 3 Follow the display messages. Refer to the instructions of the desired pump setting for more information.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE

Select this item from the Pump Settings Menu	To do this function	Comments
AIR SENSOR	Change the air sensor sensitivity.	The pump alarms when it detects an air bubble of approximately 0.5 mL when ON; or at approximately 2 mL when 2 mL is selected. The air alarm does not sound if OFF is selected. The air alarm can also be set during programming if the current therapy allows it.
OCCL SENSOR	Change the distal occlusion sensor sensitivity.	Default Setting: High The pump alarms when distal pressure is greater than the selected occlusion sensor limit.
SOUND LEVEL	Adjust the audible alarm level.	Default Setting: Maximum

Note: The current setting flashes on the display.

Air Sensitivity

WARNING: To reduce the risk of infusing air, use an air-eliminating filter when the air sensitivity is set to OFF.

- 1 Press [OPTIONS] and select AIR SENSITIVITY. The current air sensitivity flashes on the display.

- 2 Select the desired Air Sensitivity level.

Note: If OFF is selected, the pump displays:

AIR ELIMINATING
FILTER IS RECOMMENDED.
ENTER TO CONFIRM

You must press [YES/ENTER] to confirm.

- 3 After selecting the desired Air Sensitivity, the pump displays: CHANGE AIR ALARM FROM XXXX TO XXXX:

AIR SENSITIVITY
1 ON
2 2 mL
3 OFF

CHANGE AIR ALARM
FROM: XXXX
TO: XXXX
YES OR NO

- Press [NO] to keep the current air sensitivity. The current air alarm status displays and the pump exits OPTIONS.
- Press [YES/ENTER] to accept the new setting. The new air alarm status displays and the pump exits OPTIONS.

Occlusion Sensor

- 1 Press [OPTIONS] and select PUMP SETTINGS.

- 2 Select DISTAL OCCLUSION from the Pump Settings Menu.

- 3 Select the desired setting.

- 4 The pump displays MAKE XXXXX DEFAULT?

DISTAL OCCLUSION
1 HIGH
2 MEDIUM
3 LOW

- Press [NO] to temporarily change the distal occlusion setting. The new setting will be in effect until a new program is entered or the pump has been off for more than 24 hours.
- Press [YES/ENTER] to permanently change the distal occlusion setting. The new setting will remain in effect until the setting is changed again through the Options Menu.

10. Options Menu

Sound Level

- 1 Press [OPTIONS] and select PUMP SETTINGS.
- 2 Select SOUND LEVEL from the Pump Settings Menu.
- 3 When the pump displays the Sound Level function, it sounds intermittent beeps at the current sound level. The blocks on the top line of the display indicate the current sound level.
 - Press  to increase the sound level.
 - Press  to decrease the sound level.
 - Press [CHANGE] to restore the sound level to the default setting.
- 4 Press [YES/ENTER] when done.

- [REVERSE] +
USE   TO ADJUST
THE SOUND LEVEL
ENTER WHEN DONE

Set Clock

- 1 Press [OPTIONS] and select SET CLOCK.
- 2 Select 12-Hour or 24-Hour format (24-Hour format is shown in example).
- 3 Use the arrow keys to set the weekday and month.
- 4 Use the number keys to enter the current day, year, and time.

Note: When entering the time in 12-Hour format, press  or  to toggle between AM and PM.

- 5 Press [YES/ENTER] to confirm.

Tip: To view the clock from the RUN or STOP mode, press [YES/ENTER], then immediately press and HOLD [4]. The clock displays until [4] is released.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE

SET CLOCK TYPE
1 SELECT 12 HOUR
2 SELECT 24 HOUR

TIME IS 10:35
MON, JUL 26, 94
ENTER TO CONFIRM

Print

WARNING: Ensure a cassette is installed in the pump before connecting to a printer or computer.

If the pump is connected to a patient: **BEFORE** connecting the pump to a computer or printer, press [STOP], close the slide clamp, and eject the cassette. Leave the pump set connected to the patient. Install another cassette in the pump before connecting to a printer or computer.

Note: Ensure the printer is connected to the pump BEFORE selecting a print function. Refer to Connecting to the Seiko Printer on page 115 for more information.

- 1 Press [OPTIONS] and select PRINT.
- 2 Press the number key for the desired print function.
- 3 When printing is complete, the pump exits OPTIONS.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE

Select this item from the Print Menu	To do this function
PRINT HISTORY	Print the current program, shift and volume totals, and history event log.
PRINT PROGRAM	Print the current program.
PRINT PROTOCOL	Print a specified speed protocol or a complete list of all speed protocols.
PRINT BOLUS Hx	Print the current program, shift and volume totals, and the bolus history (Hx) for the latest 48 hours.
STOP PRINT	Cancel printing. This function is only available if the pump is printing.

10. Options Menu

Speed Protocol

- 1 Press [OPTIONS] and select SPEED PROTOCOL.
- 2 Press the number key for the desired protocol function.
- 3 Follow the display messages. Refer to the instructions for Speed Protocol functions for more information.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	REPLACE PROTOCOL
8	NEXT DOSE

Select this item from the Speed Protocol Menu	To do this function:
ASSIGN	Store a program in pump memory for later use.
REVIEW/RECALL	Review assigned protocols from the RUN or STOP mode. Recall an assigned protocol from the STOP mode.
DELETE	Delete protocols from pump memory.

Assign

- 1 Press [OPTIONS] and select SPEED PROTOCOL.
- 2 Select ASSIGN from the Speed Protocol menu.
- 3 Select [1]–[9] to assign to the new protocol.

If the selected speed protocol is already assigned, the pump displays REPLACE PROTOCOL ASSIGNMENT?

- Press [YES/ENTER] to replace the protocol. The pump displays the new speed protocol and exits OPTIONS.
- Press [NO] to keep the existing protocol and exit OPTIONS.

10. Options Menu

Review/Recall

- 1 Press [OPTIONS] and select SPEED PROTOCOL.
- 2 Select REVIEW/RECALL from the Speed Protocol menu.
- 3 Select the desired protocol to review/recall.
- 4 Review the program. Press  to begin review.
- 5 Press [YES/ENTER]:
 - If in the STOP mode, press [1] to recall the protocol; or press [2] to review another protocol.
 - If in the RUN mode, the pump exits OPTIONS.

Delete

- 1 Press [OPTIONS] and select SPEED PROTOCOL.
- 2 Select DELETE from the Speed Protocol Menu.
- 3 Select the desired protocol to delete.
- 4 Press [YES/ENTER] to confirm deletion or [NO] to exit OPTIONS without deleting.

Next Dose

This function is only available when the current therapy is Intermittent and the pump is not delivering a dose.

To set the time for the next dose:

- 1 Press [OPTIONS] and select NEXT DOSE.
- 2 Use the number keys to enter the start time for the next dose. Press [YES/ENTER] to accept the change and exit OPTIONS.

1	REVIEW PROGRAM
2	HISTORIES
3	KEYPAD LOCK
4	PUMP SETTINGS
5	SET CLOCK
6	PRINT
7	SPEED PROTOCOL
8	NEXT DOSE



Tip: When entering the time in 12-hour format, press  or  to toggle between AM and PM.

10. Options Menu

Connecting to the Data Port

Use the Data Port to connect the GemStar Pump to a Seiko® DPU 414 or compatible serial printer to print the Operation Test results or the files available from the Print Menu. Using communication software included with most computer operating systems (or available from local sources), these files can also be transferred to a computer.

WARNING: Ensure a cassette is installed in the pump before connecting to a printer or computer.

If the pump is connected to a patient: **BEFORE** connecting the pump to a computer or printer, press [STOP], close the slide clamp, and eject the cassette. Leave the pump set connected to the patient. Install another cassette in the pump before connecting to a printer or computer.

Equipment Needed

	Connect to Printer	Connect to Computer
GemStar Pump	*	*
Hospira Serial Printer Cable List 13073-01	*	*
Seiko DPU 414 or Compatible Serial Printer	*	
Serial DB9 Female-to-Female Gender Changer*		*
Null Modem*		*
Computer*		*

* Available from local sources.

10. Options Menu

Configuration Settings

Refer to the printer or computer documentation for information on configuration settings. Configure the printer or computer as follows to use the Data Port:

Baud Rate	=	2400
Data Bits	=	8
Stop Bits	=	1
Parity	=	None
Carriage Return (+Line Feed)	=	CR+LF
Flow Control	=	None

 *Note: When using the Seiko DPU 414 printer, the switch settings **MUST** be configured to work with the GemStar Pump.*

To configure the switch settings on the Seiko DPU 414:

- 1** Press and hold the ONLINE button while sliding the POWER switch to ON.
- 2** When the list of the current settings start printing, release the ONLINE button.
- 3** Set the switches according to the following table.
Press the ONLINE button to continue setting switches.
 - Press ONLINE to set switch ON.
 - Press FEED to set switch OFF.
- 4** Continue setting switches until "Dip SW setting complete!!" is printed.

10. Options Menu

Switch Settings for Sanko DPU-4114 Printer					
Switch	Function		Setting	Button to Press	
SW1-1	Input Method:	Serial	Off		FEED
SW1-2	Printing Speed:	High	On	ONLINE	
SW1-3	Auto Loading:	On	On	ONLINE	
SW1-4	CR function:	CR and LF	On	ONLINE	
SW1-5	DIP SW setting command:	Enabled	On	ONLINE	
SW1-6	Printing density:	100%	Off		FEED
SW1-7	Printing density:	100%	On	ONLINE	
SW1-8	Printing density:	100%	On	ONLINE	
SW2-1	Print Mode:	40 columns	On	ONLINE	
SW2-2	User defined char:	On	On	ONLINE	
SW2-3	Ordinary Characters:	On	On	ONLINE	
SW2-4	Zero Font:	No slash	On	ONLINE	
SW2-5	Character Set:	American	On	ONLINE	
SW2-6	Character Set:	American	On	ONLINE	
SW2-7	Character Set:	American	On	ONLINE	
SW2-8	Character Set:	American	Off		FEED
SW3-1	Data bit length:	8 bits	On	ONLINE	
SW3-2	Parity:	No parity	On	ONLINE	
SW3-3	Parity condition:	Odd	On	ONLINE	
SW3-4	Flow control:	H/S busy	On	ONLINE	
SW3-5	Baud rate:	2400 bps	On	ONLINE	
SW3-6	Baud rate:	2400 bps	Off		FEED
SW3-7	Baud rate:	2400 bps	Off		FEED
SW3-8	Baud rate:	2400 bps	On	ONLINE	

Connecting to the Seiko Printer


WARNING: Ensure a cassette is installed in the pump before connecting to a printer or computer.

If the pump is connected to a patient: **BEFORE** connecting the pump to a computer or printer, press [STOP], close the slide clamp, and eject the cassette. Leave the pump set connected to the patient. Install another cassette in the pump before connecting to a printer or computer.

CAUTION: Do not connect the GemStar AC Mains Adaptor to a printer.

Note: For maximum battery life, operate the pump on AC mains when connected to a printer or computer.

- 1 Install a cassette into the pump.
- 2 Connect the GemStar serial cable to the printer and to the data port on the bottom of the pump.
- 3 Ensure the printer is ONLINE and ready to receive data. Refer to printer documentation for more information.
- 4 Ensure paper is loaded in the printer.
- 5 Select the desired print function from the pump menu.

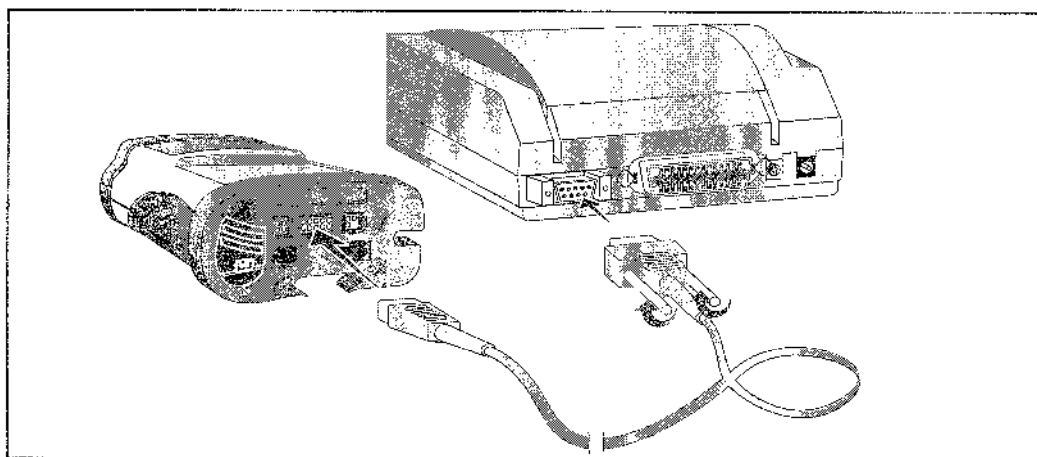

10. Options Menu

Connecting to a Computer


WARNING: Ensure a cassette is installed in the pump before connecting to a printer or computer.

If the pump is connected to a patient: **BEFORE** connecting the pump to a computer or printer, press [STOP], close the slide clamp, and eject the cassette. Leave the pump set connected to the patient. Install another cassette in the pump before connecting to a computer or printer.

Note: For maximum battery life, operate the pump on AC mains when connected to a printer or computer.

- 1 Install a cassette into the pump.
- 2 Connect the GemStar serial cable to the null modem.
- 3 Connect the null modem to the serial gender changer.
- 4 Connect the serial gender changer to the serial port of the computer.
- 5 Connect the GemStar serial cable to the data port on the bottom of the pump.
- 6 Configure the communication software on the computer. Refer to *Configuration Settings* on page 113.

Tip: After the data transfer is complete, use any text editor to view or print the file.

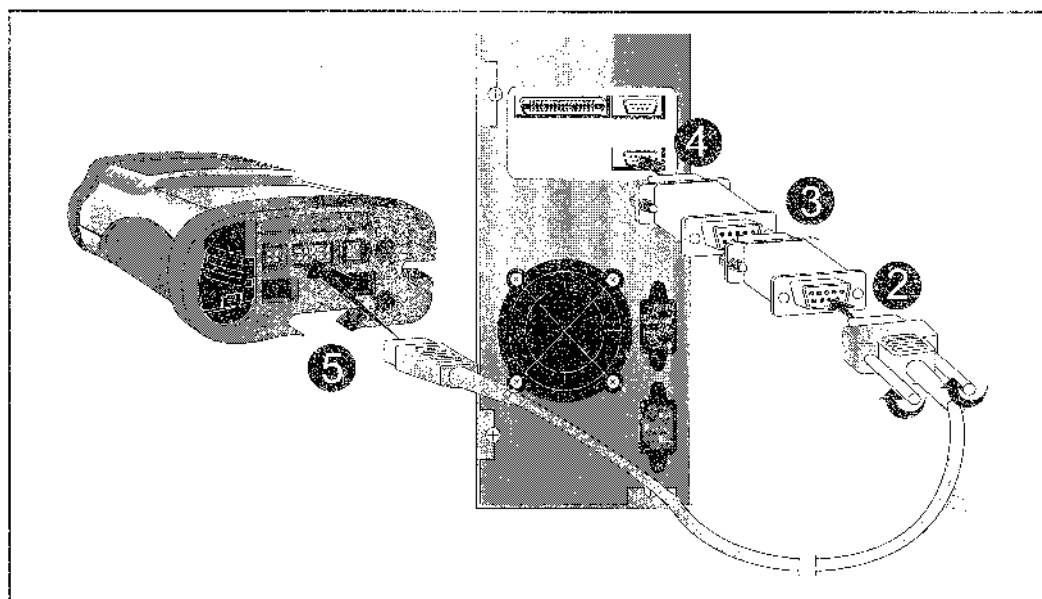

Transferring Data Using a Computer

The GemStar Pump can be connected to a Windows[®]-compatible computer. The following instructions apply to computers using Microsoft[®] Windows[®] 95, Windows[®] 98, or Windows NT[®] 4.0 or higher.

10. Options Menu

Creating a New Connection



Note: The following procedure is only needed to create a “New Connection” the first time the pump is connected to the computer.

To create a new connection:

- 1** Follow steps 1–4 in *Using an Existing Connection* on the following page.
- 2** If a message appears asking if a modem is necessary for a connection, click NO.
- 3** Enter a name for the new connection (e.g., GemStar Print Files) and choose an icon.
- 4** Click OK.
- 5** From the **Connect** drop-down menu, select the appropriate COM port.
- 6** Click OK.
- 7** Verify that the following settings are correct:
 - Baud Rate = 2400
 - Data Bits = 8
 - Stop Bits = 1
 - Parity = None
 - Flow Control = None
- 8** Click OK.
- 9** Select the **Properties** icon.
- 10** Select the **Settings** tab.
- 11** Select the **ASCII** setup button.
- 12** Ensure ‘Append line feeds to incoming line ends’ is checked.
- 13** Click **OK** in the ASCII setup window.
- 14** Click **OK** in the Properties window.
- 15** From the menu, select **File > Save**.
- 16** Exit HyperTerminal.

10. Options Menu

Using an Existing Connection

To transfer data from the GemStar Pump to a computer:

- 1** From the Windows[®] desktop, select the Start Menu.
- 2** Select **Programs**.
- 3** Select **Accessories**.
- 4** Select **HyperTerminal**.
- 5** Double click on the icon with the name entered in step 3, *Creating a New Connection*.
- 6** From the menu, select **Call > Connect**.
- 7** From the menu, select **Transfer > Capture Text**.
- 8** Type the path and file name where the printout should be saved.

Ensure the file extension is “.txt”.

The default path is: c:\windows\accessories\hyperterminal.

Note: If the file name already exists, the new printout will be added to the end of the existing file. For example, use a patient's name or id number as the file name. Then select that file for all of that patient's printouts.

- 9** Click **Start**.
- 10** Select the desired print function from the pump. Refer to the table of printing functions on page 109.
- 11** Wait for the transmission to complete.
- 12** From the menu, select **Transfer > Capture Text > Stop**.
- 13** Exit HyperTerminal.
- 14** Use any text editor to view or print the file.

Notes

Clinician Instructions

The instructions in this section apply to special features available to clinicians. **Remove this section from the manual to prevent unauthorized access to these features.**

The GemStar pump provides quick access to commonly used functions. Instead of accessing these functions by navigating through the menu system, you can enter the key combinations shown in this table:

Function	Quick Access Key Combination
Unlock Keypad	Place the pump in STOP mode and then enter the lock sequence to unlock a Full Container or Rate Lock. You cannot use quick access keys to unlock a Therapy Lock.
Displays Shift Totals	Place the pump in STOP mode and then press [HELP] to display Shift Totals.
Display Clock	Place the pump in STOP mode, press [YES], and then press and hold [1]. The clock displays until you release [1].
Program a Loading Dose	From the Pain Management STOP MODE, press [YES/ENTER], and then press [0] to program a Loading Dose.



11. Clinician Instructions

Clinician Locking Sequence Instructions

Locking the Pump

To lock the pump from the RUN or STOP mode:

- 1 Press [OPTIONS] to display the Options Menu.
- 2 Select **KEYPAD LOCK** from the Options Menu.

Note: If a lock level is already applied to the pump, the Lock Sequence Number must be entered before a second lock level can be applied. Refer to Adding a Second Lock Level on page 123 for more information.

- 3 Select the desired lock level. For THERAPY, CONTAINER and FULL Lock, go to step 5.
- 4 For RATE CHANGE Lock, set the minimum and maximum programmable ranges for the continuous delivery rate and bolus dose in a Pain Management program.

Tip: Press  to toggle between the MIN and MAX values.

- MIN flashes on the display. Enter the minimum value.
- MAX flashes on the display. Enter the maximum value.

- 5 Press [YES/ENTER] when done.
- 6 Repeat steps 4 and 5 for a bolus dose if available in the current program.
- 7 Enter the lock sequence number. Asterisks display as the numbers are entered.
 - For THERAPY Lock, enter [0][0][9][5][1].
 - For UNLOCKED, FULL, CONTAINER, and RATE CHANGE lock levels, enter [1][3][0][0][6].

Note: Default lock sequence numbers are shown.

11. Clinician Instructions

Adding a Second Lock Level

- 1** Press [OPTIONS] to display the Options Menu.
- 2** Select **KEYPAD LOCK** from the Options Menu.
- 3** Enter the Lock Sequence Number.
- 4** For **RATE CHANGE** Lock, set the minimum and maximum programmable ranges for the continuous delivery rate and bolus dose in a Pain Management program.

Tip: Press  to toggle between the MIN and MAX values.

- MIN flashes on the display. Enter the minimum value.
- MAX flashes on the display. Enter the maximum value.

- 5** Press [YES/ENTER] when done.
- 6** Repeat steps 4 and 5 for a bolus dose if available in the current program.

Note: When two lock levels are applied, if a function is available under one lock level but not the other, the more stringent restrictions apply.

- 7** Enter the lock sequence number. Asterisks display as the numbers are entered.
 - For **THERAPY** Lock, enter [0][0][9][3][1].
 - For **UNLOCKED**, **FULL**, **CONTAINER**, and **RATE CHANGE** lock levels, enter [1][3][9][9][0].

Note: Default lock sequence numbers are shown.

11. Clinician Instructions

Unlocking the Pump

- 1 Press [OPTIONS] to display the Options Menu.
- 2 Select KEYPAD LOCK from the Options Menu.
- 3 Select the lock level to be unlocked.
- 4 Enter the lock sequence number for that level.

Note: If the pump is locked under two levels, each level must be unlocked separately. Unlocking one level has no effect on the other level. For example, if both THERAPY and RATE CHANGE lock levels are applied, the clinician must unlock both the RATE CHANGE and THERAPY locks.

- 5 The pump displays the current lock levels. Press [YES/ENTER] to return to the mode the pump was in when [OPTIONS] was pressed.

Program and Deliver a Clinician Activated Loading Dose

This feature is available only when the keypad is UNLOCKED or when THERAPY is the only lock level. The default loading dose rate is 125 mL/hr.

- 1 From the stop mode, press [YES/ENTER] [0].
- 2 The pump displays: PROGRAM A LOADING DOSE?
Press [YES/ENTER].
- 3 Enter the loading dose. Press [YES/ENTER].
- 4 The pump displays: DELIVER THE LOADING DOSE NOW?

To deliver the loading dose immediately, press [YES/ENTER].

Primary infusion DOES NOT begin automatically when the clinician activated loading dose is complete. Press [START] to continue primary infusion.

To deliver the loading dose later, press [NO].

The pump returns to the stop mode. The next time [START] is pressed, the pump again displays: DELIVER LOADING DOSE NOW?

- Press [YES/ENTER] to deliver the loading dose. Primary infusion continues automatically when the loading dose is complete.
- Press [NO] to cancel the loading dose and continue primary infusion.

Shift Totals and History

Shift Totals

- 1** Press [OPTIONS] and select HISTORIES.
- 2** Select SHIFT TOTALS from the Histories Menu.
- 3** The shift totals display. Press  and  to scroll the display. When done, the pump displays CLEAR TOTALS? YES OR NO.
 - Press [YES/ENTER] to clear the shift totals.
 - Press [NO] to exit OPTIONS without clearing the shift totals.

END OF
SHIFT TOTALS
CLEAR TOTALS?
YES OR NO 



Tip: Press [HELP] to display shift totals when the pump is not in PROGRAMMING mode or an alarm condition.

Clear History Event Log and Bolus History

- 1** Press [OPTIONS] and select HISTORIES.
- 2** Select CLEAR HISTORY.
- 3** CLEARING HISTORY displays and the pump exits OPTIONS.



Notes



Chapter 12

Troubleshooting

This section contains information on audible and visual alarms that may occur with the GemStar Pump.

CAUTION: If the pump does not perform as stated in this manual, stop using it **IMMEDIATELY**.

For consultation and technical support contact either an authorized Hospira representative or the Hospira Technical Support Operations hotline (available 24 hours a day).

To return the pump for service, first contact Hospira Customer Support between 7:00 AM and 4:00 PM Pacific time. Have the pump serial number available.

In the USA

Hospira Technical Support Operations
1-800-241-4002

Hospira, Inc.
Technical Support Operations
755 Jarvis Drive
Morgan Hill, CA 95037

Outside the USA:

Contact your local Hospira sales office.



Tip: Press [HELP] for additional information during programming and alarm conditions.

12. Troubleshooting

Alarm Messages

Alarm Display Message	Is the pump infusing?	Alarm	Cause/Condition	Remedy
Y = Yes		A = Audible alarm		
N = No		V = Visual alarm		
AIR IN LINE • Pump automatically enters STOP mode.	N	A/V	The specified volume of air has been detected in the line distal to the cassette while infusing.	1 Press [SILENCE] to quiet the alarm for 1 minute. 2 Press [STOP] to clear the alarm. 3 Disconnect the patient from the pump set. 4 Press and hold [PURGE]. <i>Note: Purge rate is 250 mL/hr. The pump purges for up to 2 minutes each time you press and hold [PURGE]. Purge volume may be limited in a Pain Management program.</i> 5 Reconnect the pump set to the patient-access device. 6 Press [START] to begin infusion.
CHANGE BATTERIES	N	A/V	The pump has entered the STOP mode due to battery voltage below acceptable levels.	1 Replace both disposable batteries with fresh AA batteries and press  to restart the pump; OR connect the pump to the AC mains adaptor and the pump will automatically restart. Replace both batteries as soon as possible. 2 After power-on, resume the current program. Refer to <i>Accessing the Programming Menu</i> in Chapter 2 Basic Operation for more information.

12. Troubleshooting

Alarm Display Message	is the pump infusing?	Alarm	Cause/Condition	Remedy
Y = Yes N = No			A = Audible alarm V = Visual alarm	
CHECK CASSETTE-X	N	A/V	The cassette is improperly installed or not installed. X indicates where to look for the cause of the alarm: A = Air Sensor D = Distal P = Proximal	1 Press [SILENCE] to quiet the alarm for 1 minute. 2 Press [STOP] to place the pump in the STOP mode. 3 Reinstall the cassette. CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. <i>Note: Ensure the tubing is behind the tubing guide (at the distal end). Refer to Loading the Cassette on page 32.</i> 4 Press [START] to begin infusion. <i>Note: The pump cannot be powered off during a CHECK CASSETTE alarm (except during PROGRAMMING mode).</i>
CALL 1.800.XXX.XXXX CODE: NN/MMM/TTT	N	A/V	The pump has detected a condition that may require service.	1 Write down the telephone number and code. 2 Power-off the pump. If the pump does not power-off, disconnect external power and remove the batteries. 3 Call the telephone number shown on the display message.

12. Troubleshooting

Alarm Display Message	Is the pump infusing?	Alarm	Cause/Condition	Remedy
Y = Yes N = No			A = Audible alarm V = Visual alarm	
DIST. OCCLUSION	N	A/V	The pump has detected sustained distal pressure over the user-selected limit. To avoid a bolus infusion after an occlusion, complete the following steps: 1 Press [STOP]. 2 Close the distal clamp. 3 Press the cassette release button to remove the cassette from the pump. 4 Open the flow stop. 5 Wait 10 seconds. 6 Close the flow stop. 7 Load the cassette into the pump. 8 Open the distal clamp. 9 Review the program. 10 Press [START] to start the pump.	1 Press [SILENCE] to quiet the alarm for 1 minute. 2 Press [STOP] to place the pump in the stop mode. 3 Check for sources of distal occlusion: ▶ Kinked tubing below the pump ▶ Clamped patient-access device ▶ Clogged I.V. filter ▶ Other obstructions 4 If the problem continues, remove and reinstall the cassette. CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. <i>Note: Ensure the tubing is behind the tubing guide (at the distal end). Refer to Loading the Cassette on page 32.</i> 5 Press [START] to begin infusion.

12. Troubleshooting

Alarm Display Message	Is the pump infusing?	Alarm	Cause/Condition	Remedy
Y = Yes N = No			A = Audible alarm V = Visual alarm	
LOW BATTERIES	Y	AV	The pump has detected battery voltage below acceptable levels and the batteries are the primary power source. <i>Note: The pump allows the user to clear the LOW BATTERIES alarm. If the condition has not been corrected after 15 minutes, the alarm continues.</i> <i>If the pump is infusing and a LOW BATTERIES alarm occurs, the pump continues for approximately 36 minutes* before it can no longer deliver due to insufficient power.</i> <i>*If an external power source is used, the rate is changed, or the pump is powered-on after a LOW BATTERIES alarm and those same batteries are still in the pump, this interval may be different.</i>	1 Press [SILENCE] to quiet the alarm for 15 minutes. 2 Press [STOP] to place the pump in the STOP mode. 3 Connect the pump to AC mains adaptor, rechargeable battery pack, OR replace both disposable batteries with fresh AA batteries. 4 Press [START] to begin infusion.

12. Troubleshooting

Alarm/Display Message	Is the Pump Infusing?	Alarm	Cause/Condition	Remedy
Y = Yes			A = Audible alarm	
N = No			V = Visual alarm	
POWER LOSS	N	A/V	The pump has detected a power failure due to one of the following: - Batteries have been removed while the pump was on battery source. - External power has failed or has been disconnected and no batteries are installed. - Pump has been powered on with insufficient power source.	1 Connect the pump to AC mains adaptor, rechargeable battery pack, OR replace both disposable batteries with fresh AA batteries. 2 Power-on the pump. Review the program if required. 3 Press [START] to begin infusion.

12. Troubleshooting

Alarm Display Message	Is the pump infusing?	Alarm	Cause/Condition	Remedy
Y = Yes N = No			A = Audible alarm V = Visual alarm	
PROX. OCCLUSION	N	A/V	The pump has detected sustained negative pressure on its proximal side.	1 Press [SILENCE] to quiet the alarm for 1 minute. 2 Press [STOP] to place the pump in the stop mode. WARNING: When a Proximal Occlusion alarm activates, DO NOT power-off the pump. 3 Check for sources of proximal occlusion: ▶ Kinked tubing above pump ▶ Blockage in the container and in the spike ▶ Closed clamps above pump ▶ Empty container ▶ High torque syringe ▶ Other obstructions 4 Remove obstructions. 5 Press [START] to begin infusion. 6 If the problem continues, power-off the pump, install a new primed set and begin infusion. CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. <i>Note: Ensure the tubing is behind the tubing guides (at the distal end). Refer to Loading the Cassette on page 32.</i>

12. Troubleshooting

Alarm/Display Message	Is the Pump Running?	Alarm	Cause/Condition	Remedy
Y = Yes N = No			A = Audible alarm V = Visual alarm	
USING BATTERIES	Y	V	This is a normal response when: - The pump has been powered on with only disposable batteries as a power source; or - The pump has switched to internal power (disposable batteries) due to a loss of external power or a reduction in external power below acceptable levels.	When powering on the pump: - Press [YES/ENTER] to confirm disposable batteries are the only power source. When switching from AC mains power: - Press [SILENCE] to clear the alarm.
USING EXT BATT	Y	V	This is a normal response when: - The pump has been powered on with external batteries as a power source; or - The pump has switched to external battery power due to a loss of AC mains power; or - The pump has switched to external battery power while running on internal batteries.	When powering on the pump: - Press [YES/ENTER] to confirm external battery power source. When switching from AC mains power: - Press [SILENCE] to clear the alarm.

12. Troubleshooting

Alert Messages

Alert Display Message	Is the pump infusing?	Alarm	Cause/Condition	Remedy
Y = Yes N = No			A = Audible alarm V = Visual alarm	
ALMOST EMPTY	Y	A/V	Delivery will complete in 30 minutes or less. For bolus only program, delivery will complete during the next bolus delivery.	Press [SILENCE] to quiet the alarm for 10 minutes. <i>Note: This alert is available only with Pain Management protocols.</i>
CHECK PRINTER	Y	A/V	- The printer has not responded to commands from the pump. - The printer is not properly connected to the pump.	- Press [SILENCE] to quiet the alarm for 2 minutes. - Ensure the printer is properly connected to the pump by checking the following: - Printer is powered on and in ONLINE mode. - Cable is plugged in to both the pump and the printer. - Cable part number and printer settings are correct. - If the problem continues, select STOP PRINT from the Print Menu to cancel all print jobs.

12. Troubleshooting

Alert Display Message	Is the pump beeping?	Alarm	Cause/Condition	Remedy
Y = Yes N = No			A = Audible alarm V = Visual alarm	
EMPTY CONTAINER	N	A/V	The pump has delivered the programmed container size or the pump has delivered the programmed VTBI AND Auto KVO is delivering (if available in current program).	1 Press [SILENCE] to quiet the alarm for 2 minutes. 2 Press [STOP] to place the pump in the STOP mode and clear the alarm. 3 Use the NEW CONTAINER function to repeat the current program. Refer to <i>Starting a New Container</i> in Chapter 2 Basic Operation.
END OF INFUSION • The pump delivers KVO if selected.	N	A/V	The pump has delivered the entire programmed protocol.	If no KVO is selected: 1 Press [SILENCE] to quiet the alarm for 2 minutes. 2 Press [STOP] to place the pump in the STOP mode and clear the alarm. If KVO is selected: 1 The pump sounds 6 beeps approximately every 15 minutes. 2 Press [STOP] to place the pump in STOP mode and clear the alarm.

12. Troubleshooting

Alert Display Message	Is the pump infusing?	Alarm	Cause/Condition	Remedy
Y = Yes N = No		A = Audible alarm V = Visual alarm		
FLASHING DISPLAY • Line four of the display flashes.	N	A/V	A Non-programming Callback has occurred because the pump has been waiting for more than three minutes for a keypress on a non-programming input or a confirmation display.	Follow the instructions on the display to clear the alarm.
NO INTERNAL BATT NO POWER LOSS ALARM AVAILABLE ENTER TO CONFIRM	N	A/V	No disposable batteries or dead disposable batteries; AND external power.	1 Press [YES/ENTER] to confirm. 2 Power-off the pump and install two fresh, disposable AA batteries as a backup power source.
PROGRAMMING INCOMPLETE	N	A/V	Programming has not been completed and no keys have been pressed within three minutes.	1 Press any key to clear the alarm. 2 Complete the program entries.
START • START flashes on the display. Other parts of the display are not affected.	N	A/V	A Stop Mode Callback has occurred because the programmed pump has remained in the STOP mode for more than three minutes.	1 Press [SILENCE] to quiet the alarm for 3 minutes. 2 Follow the instructions on the display to clear the alarm.

12. Troubleshooting

Other Display Messages

Other Display Message	Cause/Condition	Remedy
KEYPAD LOCKED	User has attempted to perform a function that has been locked out.	Unlock the keypad to access the desired function.
NOT ALLOWED DURING INFUSION PRESS STOP TO HALT DELIVERY	User has attempted to perform a function that is not available in RUN mode.	- Press [STOP] to access the desired function. - Press [YES/ENTER] to return to the RUN mode display.
ROUNDING The pump sounds four quick beeps and displays the rounded number.	The pump has rounded a number that was entered during programming.	Press [YES/ENTER] to accept the rounded value and continue to the next step.
DOSE IN PROGRESS NEW CONTAINER NOT ALLOWED	User has attempted to access <i>NEW CONTAINER</i> function when it is not allowed.	1 Press [STOP]. 2 Increase the container size to complete the dose. 3 Set <i>NEW CONTAINER</i> through the Change Menu.
INFUSED XXXX.X mL NO START DELAY PRESS START TO INFUSE	The pump has not delivered any doses and there is no start delay.	Press [START] to infuse dose.
CANNOT CHANGE CLOCK WHILE TIME BASED DELIVERY IS IN PROGRESS	User has attempted to change the clock during a time-based delivery.	1 Press [STOP]. 2 Clear the program through the Change Menu. 3 Change the clock through the Options Menu. 4 Reprogram the pump.
PROGRAM EXCEEDS MAX CONTAINER SIZE. CHANGE PROGRAM VALUES.	The V(BI) + FVO volume resulted in a container size greater than the allowed value.	1 Press [HELP] to display the allowed values. 2 Change the program values.
THIS OPTION NOT AVAILABLE WHILE PROGRAMMING THE PUMP	User has attempted to access a function that is not available from the first programming display.	1 Press [YES/ENTER] to return to the first programming display. 2 Complete programming. 3 Access desired function.

Chapter 13

Maintenance

Pump Storage

Store the GemStar Pump in a cool, dry place. Remove the disposable batteries and the battery pack before storing the pump.

Program and event history are protected in the software memory for at least one year when power is removed from the pump.

Cleaning and Disinfecting

Establish a routine schedule for cleaning the pump to keep the case exterior, cassette pocket, and tubing channel free of contamination.

CAUTION: To avoid mechanical or electronic damage, do not immerse pump in any fluids or cleaning solutions.

Some cleaning and sanitizing compounds may slowly degrade components made from some plastic materials. Do not use compounds containing combinations of isopropyl alcohol and dimethyl benzyl ammonium chloride.

Do not sterilize by heat, steam, ETO, or radiation. Apply disinfectants to the outside surface of the pump only. Do not use abrasive cleaners or materials on the pump. Using abrasive cleaners or cleaning solutions not recommended by Hospira may result in product damage.

To avoid pump damage, cleaning solutions should be used only as directed in the table on the following page. The disinfecting properties of cleaning solutions vary; consult the manufacturer for specific information.

Never use sharp objects (such as pens, pencils, fingernails, paper clips, needles, etc.) to clean the pump.

13. Maintenance

Recommended Cleaning Solutions

The pump is not adversely affected by the recommended cleaning solutions in the table below:

Cleaning Solution	Manufacturer	Preparation
Super Edisonite®	S. W. Edison Chemical Co.	Per manufacturer's recommendation
Vesphene II ssc®	Calgon Vestal Laboratories	Per manufacturer's recommendation
Manu-Klenz®	Calgon Vestal Laboratories	Per manufacturer's recommendation
Formula C™	Diversey Corp.	Per manufacturer's recommendation
Household bleach (5.25% sodium hypochlorite)	Various	Per healthcare facility procedures; DO NOT exceed one part bleach in ten parts water



Note: NOT all cleaning solutions are disinfectants. Check product labeling.

Cleaning the Pump Case


- 1 Power-off the pump.
- 2 Disconnect the pump from AC mains power.
- 3 Clean the exposed surfaces of the pump with a soft, lint-free cloth dampened with the appropriate cleaning solution.
- 4 Wipe the solution from the pump surface with a water-moistened, lint-free cloth.
- 5 Use a dry, lint-free cloth to dry the pump after cleaning.

Cleaning the Cassette Pocket and Tubing Channel

The area containing the ultrasonic sensors, located in the cassette pocket, should be cleaned on a regular basis.

CAUTION: Use care not to damage the silicone seals around the sensor bodies.

- 1** Remove the cassette, if installed.
- 2** Use a cotton swab, moistened with an approved cleaning solution, to clean the sensor faces (A), tubing channel (B), plunger tip (C), and sensor pins (D).
- 3** Dry the sensor faces and tubing channel and ensure the sensor faces are free of detergent film and/or debris.

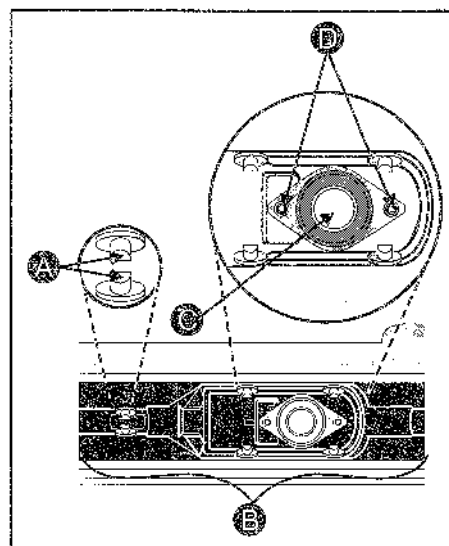

Repair

- ✦ The GemStar Pump has no user-serviceable components, with the exception of disposable batteries.
- ✦ If the display clock stops working, contact Hospira Technical Support Operations to arrange for internal lithium battery replacement.
- ✦ **Homecare Customers:** Call your healthcare professional or homecare company regarding any required service or repairs. **DO NOT** attempt to repair the pump for any reason.
- ✦ The GemStar Pump is covered by a manufacturer's warranty for one year after purchase. During this time, opening the pump case for any reason voids this warranty. Refer to Chapter 18 Warranty/Technical Assistance, page 189 for details.
- ✦ Refer all service to Hospira certified technicians only. A *GemStar Technical Field Service Manual* is available to Hospira certified technicians.

Notes

Operation Test

This test is designed to ensure the GemStar Pump is operating properly. Hospira recommends performing this test a minimum of once every 12 months. Refer to facility requirements to determine additional testing needs.

The test can be performed when the pump is in the STOP mode and requires approximately 15 minutes to complete. The pump provides step-by-step instructions through each section of the test.

Note: This test is available only when the pump is UNLOCKED.

WARNING: Disconnect the pump from the patient BEFORE performing the Operation Test.

Equipment Required

- GemStar Pump
- Two GemStar Microbore Pump Sets without a filter
- Reservoir with at least 50 mL of water
- 20+ mL or larger graduated cylinder (graduated in 0.1 mL)
- Two fresh, disposable AA batteries
- Appropriate GemStar AC Mains Adaptor or other Hospira approved external power source

14. Operation Test

Test Setup



Note: The pump clears any program and shift totals stored in pump memory before starting the Operation Test.

- 1 Connect the external power source and install two fresh, disposable AA batteries.
- 2 Power-on the pump.
- 3 Press [OPTIONS].
- 4 Select HISTORIES from the Options Menu.
- 5 Select OPERATION TEST from the Histories Menu.
- 6 Follow the instructions on the display.

Performing the Operation Test




Tip: The checklist at the end of this section may be copied to record test results.

Case Inspection

- 1 Visually inspect the case and keypad for damage. Check the case for missing screws, or labels.
- 2 Enter the results:
 - If the pump is free of damage and all the screws and labels are in place, press [YES/ENTER].
 - If damage is found on the case or keypad, or any screws or labels are missing, press [NO].

Cassette Pocket and Latch Inspection

- 1** Visually inspect the cassette pocket and the sensors (A). Verify that they are free of dirt and residue.
- 2** Verify that the four cassette latches (B) are present. Verify that they open and close when the cassette release button is pressed.
- 3** Verify that the plunger seal (C) and sensor seals (D) are present and intact.
- 4** Verify that the plunger (E) is not scraped or damaged.
- 5** Verify that the sensor pins (F) are not damaged.
- 6** Install the cassette. Verify the cassette fits smoothly and securely.
- 7** Enter the results:
 - If the cassette pocket area and latches are free of damage, and ALL of the above steps are verified, press [YES/ENTER].
 - If damage is found, or ANY of the above steps cannot be verified, press [NO].

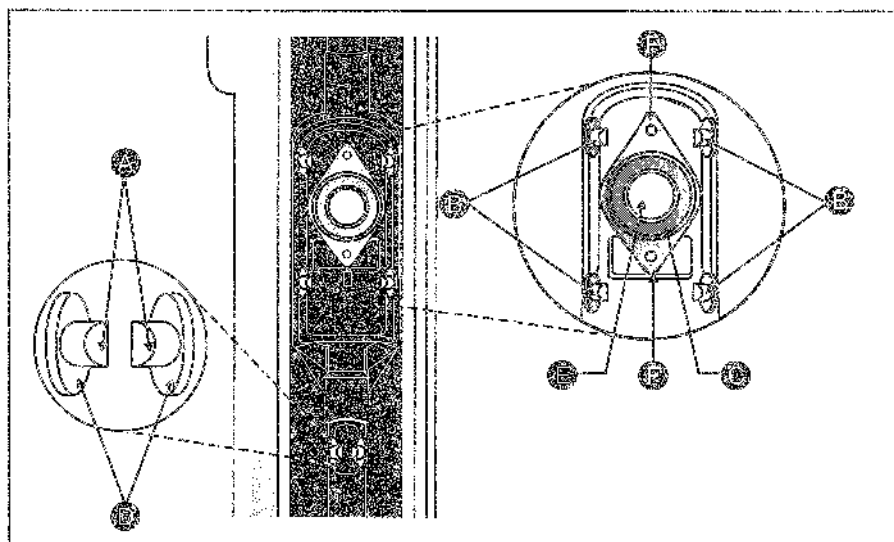

14. Operation Test

Power Check



Tip: Opening the battery compartment door will disconnect the batteries enough for the pump to recognize external power. Then the door can be tightened and testing can continue. This can save time when performing multiple tests on a group of pumps.

- 1 The display indicates that an external power source and batteries are connected. Press .
- 2 Remove battery door and press . The display indicates that only external power is connected. Press .
- 3 Return the battery door to the closed position.
- 4 Disconnect the external power source and press .
- 5 The display indicates that only batteries are connected and displays the test results. Press .
- 6 If the pump fails, reconnect the external power source and press  to rerun the test, or press [YES/ENTER] to confirm failure.

Keypad Test

- 1 Press each key, including the Bolus Button. Press [YES/ENTER] last.
- 2 Enter the results:
 - If each key makes an audible tone, press [YES/ENTER].
 - If any key does not make an audible tone, press [NO].

Display Test

- 1 Press [YES/ENTER] when prompted. The display fills with solid boxes for approximately 10 seconds. To view the solid boxes for another 10 seconds after the display times-out, press .
- 2 Enter the results:
 - If the display fills with solid boxes, press [YES/ENTER].
 - If the display does not fill with solid boxes, press [NO].

LED Test

- 1 Verify that both the Power (green) and Alarm (red) LEDs are lit. Refer to the *GemStar Pump Layout* in the "Introduction" section for the location of the LEDs.
- 2 Enter the results:
 - If both LEDs are lit, press [YES/ENTER].
 - If either of the LEDs does not light, press [NO].



Tip: Use two microbore pump sets to save time during the next four tests. The first set should be fully primed BEFORE beginning the volume accuracy test. Use a second, DRY set for the air-in-line test.

When the proximal occlusion test is complete, remove the WET set from the pump and insert the DRY set. Press [START] to begin the air-in-line test. The WET set can then be placed in the next unit ready for testing.

Volume Accuracy Test

The pump should deliver between 19.0 to 21.0 mL within approximately 3 minutes to pass this test.

- 1 Load a primed pump set into the pump. Press .
- 2 Place the proximal end in the reservoir of water. Press .
- 3 Place the distal end in a 20+ mL graduated cylinder. Press .
- 4 Press [START].
- 5 Enter the results:
 - If the pump delivers between 19.0 to 21.0 mL, press [YES/ENTER].
 - If the pump does not deliver between 19.0 to 21.0 mL, press [NO].

14. Operation Test

Distal Occlusion Test

The occlusion should occur within approximately 30 seconds to pass this test.

- 1 Clamp the pump set 10 to 12 inches (25.4 to 30.5 cm) below the cassette.
- 2 Press [START].
- 3 The test result displays. Release the clamp and press .
 - If the test is passed, the pump advances to the next test.
 - If the test is failed, clean the sensors and press  to rerun the test or press [YES/ENTER] to confirm failure.

 *Note: Clean the sensor pins before repeating the distal occlusion test. Refer to Chapter 13 Maintenance for instructions.*

Proximal Occlusion Test

The occlusion should occur within approximately 30 seconds to pass this test.

- 1 Clamp the pump set above the cassette.
- 2 Press [START].
- 3 The test result displays. Press .
 - If the test is passed, the pump advances to the next test.
 - If the test is failed, clean the sensors and press  to rerun the test or press [YES/ENTER] to confirm failure.

 *Note: Clean the sensor pins before repeating the proximal occlusion test. Refer to Chapter 13 Maintenance for instructions.*

Air-in-line Test

The air-in-line alarm should occur within approximately 1 minute to pass this test.

- 1** Release the clamp.
- 2** Remove the proximal end of the pump set from the reservoir. Press .
- 3** Press [START].
- 4** The test result displays. Press .
 - If the test is passed, the pump advances to the next test.
 - If the test is failed, clean the sensors and press  to rerun the test or press [YES/ENTER] to confirm failure.



Note: Clean the sensor faces before repeating the air-in-line-test. Refer to Chapter 13 Maintenance for instructions.

Printing Test Results

When the test is complete, the pump displays an option to print the test results. If desired, follow the steps below to print the test results. An example of the printed test results is shown.



Note: For maximum battery life, operate the pump on AC mains power when connected to a printer or computer.

- 1** Connect the serial printer cable to the pump. Refer to *Connecting to the Seiko Printer* in Chapter 10 Options Menu, on page 115.
- 2** Press [YES/ENTER] to print the test results.
- 3** The pump displays the Programming Menu when the results have finished transmitting or if [NO] is pressed.

14. Operation Test

```
*****
*          HOSPIRA GemStar          *
*          OPERATION TEST            *
*          RESULTS                   *
*****

PUMP SERIAL NUMBER: XXXXXXXXXX

          10:45AM JUL  4, 04
          VERSION X.X.X

1. Case Inspection          PASS
2. Cassette Pocket & Latch PASS
3. Power Check              PASS
4. Keypad Test              PASS
5. Display Test             FAIL
6. LED Test                 PASS
7. Volume Accuracy Test     PASS
8. Distal Occlusion Test    PASS
9. Proximal Occlusion Test  PASS
10. Air-in-Line Test        PASS

*****
*          VERIFIED BY:             *
*****
*          END OF                   *
*          OPERATION TEST           *
*          RESULTS                   *
*****
1234567890123456789012345678901234567890
```

14. Operation Test

GemStar® Operation Test Checklist

	Pass	Fail
Case Inspection	☐	☐
Cassette Pocket and Latch Inspection	☐	☐
Power Check	☐	☐
Keypad Check	☐	☐
Display Test	☐	☐
LED Test	☐	☐
Volume Accuracy Test	☐	☐
Distal Occlusion Test	☐	☐
Proximal Occlusion Test	☐	☐
Air-in-line Test	☐	☐
Printed Test Results	☐	☐
Serial Number:		
Verified By: _____ Date: _____		

Notes

Optional System Components

All GemStar Optional System Components are composed of latex-free components only.

Product Description	List #	Description
GemStar AC Mains Adaptor (110 Volts - USA)	13026-01	Wall-mount adaptor plugs into a 110 VAC mains supply to power the GemStar Pump. 6 ft (1.8 m) cord.
GemStar AC Mains Adaptor (100 – 240 Volts - Universal)	13072-01	Tabletop adaptor plugs into a 100 – 240 VAC mains supply to power the GemStar Pump.
GemStar Bolus Cord	13027-01	Allows bolus requests to be made up to 6 feet (1.8 m) from the pump.
GemStar Pole Clamp	13076-01	Attaches the GemStar Pump to an I.V. pole. The clamp is mounted to the track on the back of the pump. Press the release button to remove the pump.
GemStar Battery Pack	13073-01	The rechargeable nickel metal hydride batteries power the GemStar Pump during periods of transport or when use of AC power is not desirable. The battery pack is recharged using a GemStar AC Mains Adaptor.
GemStar Universal Lockbox	12499-01	Locks the GemStar Pump and up to a 250 mL flexible I.V. container or 30 mL PCA vial in a plastic enclosure. The lockbox can be attached to an I.V. pole using the GemStar Pole Clamp (available separately).
GemStar 250 mL Lockbox	13024-01	Locks the GemStar Pump and up to a 250 mL flexible I.V. container in a plastic enclosure. The lockbox can be attached to an I.V. pole using the GemStar Pole Clamp (available separately).

15. Optional System Components

Product Description	List #	Description
GemStar 30 mL PCA Vial Lockbox	13077-01	Locks the GemStar Pump and a 30 mL PCA vial in a plastic enclosure. The lockbox can be attached to an I.V. pole using the GemStar Pole Clamp (available separately).
Hospira Carrying Case - 250 mL	13534-01	Soft-sided carrying case can be worn around the waist to hold the GemStar Pump and up to a 250 mL flexible container.
Hospira Carrying Case - 500 mL	13079-01	Soft-sided carrying case can be worn around the waist or over the shoulder to hold the GemStar Pump and up to a 500 mL flexible container.
Hospira Carrying Case - 1000 mL	13980-01	Soft-sided carrying case can be worn around the waist or over the shoulder to hold the GemStar Pump and up to a 1000 mL flexible container.
Hospira Carrying Case - 3 Liter	13593-01	Soft-sided carrying case can be worn around the waist or as a backpack to hold a 3 liter or smaller flexible container and the GemStar Pump with batteries or battery pack. Includes List 13590-01.
Hospira 3 Liter Carrier	13590-01	Interior support for the Hospira Carrying Case - 3 Liter.
GemStar Docking Station (100 - 240 Volts)	13075-01	Plugs into a 100 - 240 VAC Volt mains supply to power the GemStar Pump while mounted on an I.V. pole. During periods of transport, an internal rechargeable battery pack powers the pump for up to 16 hours at a rate of 125 mL/hr or less. When the AC mains supply is attached, the internal batteries are recharged in 4 to 6 hours. 12 ft. (3.8 m) cord.
GemStar Serial Cable	13078-01	Connects the GemStar Pump to a compatible serial printer or computer to download the event history, program, or speed protocol information.

15. Optional System Components


Product Description	List #	Description
GemStar Connect Clinician Kit	13048-01	Allows remote access from a clinician's computer to a GemStar Pump. The Kit includes a modem, phone cable, serial cable, and GemStar Connect software.
GemStar Connect Patient Kit	13049-01	Use with the GemStar Connect Clinician Kit to allow remote access to a GemStar Pump. The Patient Kit includes a modem, phone cable, and remote pump cable.
GemStar Triple Pole Clamp	13022-01	Not for use with a lockbox.
GemStar Technical Field Service Manual	13000-38	Provides advance instructions for Hospira Certified Technicians.
GemStar Configuration Guide	13000-97	Provides advance instructions for biomed technicians.
GemStar System Operating Manual	13981-01	Provides operating instructions for patients and caregivers.

15. Optional System Components

Connecting the Bolus Cord

(List 13027-01)

Connect the bolus cord to the port labeled "Bolus" on the bottom of the pump.

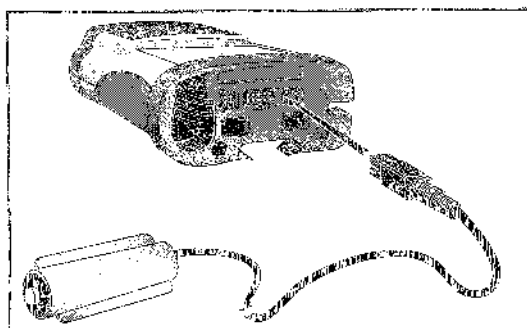

CAUTION: When removing the cord from the bottom of the pump, firmly grasp the connector and pull straight out. **DO NOT** twist or bend the cord or connector.

Using the GemStar Lockboxes


The GemStar Lockbox is an enclosure that can be locked with a key to secure the pump, and a flexible container or vial. When the pump is secured in the lockbox, the user has access to the AC mains connector, data and bolus ports, keypad, and bolus button.

Note: For the GemStar Universal Lockbox (List 12499-01) refer to the lockbox for instructions.


15. Optional System Components

30 mL PCA Vial Lockbox (List 13077-01)


- 1 Open the lockbox (A).
- 2 Flip the top door up (B).

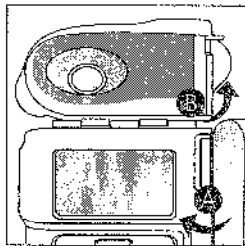


- 3 Align the pump with the lockbox (C).
- 4 Slide the pump down into the lockbox until fully seated.

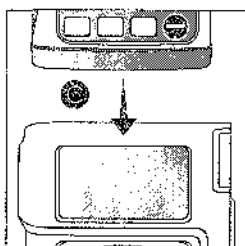


- 5 Flip the top door down.
- 6 Connect the pump set to the vial and prime.
- 7 Place the vial inside the lockbox (D).

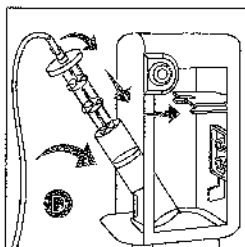


- 8 Load the cassette (E).

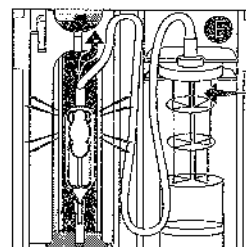

CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

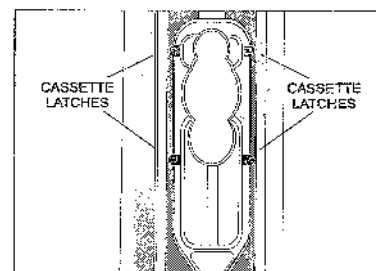

Note: Ensure the tubing is in the tubing guide (at the distal end).

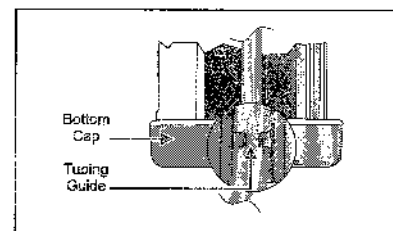


- 9 Ensure the tubing will not be pinched when closing the lockbox.
- 10 Close and lock the lockbox.
- 11 Remove the key.


15. Optional System Components

250 mL Lockbox (L11 13020-01)


- 1 Open the lockbox (A).
- 2 Flip the top door up (B).

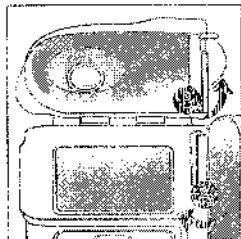


- 3 Align the pump with the lockbox (C).
- 4 Slide the pump down into the lockbox until fully seated.

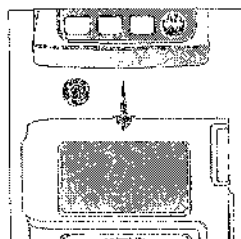


- 5 Flip the top door down.
- 6 Connect the pump set to the bag and prime.
- 7 Place the flexible container inside the lockbox with the administrative port down (D).

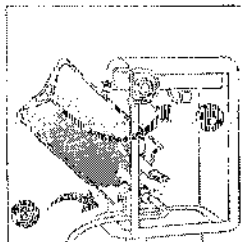


- 8 Hang the bag on the hook inside the lockbox (E).
- 9 Load the cassette (F).

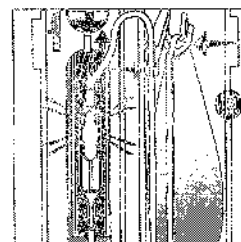

CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 31.

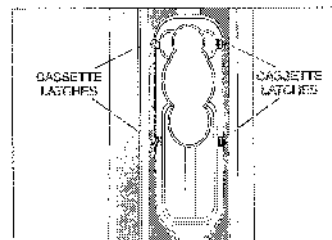

Note: Ensure the tubing is in the tubing guide (at the distal end).

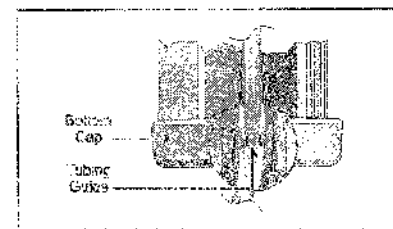

- 10 Ensure the tubing will not be pinched when closing the lockbox.
- 11 Close and lock the lockbox.
- 12 Remove the key.

15. Optional System Components

Removing the Lockbox

- 1 Remove the lockbox from the pole.
- 2 Unlock the lockbox.
- 3 Open the lockbox and flip the top door up.
- 4 Press the cassette release button to remove the cassette.
- 5 Lift the pump out of the lockbox.


Using the Hospira Carrying Cases

Ensure the administration set is fully primed and the cassette is properly loaded before placing the pump in the carrying case. Refer to Chapter 2 Basic Operation for more information.


Carrying Case - 250 mL (List 13594-01)

This soft-sided carrying case can be worn around the waist. It is designed to hold the GemStar Pump and up to a 250 mL flexible fluid bag.

- 1 Place the fluid bag in the zippered compartment (A) on the backside of the case. Ensure the bag spike is positioned so that the fluid flow is not restricted.

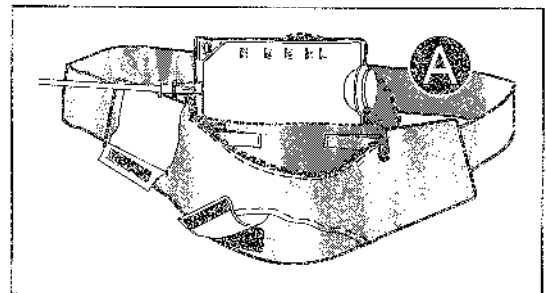


- 2 Place the pump inside the front compartment (B) of the case.

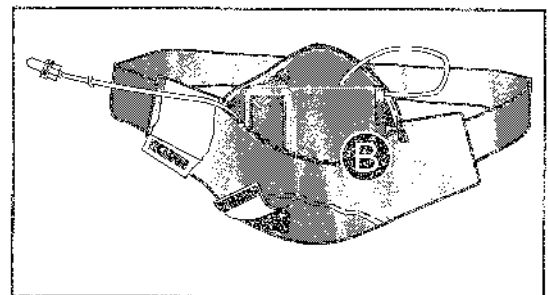


15. Optional System Components

- 3 Secure the tubing to the inside of the pack with the Velcro® tabs (C). Ensure the tubing is not kinked.
- 4 Use the large Velcro flap (D) on the upper left side of the case to secure the tubing on the outside of the case.

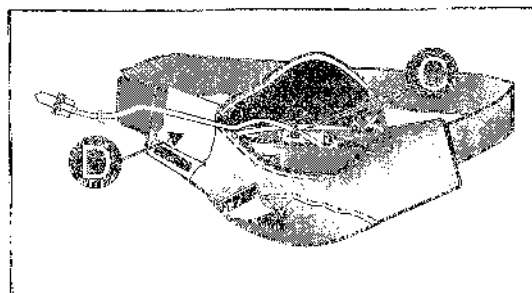

- 5 Close the zipper. Ensure the tubing exits the side of the case near the zipper tab when the zipper is closed (E).
- 6 Place the case around the waist. Snap the buckle together to secure.
- 7 Adjust the belt for a snug, comfortable fit.

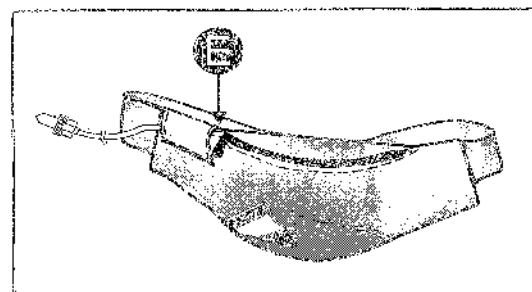


Carrying Case - 500 mL (List 13079-01) or 1000 mL (List 13980-01)


These soft-sided carrying cases can be worn around the waist or over the shoulder. They are designed to hold the GemStar Pump and up to a 500 mL (List 13079-01) or a 1000 mL (List 13980-01) flexible fluid bag.

15. Optional System Components

- 1 Place the pump on the inside front window (A). Ensure the keypad is facing the “window.”

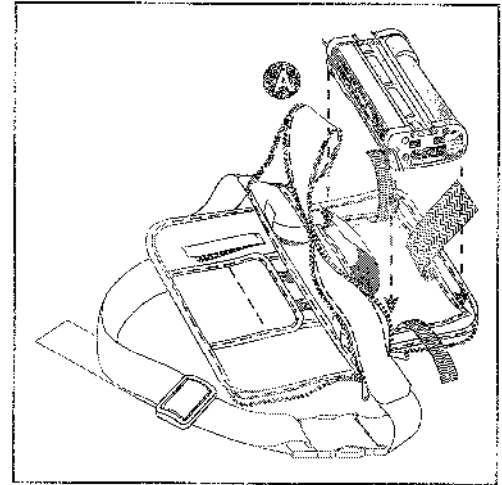


- 2 Secure the pump with the large side straps (B) and (C). Then secure the top (D) and bottom (E) straps.
- 3 Place the fluid bag on the backside of the compartment (F).
- 4 Secure the bag with flaps (G) and (H).
- 5 Secure the bag spike with strap (I). Ensure the fluid flow is not restricted.
- 6 Use a “figure-8” pattern to wrap the tubing. Secure the tubing to the inside of the case with tabs (J) and (K). Ensure the tubing is not kinked.
- 7 Zip the case closed. Ensure the tubing exits through one of the sides near the zipper tab when the case is closed.

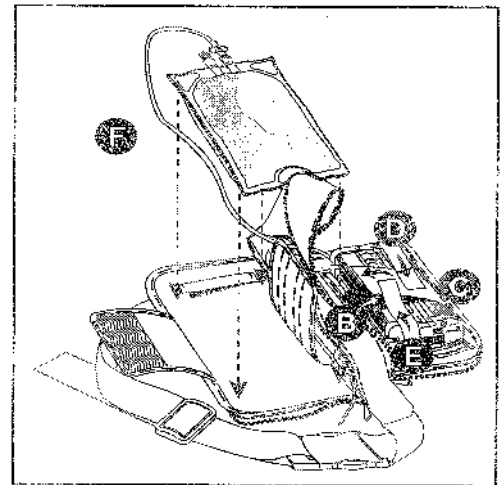

15. Optional System Components

- 2 Place the case around the waist with the front "window" facing away from the body. Snap the buckle together to secure. Adjust the belt for a snug, comfortable fit.

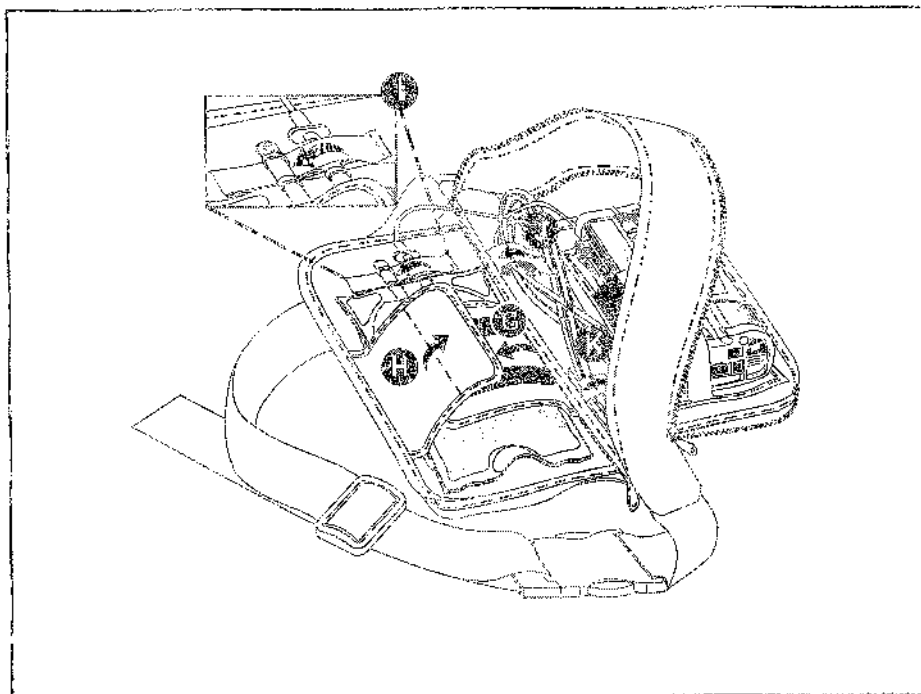

Carrying Case - 3 Liter (List 13590-01)

This soft-sided carrying case can be worn around the waist or as a backpack. It is designed to hold the GemStar Pump with batteries or a battery pack, and up to a 3 liter flexible fluid bag.

Note: This carrying case includes an interior support (Hospira 3 Liter Carrier List 13590-01) that may be removed from the case and used separately.


15. Optional System Components

- 1** Unzip the carrying case and release all Velcro straps and flaps.
- 2** Place the fluid bag on the carrying case as shown (A).
- 3** Hook top strap (B) through the hole in the top of the fluid bag and secure.
- 4** Fold small flaps (C) and (D) over the bag as shown and secure.
- 5** Secure straps (E) and (F) across the length of the fluid bag and over the small flaps.
- 6** Secure the spike connection with strap (G).
- 7** Fold large flaps (H) and (I) over as shown and secure.
- 8** Place the pump, with the keypad facing upward, inside the pocket (J) on the flap. The pump's display should be visible through the opening on the front of the pocket.
- 9** Secure the pump with straps (K) and (L).
- 10** Secure the tubing with strap (M).
- 11** Place the secured pump (N) inside the large zipper case or other carrying case.
- 12** Zip the large case closed.
- 13** If desired, route the tubing through the loops on the bottom of the carrying case. Ensure the tubing is not kinked.
- 14** Adjust the shoulder straps and waist belt for a snug, comfortable fit.



Tip: When not using the waist belt, secure it under the large flap on the back of the case.

15. Optional System Components

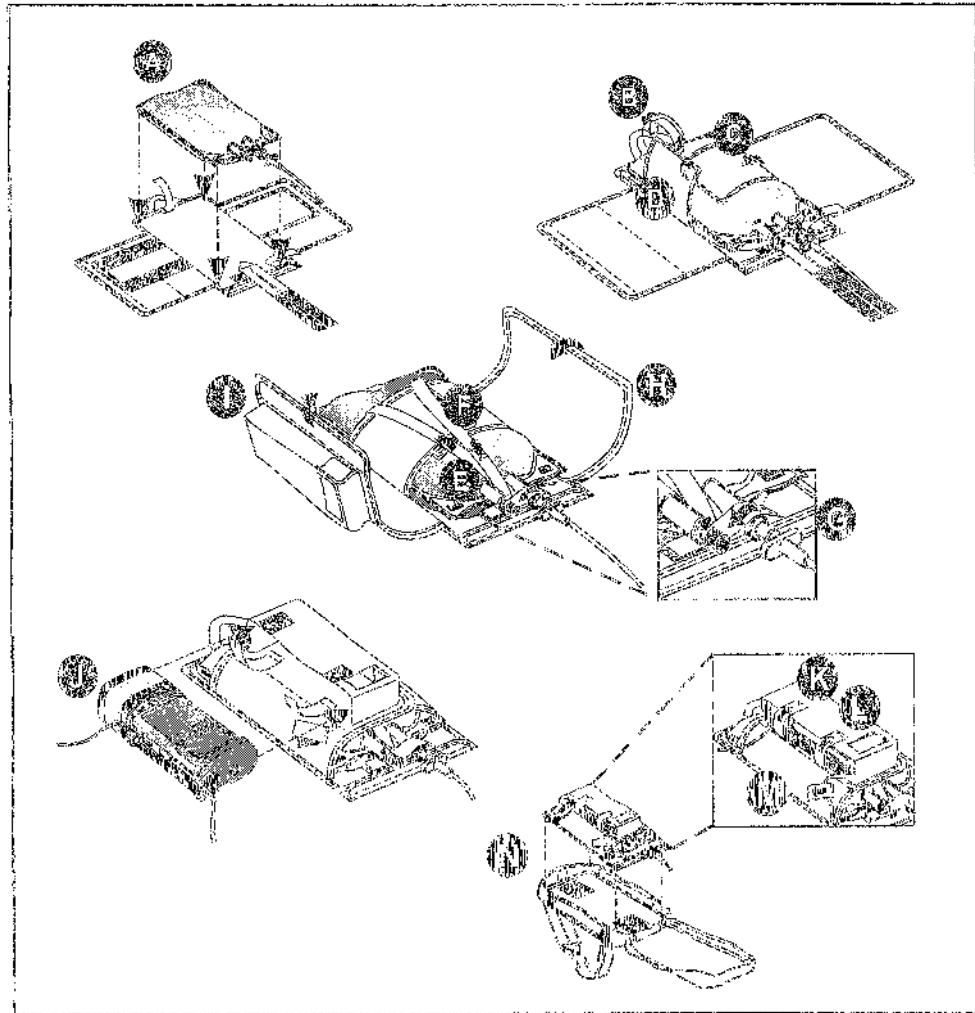


15. Optional System Components

Using the GemStar Docking Station

(List 13075-01)

Connecting the Pump

- 1** Mount the docking station to the I.V. pole and tighten the knob.
- 2** Align the connectors on the bottom of the pump with the connectors on the docking station.
- 3** Slide the pump into the docking station until it snaps into place.
- 4** Connect the docking station to AC mains power.

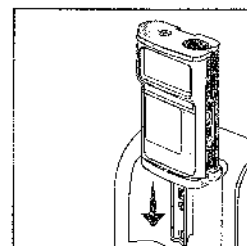

Note: The docking station light is yellow while the docking station is charging. When the docking station is fully charged, the light is green.


15. Optional System Components

Connecting the Bolus Cord and Printer Cable


WARNING: Ensure a cassette is installed in the pump before connecting to a printer or computer.

If the pump is connected to a patient: **BEFORE** connecting the pump to a computer or printer, press [STOP], close the slide clamp, and eject the cassette. Leave the pump set connected to the patient. Install another cassette in the pump before connecting to a printer or computer.

- 1 Plug the bolus cord into the bolus port on the rear, right side of the docking station.

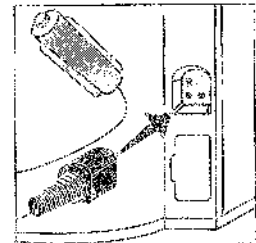


- 2 Plug the printer cable into the printer port on the rear, right side of the docking station.

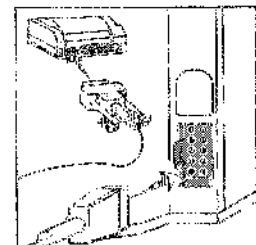

Disconnecting the Pump

To remove the pump:

- 1 Press the release button on the top of the docking station.
- 2 Slide the pump out.

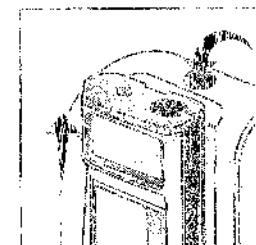

Chapter 16

Pump Sets

CAUTION: USE ONLY GemStar Pump Sets with the GemStar Pump. Use of unauthorized sets may result in injury to the patient or damage to the pump.

Contact a Hospira Representative for appropriate pump set configurations.

Fluid path and areas beneath undisturbed protective set covers are sterile and nonpyrogenic in the intact unit package.

STERILE ②

The GemStar Pump Sets, listed in the table on the following page are available for use with the GemStar Pump.

Note: The disposable pump sets should be changed per the Centers for Disease Control (CDC) guidelines or health care provider policy. Discard after use. All GemStar Pump Sets are composed of latex-free components only.

GemStar Pump Sets contain an anti-siphon valve to reduce the potential for gravity flow at a typical 36-inch head height.

16. Pump Sets

The following sets are a representation of the GemStar® Pump Sets currently available:

List No	Description
Ambulatory Sets	
13263	Latex-Free GemStar® Pump Set and Extension Set with integral pressure-activated anti-siphon valve, Secure Lock, Nonvented, Microbore tubing, 96 in (244 cm)
13264	Latex-Free, LifeShield®, GemStar® Pump Set with 0.2 Micron Filter, Secure Lock, Prepierced reseal injection site, Nonvented, Pre-attached pressure-activated anti-siphon valve, Microbore tubing, 96 in (244 cm)
13265	Latex-Free, LifeShield®, GemStar® Pump Set with 1.2 Micron filter, Secure lock, Prepierced reseal injection site, Nonvented, Pre-attached pressure-activated anti-siphon valve, NonPhthalate (Non-DEHP) microbore tubing, 96 in (244 cm)
13273	Latex-Free GemStar® Pump Set, Secure lock, Non vented, Pre-attached pressure-activated anti-siphon valve, Microbore tubing, 96 in (244 cm)
General Infusion Sets	
13244	Latex-Free, GemStar® REMA® II, Y-Type Blood Set, 210 Micron Filter, Lower CLAVE® Connector, Macrobore tubing, 110 in (279 cm)
13267	Latex-Free, LifeShield®, GemStar® Primary Piggyback Pump Set, Secure Lock, Convertible, Macrobore Tubing Piercing Pin, 2 CLAVE® connectors, Pre-attached pressure-activated anti-siphon valve, Macrobore tubing, 110 in (279 cm)
13272	Latex-Free, LifeShield®, GemStar® 150 x 60 micro drip SoluSet®, Option Lock, Convertible Piercing Pin, 2 CLAVE® connectors, Macrobore tubing, 110 in (279 cm)
13274	Latex-Free, LifeShield®, GemStar® Piggyback Pump Set, with 0.2 Micron Filter, Secure Lock, Convertible Piercing Pin, 3 CLAVE® Connectors, Pre-attached pressure-activated anti-siphon valve, Macrobore tubing, 110 in (279 cm)

16. Pump Sets

List No	Description
Pain Management Sets	
13045	Latex-Free, GemStar® Pump Set, Integral pressure-activated anti-siphon valve Secure Lock, Nonvented with PCA injector and female Luer Lock™ adapter, Microbore tubing, 96 in (244 cm)
13047	Latex-Free, GemStar® Pump Set with PCA Extension, Integral pressure-activated anti-siphon valve, Secure Lock, Nonvented with PCA injector and female Luer Lock™ adapter, Microbore tubing, 96 in (244 cm)
13261	Latex-Free, GemStar® Pump Set and Extension Set with integral pressure-activated anti-siphon valve, Secure Lock, Non-vented, Microbore tubing, 110 in (279 cm), Yellow Striped Tubing

Notes

GemStar Pump Specifications

Physical Dimensions: Weight:	The GemStar Pump is composed of latex-free components only. Height: 5.5 in. (14.0 cm) Width: 3.8 in. (9.7 cm) Depth: 2.0 in. (5.1 cm) Approximately 17 oz. (482 grams) excluding batteries
Transport and Storage Environment Ambient Temperatures: Relative Humidity: Atmospheric Pressure:	Store in cool and dry environment. -4° F to +140° F (-20° C to +60° C) 10% to 90% 0–10,000 ft. (0–3,000 m) equivalent pressure
Operating Environment Ambient Temperatures:	+41° F to +104° F (+5° C to +40° C)
Power Sources AC Mains: List 13026-01 (USA) complies with UL 544 List 13072-01 (Universal) complies with CSA 22.2 No. 601.1 and UL 60601-1 Battery: Battery Pack: Docking Station:	Wall-mount AC mains adaptor; 6 ft. (1.8 m) cord molded plug Input: 110 VAC Output: 3.0 VDC Tabletop AC mains adaptor; molded plug input: 100–240 VAC Output: 3.0 VDC Two disposable AA alkaline or lithium batteries Rechargeable using GemStar AC Mains Adaptor Mounts on an I.V. pole; 12 ft. (3.6 m) cord Input: 100–240 VAC Output: 3.0 VDC Internal rechargeable battery

17. GemStar Pump Specifications

Power Capacity	Using two fresh, disposable AA batteries or a charged battery pack, at temperatures between 59° F and 95° F (15° C and 35° C), the pump is capable of delivering approximately: 96 hours at rates below 5 mL/hr 48 hours at rates at or above 5 mL/hr but below 25 mL/hr 24 hours at rates at or above 25 mL/hr but below 125 mL/hr 3000 mL at a rate of 125 mL/hr or higher
Pump Mechanism	Volumetric, Piston Driven
Memory Protection	Current program and 400-event (history) log protected by internal lithium battery-backed memory for at least one year after power is removed
Operating Controls	One 23-key keypad; bolus, data port, and AC mains jacks; bolus button
External Power LED	Green LED marked with plug icon lights continuously when pump is connected to AC mains power Green LED flashes when pump is connected to external batteries
Alarm LED	Red LED marked with alarm icon lights during alarm conditions
Audible Alarm	The audible alarm is user-adjustable from the maximum volume down to silent The alarm automatically reverts to the maximum volume if the user does not respond within one minute
Display Backlight on AC: Backlight on Batteries:	Backlit, four-line-by-sixteen-character alphanumeric graphics display Continuous Continuous during programming, program review, and history display; otherwise activated by keypress or alarm; not activated by bolus request
Real Time Clock	Accuracy of ± 5 minutes per month or better at temperatures between 59° F and 95° F (15° C and 35° C)
Print Function Port & Interface: Printers:	RS-232 serial interface port (data port), minimum baud rate of 2400, isolated circuit Seiko DPU 414 or compatible serial printer
System Accuracy	$\pm 10\%$ for rates of 0.1 mL to less than 5 mL/hr $\pm 5\%$ for rates of 5 mL to 1000 mL/hr

17. GemStar Pump Specifications


Air Sensitivity	ON: Pump alarms at approximately 0.5 mL of air Alarms for any bubble greater than 500 microliters with a tolerance of 200 microliters 2 mL: Pump alarms at approximately 2 mL of air Alarms when pump detects 2.0 +1.0/-0.2 mL of air in 6 mL of total volume delivered OFF: Alarm is not activated Recommend using air-eliminating filter when Air Sensitivity is set to OFF
Occlusion Sensitivity	Distal Occlusion: Alarms when pump detects distal pressure within a selectable range: Low: 7 psi (48 kPa) $\pm$ 5 psi ($\pm$34 kPa) Medium: 12 psi (83 kPa) $\pm$ 8 psi ($\pm$55 kPa) High: 26 psi (179 kPa) $\pm$ 14 psi ($\pm$96 kPa) Proximal Occlusion: Alarms when pump detects proximal pressure less than or equal to -4 psi (-28 kPa)
Piggyback Rate	0.1 to 300 mL per hour
Pump Self-Tests and Safety Features	Self-test performed when the power switch is activated Diagnostic routine, including motor speed and air-detection monitoring, is repeated continuously while the pump is powered on Error and alarm conditions are indicated by both audible and visual alarms; delivery in progress is stopped, if appropriate (refer to Chapter 12 Troubleshooting)

17 GemStar Pump Specifications


Occlusion Information

Avoiding Bolus Infusion After Occlusion


- 1 Press [STOP].
- 2 Close the distal slide clamp or CAIR[®] (roller) clamp.

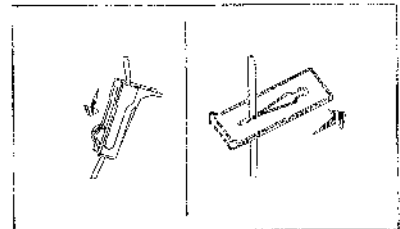


- 3 Press the cassette release button to eject the cassette from the pump.

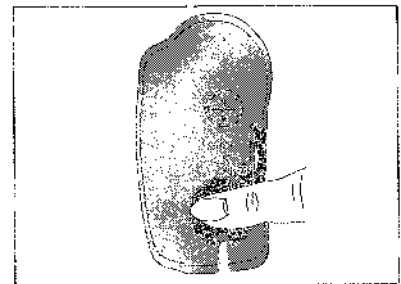

- 4 Press down on (S) to open the flow stop.
- 5 Wait 10 seconds.

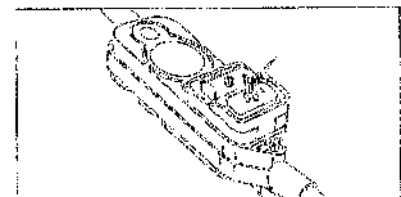



- 6 Press down on (S) to close the flow stop.

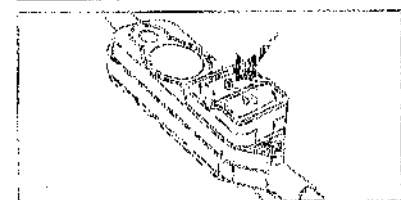


17. GemStar Pump Specifications

- 7** Load the cassette.


CAUTION: Ensure all four cassette latches are clearly visible after the cassette is installed. Refer to *Loading the Cassette* on page 32.

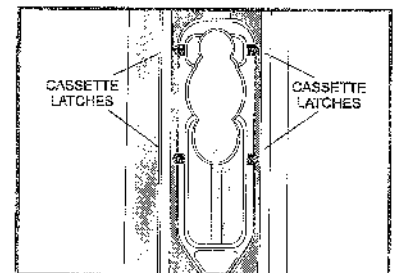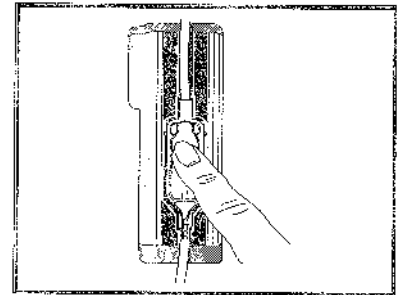

Note: Ensure the tubing is behind the tubing guide (at the distal end).

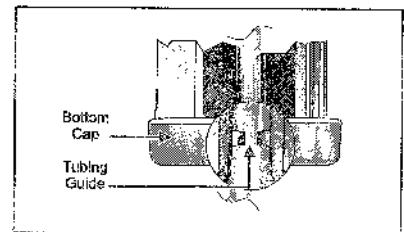

- 8** Open the distal slide clamp or CAIR[®] (roller) clamp.
- 9** Review the program. Press  to begin review.
- 10** Press [START] to start the pump.

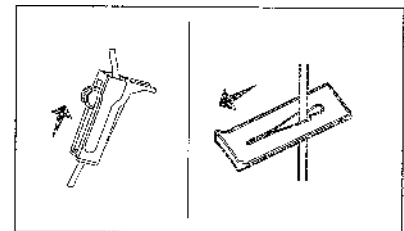

17. GemStar Pump Specifications

Stored Occlusion Volume

Delivery Rate of 25 mL/hr at	Pressure Limit		Stored Volume in mL		Time to Occlude in seconds	
Pressure Setting	Macro	Micro	Macro	Micro	Macro	Micro
Low	7 psig (48.3 kPa)	7 psig (48.3 kPa)	0.26	0.09	50	15
Medium	14 psig (96.5 kPa)	13 psig (89.6 kPa)	0.54	0.14	93	23
High	32 psig (220.6 kPa)	30 psig (206.8 kPa)	1.24	0.40	254	83

Delivery Rate Accuracy

The following trumpet curves represent the typical flow rate deviations, both positive and negative, from the programmed flow rate under test conditions defined by International Electrotechnical Commission (IEC) Standards for infusion pumps. Using these curves, the medical professional may determine if the device can be expected to perform in a manner suitable for the drug or fluid to be infused. The pump sets used in this test are comparable to GemStar Pump Set List 13273.

The typical accuracy for this device is $\pm 5\%$ of the programmed rate for most conditions, after a period of normalization. It is recommended that the medical professional refer to the following curves when making decisions regarding drug and fluid administration.

A single fault condition occurs when a single means for protection against a safety hazard in equipment is defective or a single external abnormal condition is present.

The maximum over-infusion volume from a single fault condition is based upon the calculation of the maximum time to generate an alarm times the maximum rate possible (1000 mL/hr) as a consequence of the fault. This has been determined not to exceed 2.78 mL.

CAUTION: When using pre-filled vials, the system may underdeliver at low flow rates.

Note: False occlusion alarms and fluid delivery inconsistencies may occur when pumping viscous fluids in microbore tubing at rates greater than 500 mL/hr.

Trumpet Curves

The Trumpet Curve Graphs following the Example show representative maximum and minimum percent flow rate deviation from the programmed rate over time. This information was developed in accordance with IEC 60601-2-24: 1998, Sub-Clause 50.102. Refer to this standard for detailed information.

How to read a Trumpet Curve Graph (Refer to example on the following page):

The graphs following the Example plot flow rates at 30 second intervals for the first 2 hours and for the 96th hour of delivery. The graph plots mean delivery rate error (Average of 3 pumps) for the 2nd hour and the 96th hour as a straight line. The graph also presents maximum and minimum average delivery rate error for this interval plotted by averaging delivery errors over intervals of 2, 5, 11, 19 and 31 minutes ("Trumpet Curve").


Example

From the Trumpet Curve Graph sample that follows, find the 5 minute interval (A) at the horizontal axis and read the corresponding points (B) and (C) on the vertical axis. The values are approximately +2.8% and -0.5%.

This means that at the rate of 25 mL/hr the average maximum flow rate fluctuation for any 5 minute time interval during the 2nd hour of operation was within the limits of +2.8% and -0.5% from the nominal rate. The average delivery rate error over the entire 2nd hour was +1.6% (D).

For other time intervals look at other points at the horizontal axis and determine corresponding limits as above.



17. GemStar Pump Specifications

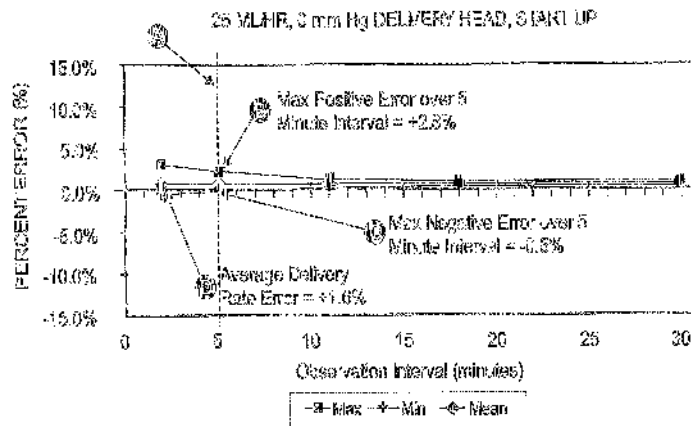



A trained professional can use the resulting graphs to select a pump with the appropriate startup and flow characteristics to suit the clinical application.

Note: As an example of how the trumpet curves can be used, consider the maximum and minimum deviations at the 5 minute average interval. The upper curve provides the maximum expected delivery rate error over a 5 minute interval, the lower curve provides the minimum expected delivery rate error over a 5 minute interval. An example would be Dopamine administered at 5 µg/kg/min. At 5 minutes, the average drug delivery error would be within the range of +2.8% and -0.5% of the expected nominal rate.



17. GemStar Pump Specifications

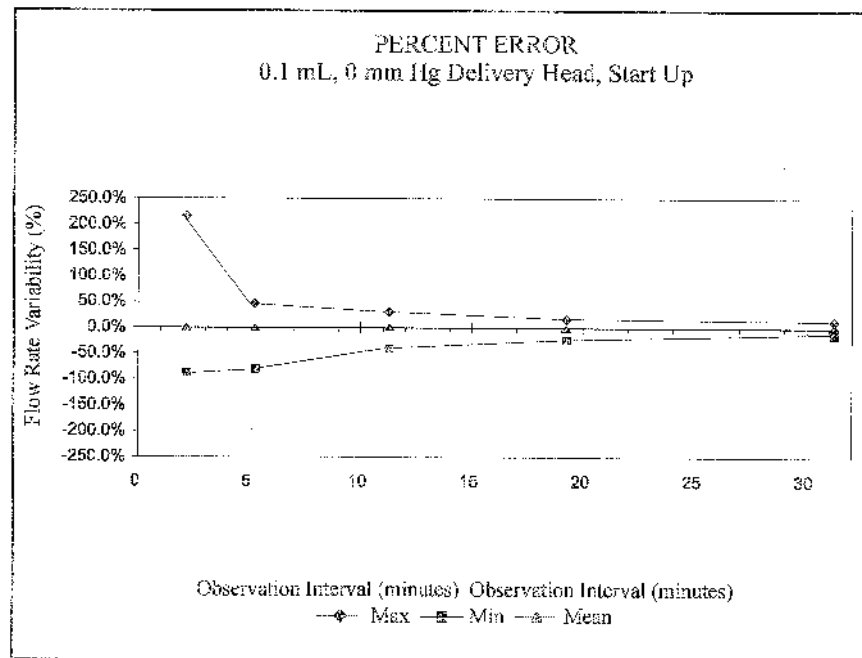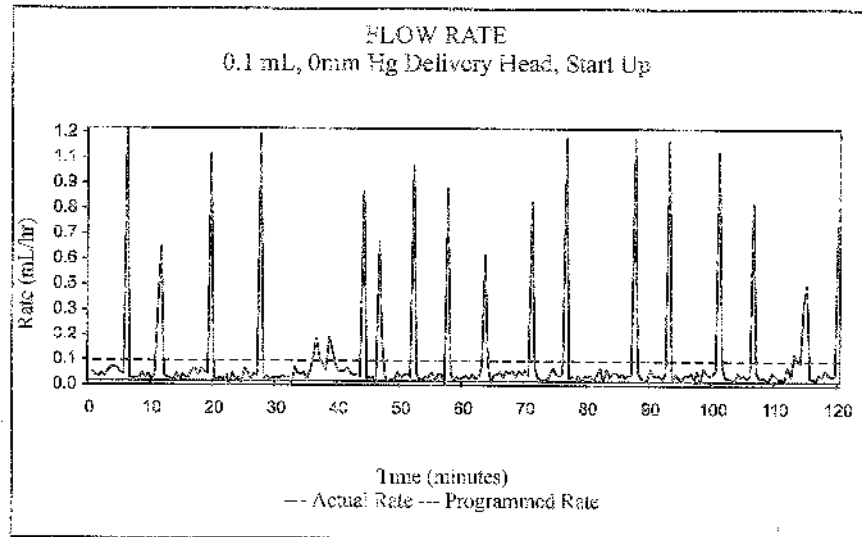



17 GemStar Pump Specifications

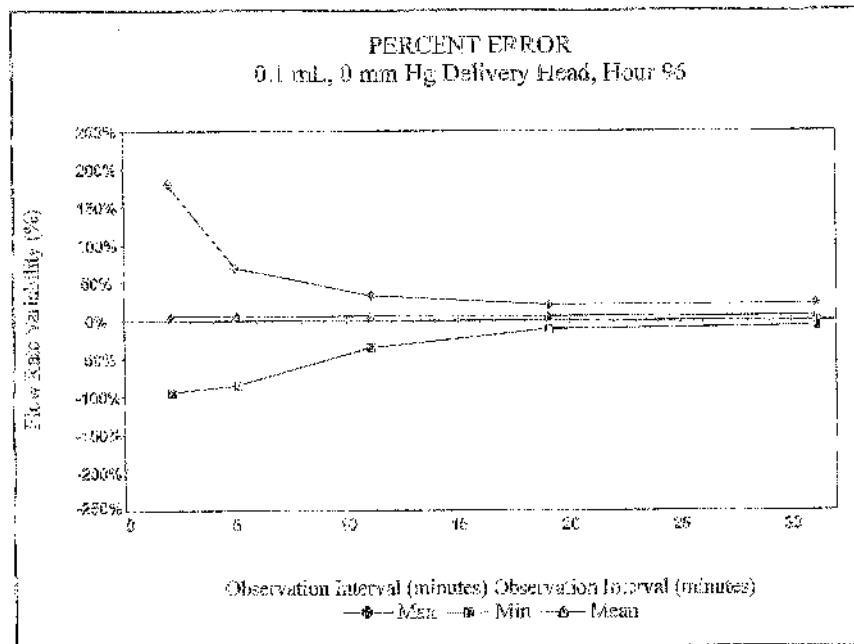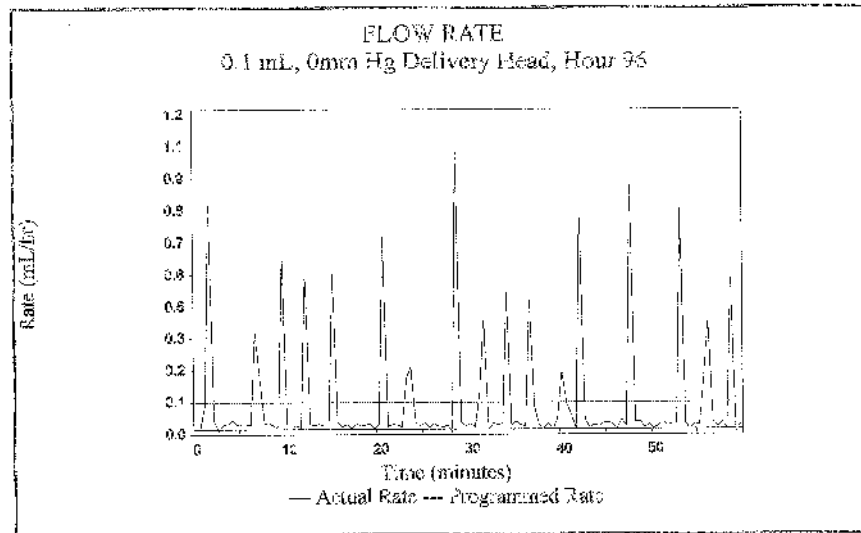



17. GemStar Pump Specifications

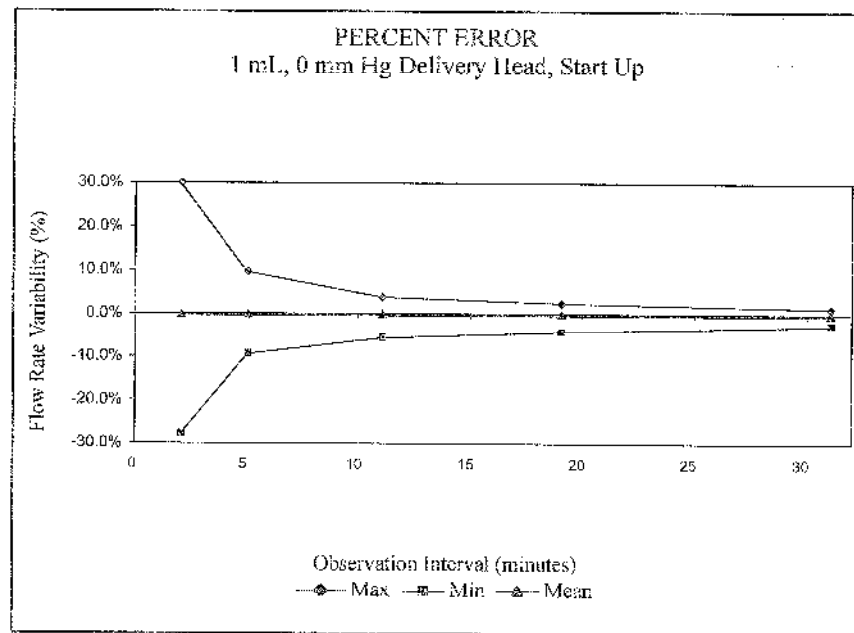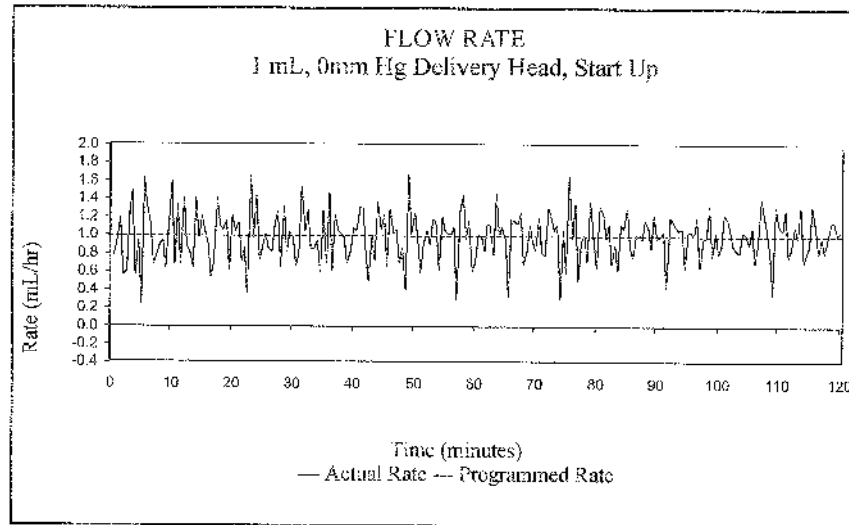



17. GemStar Pump Specifications

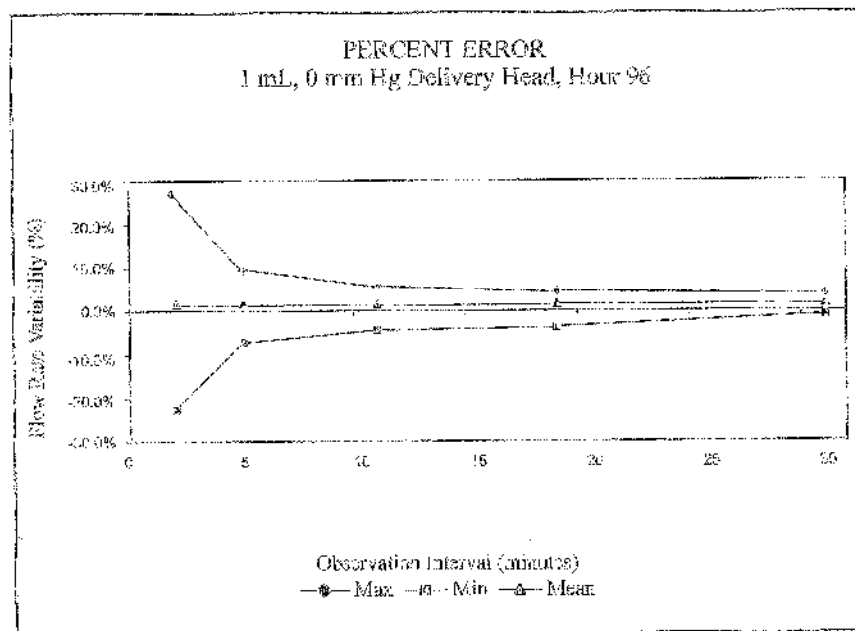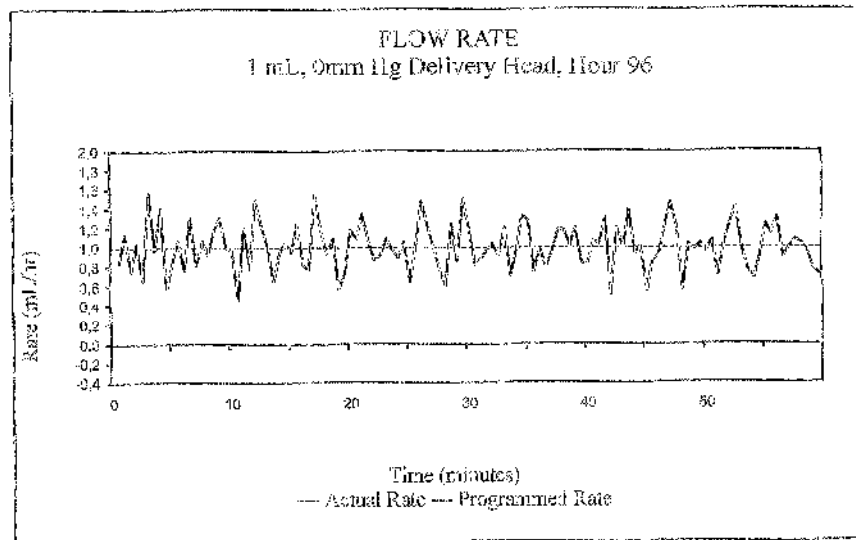



17. GemStar Pump Specifications

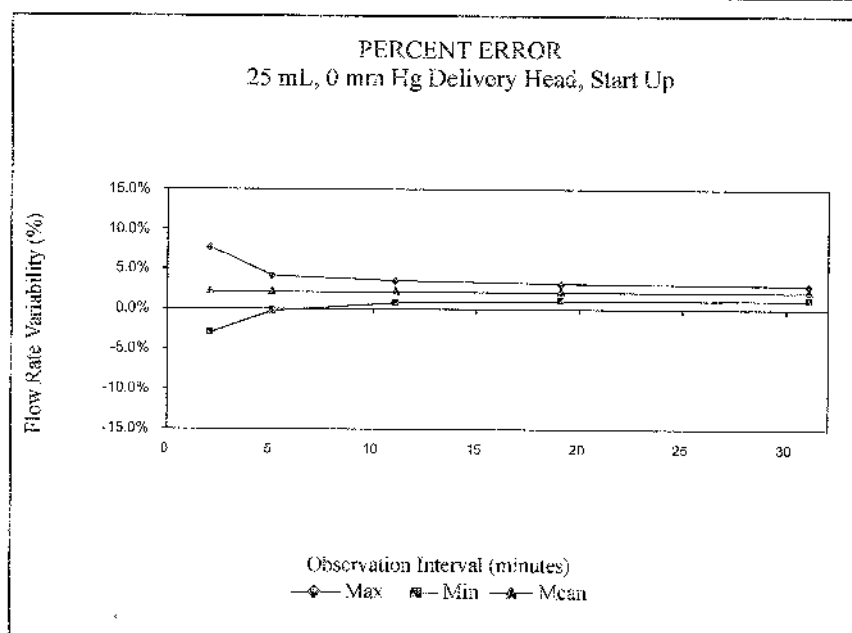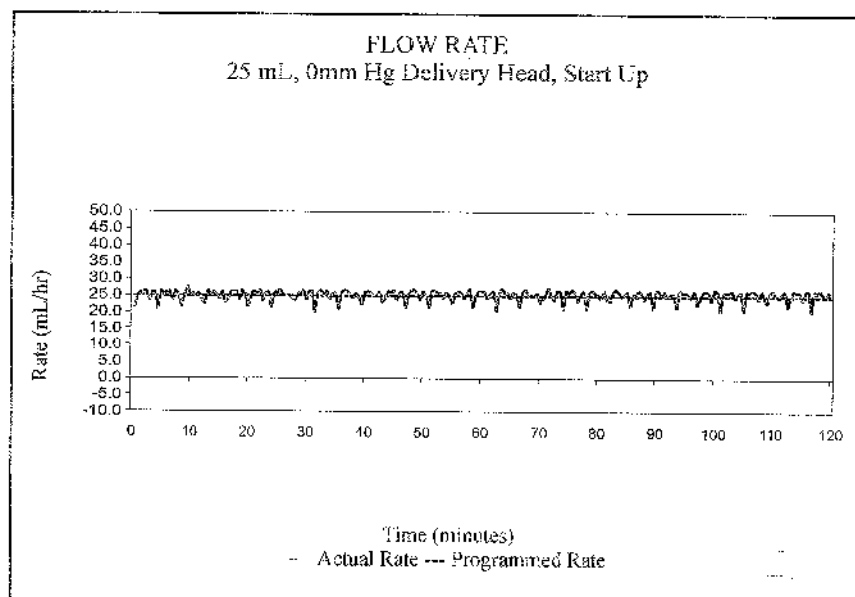



17 GemStar Pump Specifications

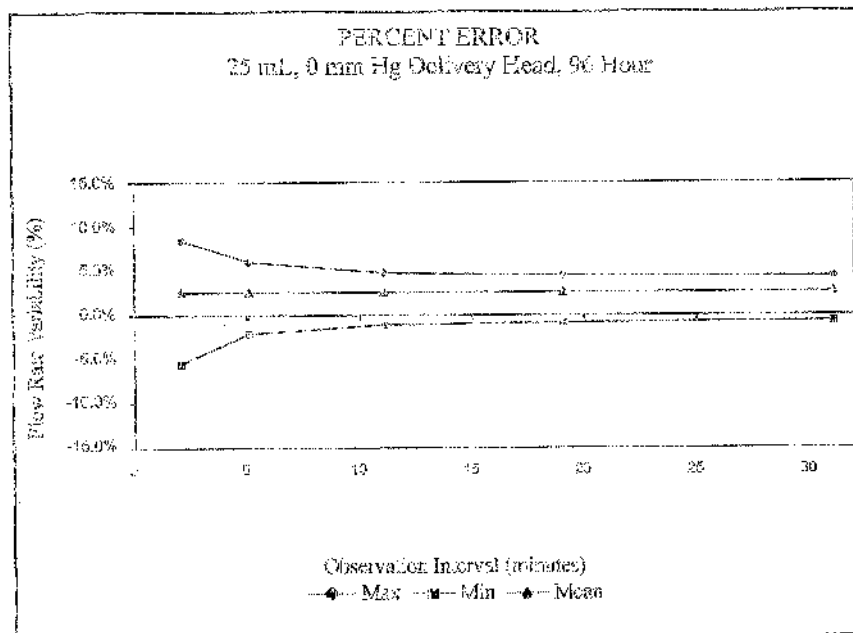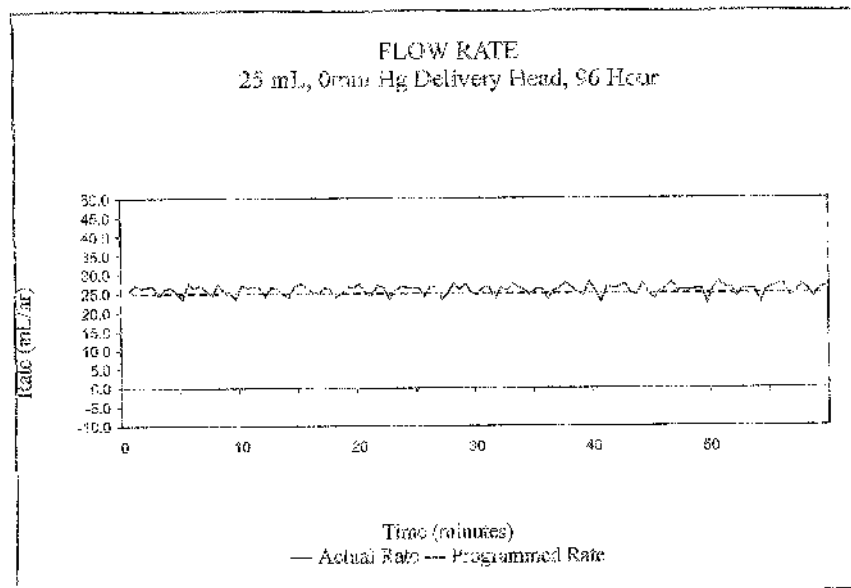



17. GemStar Pump Specifications

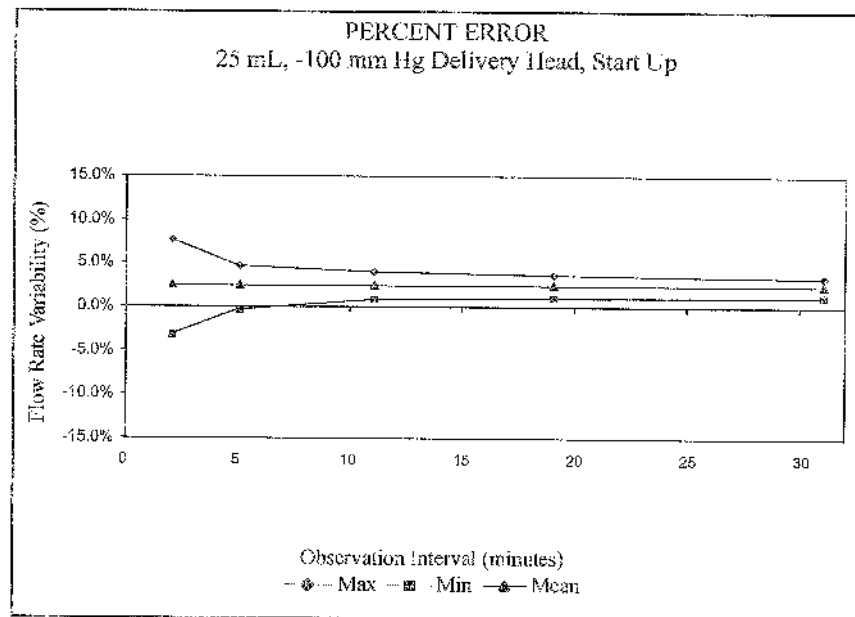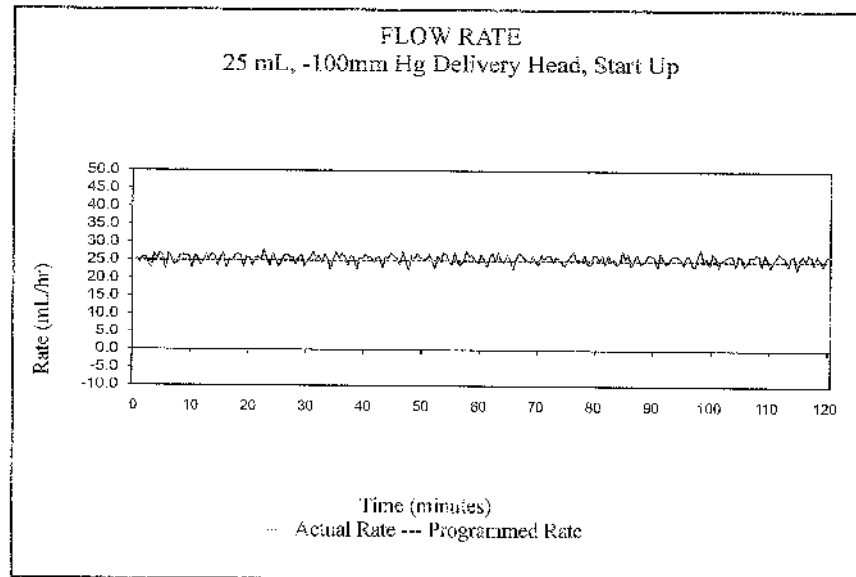

17. GemStar Pump Specifications

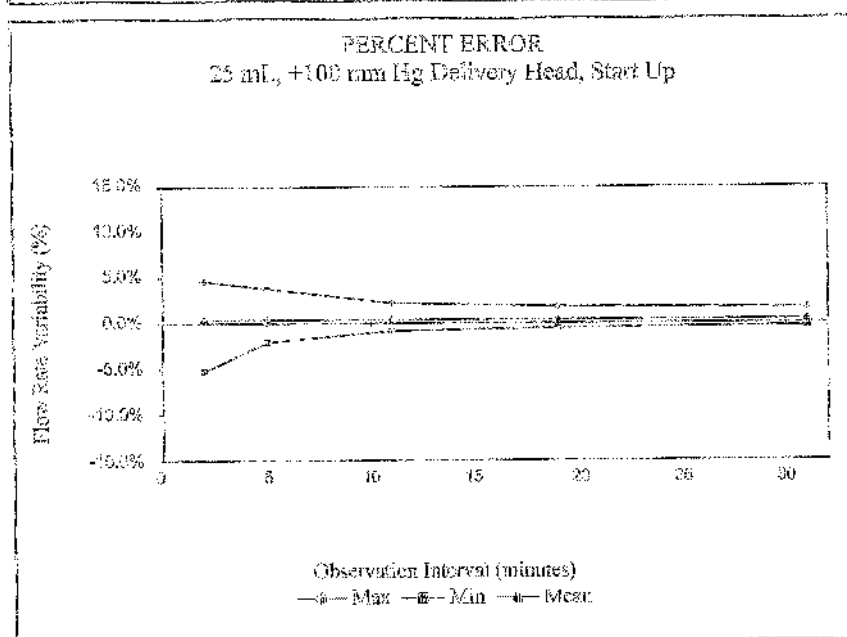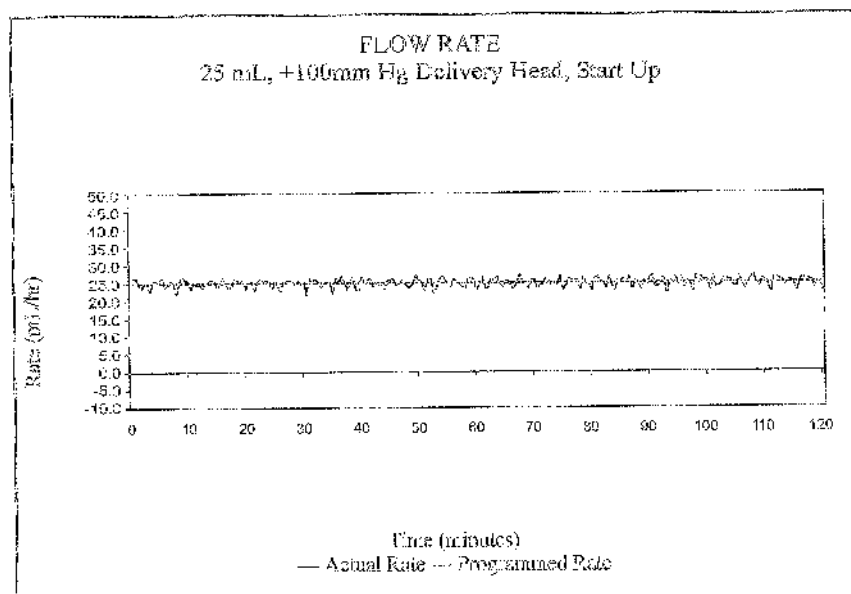

17. GemStar Pump Specifications

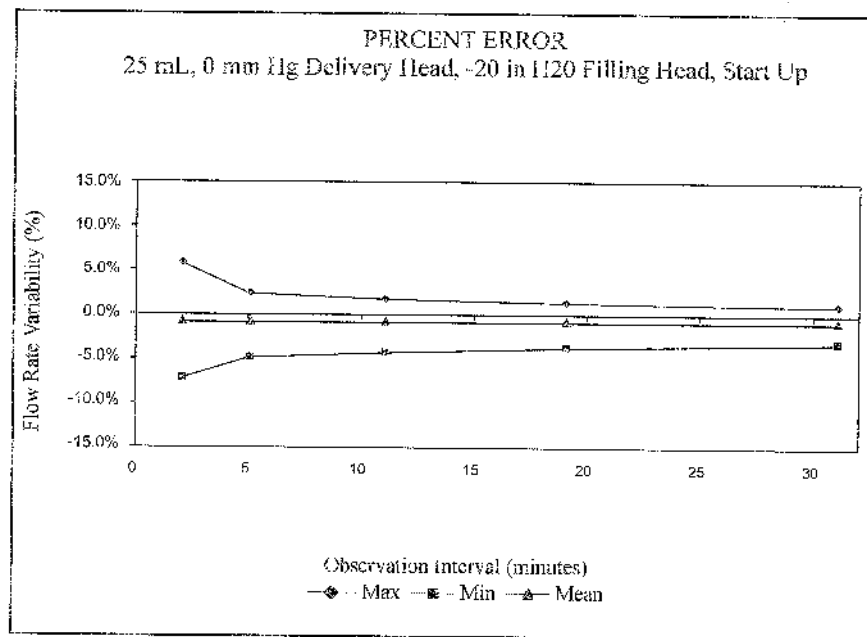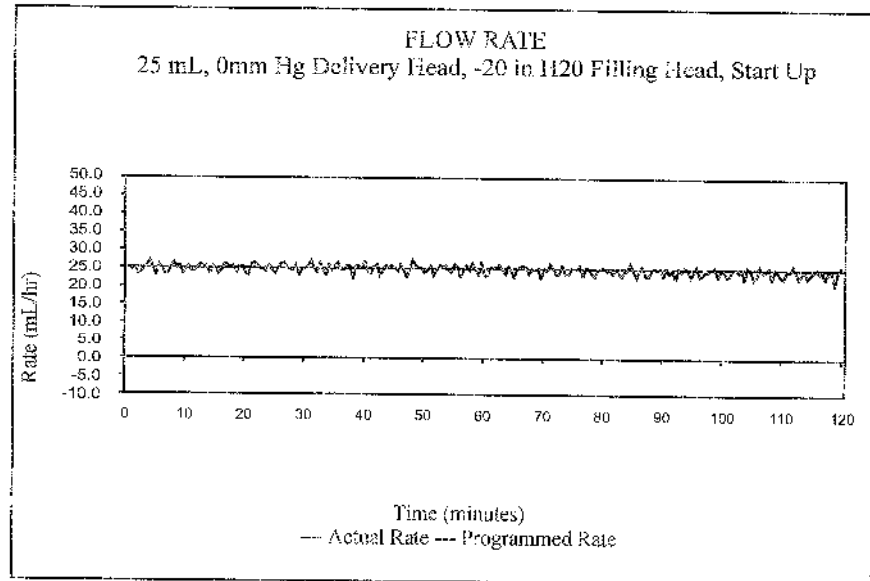

Notes

Chapter 18

Warranty/Technical Assistance

The GemStar Pump has been carefully manufactured using high-quality components. It is warranted to be free from defects in material and workmanship for a period of one year from the date of purchase under normal use and service. The warranty on the optional battery pack and other accessories is limited to 90 days.

On return to Hospira, the pump will be repaired or replaced within the terms of this warranty. Material returned to Hospira must be properly packaged and sent freight prepaid.

This warranty shall not apply if the pump has been repaired by anyone other than Hospira qualified service personnel; or altered in any way which, in Hospira's judgment, affects its stability or reliability; or if the serial number has been altered, effaced, or removed.

THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED.

18. Warranty/Technical Assistance

For customer service within the United States, contact:

1-877-9HOSPIRA (1-877-946-7747)

For technical assistance and product return authorization within the United States, contact:

1-800-241-4002

For additional technical assistance, including Technical Service Bulletins, technical training, and product information, visit the Web site at:

www.hospiraparts.com

After authorization, ship prepaid product returns to the following address:

Hospira, Inc.
Technical Support Operations
755 Jarvis Drive
Morgan Hill, CA 95037

 *Note: Outside the U.S., contact your local Hospira sales office.*

Index

Numerics

1-hour dose limit, 55

4-hour dose limit, 55

A

AC mains

power connector, 4

power indicator, 3

Pump Specifications, 171

regarding handling and maintenance, 13

regarding pump operation, 13

AC Mains Adaptor

connecting to, 20

administration set. *See* pump set

air sensitivity, 77, 83, 107

regarding air-in-line infusion, 13

specifications, 173

air-in-line, alarm, 11, 13

test of, 149

alarm messages, 128–134

alert messages, 135–137

anti-siphon valve, 12, 30, 31

Auto KVO, 5

in mL/hr Only therapy, 85

in Weight Dosed therapy, 79

Auto Taper Down

in TPN therapy, 50

See also taper down; taper up

auto tapering. *See* Auto Taper Down; taper down; taper up

B

backcheck valve, 71, 72, 86, 87

BACK-UP key, 3, 28, 97, 102

base rate, 92, 93, 94

batteries, 115

alarm messages, 128–134

battery compartment, 4

battery pack connector, 4

installing disposable batteries, 19

maintenance, 139

power check, 146

specifications, 171

using the rechargeable battery pack, 21

warnings and cautions, 13–16

baud rate, 172

configuration settings for the Data

Port, 113, 118

switch settings for Seiko printer, 114

bolus, 5, 10, 57

bolus button, 4

bolus delivery, 56–57

bolus dose interruptions, 57

clinician instructions, 122, 123

connecting the bolus cord, 156, 166

history, 98

occlusion information, 174

optional system components, 153

in Pain Management therapy, 53–57

print, 98, 109

in Variable Time therapy, 91, 94

1-hour or 4-hour dose limit, 55

Index

C

- call back alarm, 51
- carrying cases
 - using the Abbott carrying cases, 159–164
- cassette, 30
 - cassette pocket, 4
 - cassette release button, 4
 - cleaning of, 139–141
 - error message, 129
 - inspection of, 143
 - loading the cassette, 32
 - releasing the cassette, 33
 - warnings and cautions, 13
- Change History, ii
- CHANGE key, 3, 28
- Change Menu, 8
- clinician instructions, 121–125
 - adding a second lock level, 123
 - locking the pump, 122
 - program and deliver a clinician activated loading dose, 124
 - unlocking the pump, 124
- configuration settings, data port, 113
- connecting the Abbott GemStar AC Mains Adapter, 20
- connecting the bolus cord, 156
- connecting to a computer, 116
 - See also* Data Port
- connecting to the Seiko printer, 115
 - See also* Data Port
- continuous delivery
 - in Pain Management therapy, 53–58
 - in TPN therapy, 47, 50
- Continuous therapy, 5, 69–75
 - changing a Continuous Program, 74
 - Clinician instructions, 122, 123
 - features and specifications, 69
 - piggybacking in Continuous therapy, 71–74
 - programmable ranges, 75

- programming of, 69
 - Programming Worksheet, 75
 - rate titration, 71
- contraindications for use, 11

D

- Data Port, connecting to, 112–117
 - configuration settings, 113
 - connecting to a computer, 116
 - connecting to the Seiko printer, 115
 - connector, 4
- disconnecting the pump, 23
- distal end, 30, 35, 147
- distal occlusion, 6, 106
 - alarm message, 130
 - occlusion sensitivity, 173
 - occlusion sensor setting, 107
 - test of, 148
 - troubleshooting, 130
- docking station, 154
 - connecting the bolus cord and printer cable, 166
 - connecting the pump, 165
 - disconnecting the pump, 166
- dose limit, 55
- drip chamber
 - priming pump sets with a drip chamber, 38
 - priming pump sets with a drip chamber using the PURGE Key, 40

E

- event log, 6
 - history, display, 100
 - print, 98

F

- flow stop, 30, 33, 34, 36, 38, 40, 42, 45, 174
- Full Lock, 104
 - lock level functions, 105

locking the pump, 122
unlocking the pump, 124

G

GemStar pump components, 4
GemStar pump features, 6
GemStar pump layout, 2–3

H

Heading 1, 17
HELP key, 3, 9, 28, 47, 54, 61, 69, 78, 84,
91, 101, 125, 127
HELP mode, 9
history, 6, 56, 57, 98, 100, 139
display bolus history, 103
display, clear history, 102
print, bolus history, 98, 109

I

Intermittent therapy, 5, 61–67
changing an Intermittent therapy, 66
changing the start time, 98
delivery interruptions, 64–65
features and specifications, 61
Next Dose function, 111
programmable ranges, 67
programming of, 62–63
Programming Worksheet, 67

K

Keep Vein Open. *See* KVO; Auto KVO
KVO, 5
in Continuous therapy, 69
in Intermittent therapy, 61, 67
in mL/hr Only therapy, 85
in TPN therapy, 47
in Variable Time therapy, 91
in Weight Dosed therapy, 77, 79
troubleshooting, 138

L

loading dose, 5, 55
clinician activated loading dose, 124
delivery, 57–58
in Pain Management therapy, 53
loading the cassette, 32
lock, 3
adding a second lock level, 123
Full Lock, 104, 105
keypad lock, 6, 98, 104
keypad lock level functions, 105
locking the pump, 122
Luer Lock, 12
Rate Change Lock, 58, 104, 105, 122, 123
unlocking the pump, 124
lockboxes
using the Abbott GemStar
lockboxes, 156–159
locking the pump, 122
adding a second lock level, 123
Luer Lock, 12

M

maintenance, 139–141
cleaning and disinfecting, 139
cleaning the cassette pocket and tubing
channel, 141
cleaning the pump case, 140
recommended cleaning solutions, 140
repair, 141
Menu
Programming, Change, Options, 7
mL/Hr Only therapy
VTBI Titration, 86
mL/hr Only therapy, 5, 83–90
Auto KVO, 85
changing a mL/hr Only program, 89
features and specifications, 83
piggybacking, 86–89
programmable ranges, 90

Index

programming of, 83, 84
Programming Worksheet, 90
rate titration, 85

N

new container, 6
 starting a new container, 29
 troubleshooting, 136, 138
New Container Lock
 locking the pump, 122
 unlocking the pump, 124
next dose, 92, 111

O

occlusion, 106
 changing the occlusion sensor
 settings, 107
 occlusion information, 174–176
 occlusion sensitivity, 173
 See also distal occlusion; proximal occlusion
occlusion sensor
 changing the occlusion sensor
 settings, 107
 See also occlusion; distal occlusion;
 proximal occlusion
Operating Modes
 Help, Programming, Run, Stop, 9
operation test, 143–151
 air-in-line test, 149
 cassette pocket and latch inspection, 145
 display test, 146
 distal occlusion test, 143
 equipment required, 143
 keypad test, 146
 LED test, 147
 Operation Test Checklist, 151
 power check, 146
 printing the results, 149
 proximal occlusion test, 148

test setup, 144
volume accuracy test, 147
optional system components, 153–156
OPTIONS key, 3
 display shift totals, 161, 125
Options Menu, 8, 97–117
 change Pump Settings, 106
 connecting to the data port, 112–117
 Histories, 99
 Keypad Lock, 104
 Next Dose, 111
 Operation Test, 103
 Print, 109
 Review Program, 99
 Set Clock, 103
 Speed Protocol, 110
 tips for using the options menu, 97
other display messages, 138

P

Pain Management therapy, 5, 53–60
 1-Hour or 4-Hour Dose Limit, 55
 alert messages, 135
 bolus delivery, 56
 changing a Pain Management therapy, 58
 features and specifications, 53
 loading dose delivery, 57
 programmable ranges, 60
 programming of, 54–56
 Programming Worksheet, 60
piggybacking, 5, 6
 canceling a piggyback delivery, 74, 89
 changing a piggyback delivery, 73, 88
 in Continuous therapy, 71–74
 in mL/hr Only therapy, 85–89
 preparing for secondary delivery, 72, 87
 programming a secondary delivery, 73,
 88

- power indicator, 3
 - power sources
 - power check, 146
 - selecting, 18
 - powering on the pump, 24
 - prime, 3, 31
 - priming, 11
 - priming extension sets, 45
 - priming pump sets using pre-filled vials, 42
 - priming pump sets with a drip chamber, 38
 - priming pump sets with a drip chamber using the PURGE Key, 40
 - priming pump sets without a drip chamber using the PURGE Key, 36
 - pump sets without a drip chamber, 34
 - See also* prime; purge
 - Print, 109
 - programmable ranges
 - in Continuous therapy, 75
 - in Intermittent therapy, 67
 - in mL/hr Only therapy, 90
 - in Pain Management therapy, 60
 - in TPN therapy, 51
 - in Variable Time therapy, 95
 - in Weight Dosed therapy, 81
 - Programming Menu, 8
 - accessing the Programming Menu, 26
 - Programming Worksheet
 - Continuous, 75
 - Intermittent, 67
 - mL/hr Only, 90
 - Pain Management, 60
 - TPN, 51
 - Variable Time, 95
 - Weight Dosed, 81
 - proximal end, 30, 147
 - proximal occlusion
 - alarm messages, 133
 - occlusion sensitivity, 173
 - test of, 148
 - pump components, 4
 - pump features, 6
 - pump layout, 2–3
 - pump set, 12–16, 29, 167–169
 - GemStar Piggyback Pump Set, 69
 - pump components, 4
 - See also* pump settings
 - pump settings, 3, 98
 - air sensitivity, 107
 - change pump settings, 106
 - occlusion sensor, 107
 - sound level, 108
 - pump specifications, 171–173
 - PURGE key, 3
 - priming pump sets with a drip chamber using the PURGE Key, 40
 - priming pump sets without a drip chamber using the PURGE Key, 36
- R**
- Rate Change Lock, 58, 104, 122, 123
 - lock level functions, 105
 - locking the pump, 122
 - unlocking the pump, 124
 - rate titration
 - in Continuous therapy, 71
 - in mL/hr Only therapy, 85
 - in Weight Dosed therapy, 79
 - releasing the cassette, 33
 - RUN mode, 9
- S**
- SILENCE key, 3
 - Speed Protocols, 6, 98, 109
 - assign, review/recall, delete, 110–111

Index

- START key, 3
- starting a new container, 29
- STOP key, 3
- STOP mode, 9
- T**
 - taper down, 47, 48, 50
 - taper up, 47
 - technical assistance, 190
 - Therapies, 5
 - Continuous, 69–75
 - Intermittent, 61–67
 - mL/hr Only, 83–90
 - Pain Management, 53–60
 - Total Parenteral Nutrition (TPN), 47–51
 - Variable Time, 91–96
 - Weight Dosed, 77
 - Therapy Lock, 104
 - lock level functions, 105
 - Locking Sequence Instructions, 122
 - unlocking the pump, 124
 - Titrate Rates, 5
 - Titrate VTBI, 5
 - Total Parenteral Nutrition. *See* TPN therapy
 - TPN therapy, 5, 47–51
 - Auto Taper Down, 50
 - changing a TPN program, 49
 - features and specifications, 47
 - programmable ranges, 51
 - programming of, 48
 - Programming Worksheet, 51
 - troubleshooting, 127–138
 - alarm messages, 128–134
 - alert messages, 135
 - other display messages, 138
- U**
 - unlocking the pump, 124
 - UP/DOWN key, 28
- V**
 - Variable Time therapy, 5, 91–96
 - changing a Variable Time program, 94
 - delivery interruptions, 93
 - features and specifications, 91
 - programmable ranges, 95
 - programming of, 91
 - Programming Worksheet, 95
 - Volume To Be Infused. *See* VTBI
 - VTBI, 70, 78, 84
 - piggyback VTBI, 71–73, 86–88
 - VTBI Titration, 80, 86
- W**
 - warnings and cautions, 11
 - Weight Dosed therapy, 5, 77
 - Auto KVO, 79
 - changing a Weight Dosed program, 80
 - features and specifications, 77
 - programmable ranges, 81
 - programming of, 78
 - programming worksheet, 81
 - rate titration, 79
 - VTBI Titration, 80
- Y**
 - YES/ENTER key, 3, 28, 47, 54, 61, 69, 78, 84, 91

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician or other licensed practitioner.

DANGER: Possible explosion hazard exists if used in the presence of flammable anesthetics.



Authorized Representative:
Hospira
Donegal Town
Co. Donegal, Ireland



The 'C' and 'US' indicators adjacent to the CSA Mark signify that the product has been evaluated to the applicable CSA and ANSI/UL Standards, for use in Canada and the U.S., respectively. This 'US' indicator includes products eligible to bear the 'NRTL' indicator. NRTL, i.e. National Recognized Testing Laboratory, is a designation granted by the U.S. Occupational Safety and Health Administration (OSHA) to laboratories which have been recognized to perform certification to U.S. Standards.



Type BF

IPX1

Equipment providing a degree of protection higher than that for Type B equipment against electrical shock particularly regarding allowable leakage currents and having an F-Type applied part.

Protected against dripping water.



Attention: Consult accompanying documents.

GemStar, LifeShield, HEMA, and CAIR are registered trademarks of Hospira, Inc. Lucr Lock is trademark of Hospira, Inc., Super Edisonite, Seiko, Vesphene II se, Mann-Klenz, Formula C, Teflon, Velcro, CLAVE, and Microsoft, Windows, and Windows NT are not trademarks of Hospira.

This document and the subject matter disclosed herein are proprietary information. Hospira retains all the exclusive rights of dissemination, reproduction, manufacture and sale. Any party using this documentation accepts it in confidence, and agrees not to duplicate it, in whole or in part, nor disclose it to others without the written consent of Hospira.

© 2005 Hospira, Inc.
Lake Forest, Illinois 60045, USA
Printed in USA

Hospira GemStar Manual 430-00897-A02 (Rev. 02/05)

